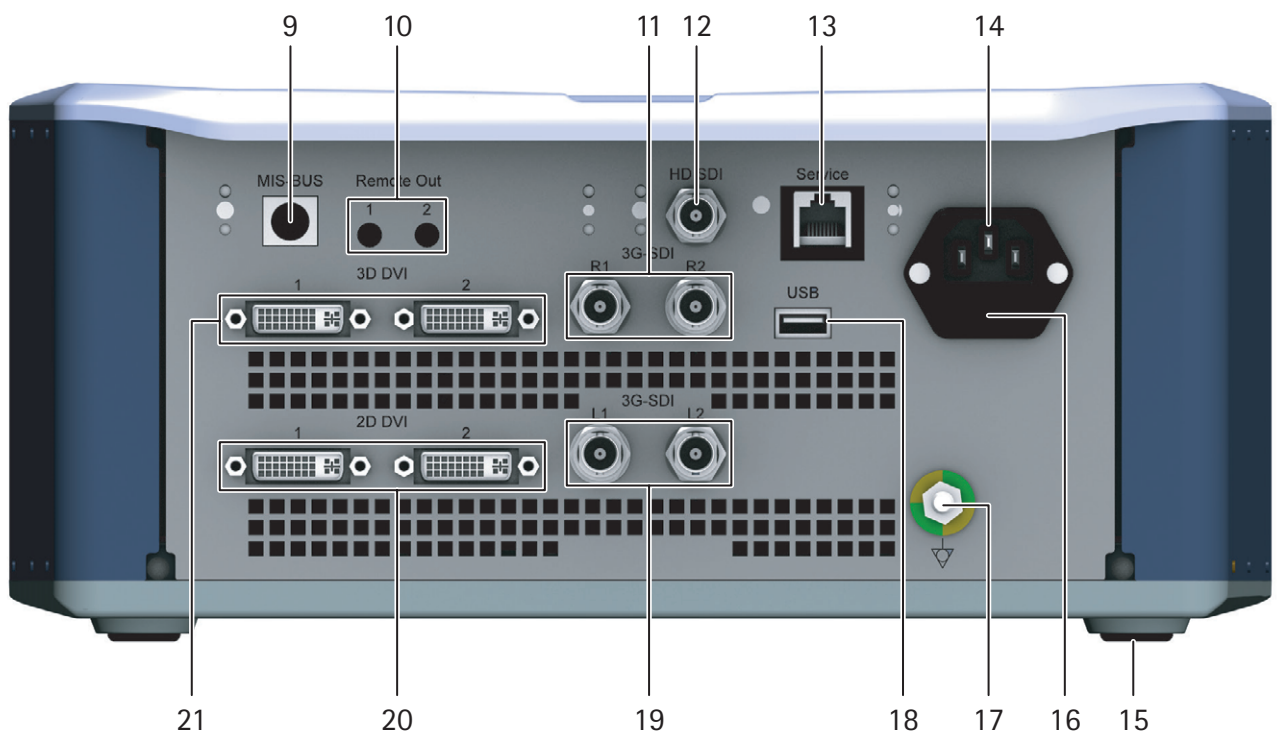
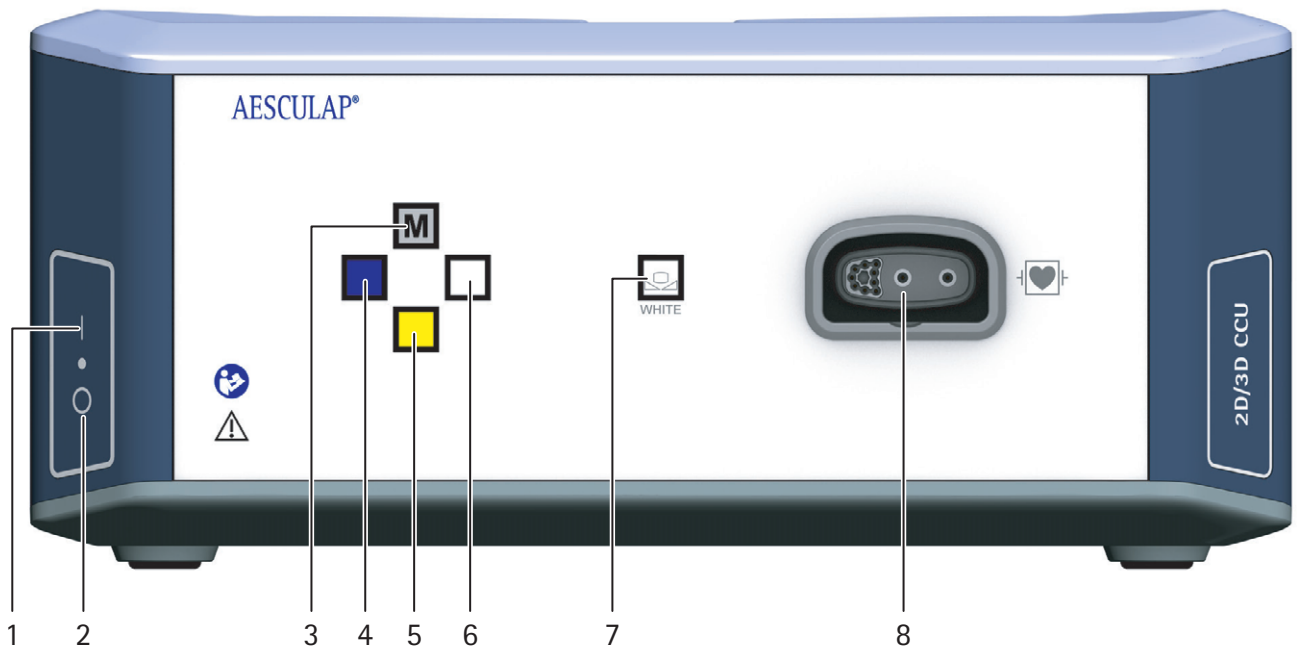
Aesculap Endoscopic Technology

- GB** Instructions for use/Technical description
USA EinsteinVision 3.0 camera control unit PV630
- D** Gebrauchsanweisung/Technische Beschreibung
EinsteinVision 3.0 Kamerakontrolleinheit PV630
- F** Mode d'emploi/Description technique
Unité de contrôle de caméra EinsteinVision 3.0 PV630
- E** Instrucciones de manejo/Descripción técnica
Unidad de control de la cámara EinsteinVision 3.0 PV630
- I** Istruzioni per l'uso/Descrizione tecnica
Unità di controllo della telecamera EinsteinVision 3.0 PV630
- P** Instruções de utilização/Descrição técnica
Unidade de controlo de câmara EinsteinVision 3.0 PV630
- NL** Gebruiksaanwijzing/Technische beschrijving
EinsteinVision 3.0 cameracontrole-unit PV630
- DK** Brugsanvisning/Teknisk beskrivelse
EinsteinVision 3.0 Kamerakontrolenhed PV630
- S** Bruksanvisning/Teknisk beskrivning
EinsteinVision 3.0 kamerakontrollenhet PV630
- FIN** Käyttöohje/Tekninen kuvaus
EinsteinVision 3.0 kameran ohjausyksikkö PV630
- RUS** Инструкция по применению/Техническое описание
блок управления камерой EinsteinVision 3.0 PV630
- CZ** Návod k použití/Technický popis
EinsteinVision 3.0 kamerová kontrolní jednotka PV630
- PL** Instrukcja użytkowania/Opis techniczny
Zespół sterujący kamery EinsteinVision 3.0 PV630
- SK** Návod na použitie/Technický opis
EinsteinVision 3.0 kamerová kontrolná jednotka PV630
- TR** Kullanım Kılavuzu/Teknik açıklama
EinsteinVision 3.0 kamera kontrol birimi PV630
- KR** 사용 설명서 / 기술 설명
EinsteinVision 3.0 카메라 컨트롤 유닛 PV630

Aesculap®

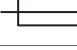
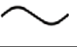
EinsteinVision 3.0 camera control unit PV630

Legend

- 1 Power ON switch
- 2 Power OFF switch
- 3 Short press: Image rotation 180° / Long press: Open menu (function in menu: up)
- 4 Short press: Image documentation / Long press: Video documentation (function in menu: left)
- 5 Short press: Light source on / Long press: Light source on/off (function in menu: down)
- 6 Short press: Zoom / Long press: White balance (function in menu: right)
- 7 White balance
- 8 Connector for camera head
- 9 Connector for MIS-Bus connection to LED light source OP950
- 10 Remote connectors (3.5 mm latches) for external documentation system
- 11 3G SDI outlets (1080p) for 2D or 3D monitors (Channel R1, R2)
- 12 HD-SDI output (1080i) for external documentation system
- 13 Connector for service PC
- 14 Power cord connector
- 15 Device foot
- 16 Fuse holder
- 17 Equipotential bonding connection
- 18 Connector Service-Dongle
- 19 3G SDI outlets (1080p) for 2D or 3D monitors (Channel L1, L2)
- 20 DVI outputs (1080p) for 2D monitors, left channel
- 21 3D DVI outlets (1080p) for 3D monitors (operating with 2D camera head: 2D DVI outlet).

Symbols on product and packages

Symbol	Explanation
	Signal word: WARNUNG Indicates a hazard. If not avoided, the hazard can result in death or serious injury.
	Signal word: CAUTION Indicates a potential hazard. If not avoided, this hazard may result in injury and/or damage to the product.
	Caution (IEC/EN/DIN 60601-1 3rd edition) / Attention, note the accompanying documents (IEC/EN/DIN 60601-1 2nd edition)
	Follow the instructions for use

	Defibrillation protected application part type CF pursuant to IEC/EN/DIN 60601-1
	Equipotential bonding
	Permissible storage temperature
	Permissible relative air humidity during storage
	Permissible atmospheric pressure during storage
	Article number
	Serial number
	Manufacturer
	Date of manufacture
	Fuse
	Separate collection for WEEE (waste of electrical and electronic equipment)
	Hot temperature
	Not permitted for use in a magnetic resonance environment
	Alternating current
	US federal law requires that the product may only be supplied to a physician or a person acting on the physician's behalf

Contents

1.	Applicable to	3
2.	Safe handling	4
3.	Product description	6
3.1	Scope of supply	6
3.2	Components required for operation	6
3.3	Intended use	6
3.4	Operating principle	7
4.	Preparation and setup	12
4.1	Electromagnetic tolerance	13
4.2	Combination with medical electrical equipment	13
4.3	Setup	13
4.4	First use	14
5.	Working with the camera control unit PV630	15
5.1	System set-up	15
5.2	Function checks	16
5.3	Configuration	16
5.4	Safe operation	17
6.	Reprocessing procedure	20
6.1	General safety notes	20
6.2	General information	20
6.3	Cleaning/disinfection	20
6.4	Inspection, maintenance and checks	21
6.5	Storage	21
7.	Repairs	21
8.	Maintenance	22
9.	Troubleshooting list	22
9.1	User maintenance	23
9.2	Repairs	23
10.	Technical Service	23
11.	Accessories/Spare parts	24
12.	Technical data	24
12.1	Ambient conditions	25
13.	Disposal	25
14.	Distributor in the US/Contact in Canada for product information and complaints	25

1. Applicable to

These instructions for use apply to the following device:

Article name: **EV3.0 camera control unit**

Art. no.: **PV630**

These instructions for use are an integral component of the device and contain all the information required by users and operators for safe and proper use.

Target group

These instructions for use are intended for physicians, medical assistants, medical technicians and employees of sterile services who are entrusted with the installation, operation, maintenance and reprocessing of the device.

Using and storing this document

These instructions for use must be stored in a defined location so that it may be accessed at all times by the target group.

In the event of the sale of this device or its relocation, this document must be handed over to the new owner.

Supplementary documents

The electromagnetic compatibility (EMC) requirements are described in the separate brochure enclosed. The rules and information contained therein must be observed during installation and operation of this device.

The instructions for use of the camera head used and instructions for use of all other devices used must be observed to ensure the safe use of the device.

- For item-specific instructions for use and information regarding material compatibility, see Aesculap Extranet at www.aesculap-extra.net

2. Safe handling

Risk to user and patients if the instructions, warnings and precautions are not followed!

- ▶ Use the product only according to these instructions for use.
- ▶ Check the product is in good working order before use.
- ▶ Carry out function checks before using the equipment.
- ▶ Do not use the product if irregularities are observed during inspections and checks.

Risk of injury due to electric shock!

- ▶ Ensure that the rear of the device and access to the mains plug is freely accessible at all times.
- ▶ When installing a medical electrical system, there is a risk of fire, short-circuit or electric shock. Installation must be performed by qualified personnel only.
- ▶ For the combination of electrical devices, follow annex 1 of IEC/EN/DIN 60601-1. Non-medical devices which adhere to the applicable IEC safety standards may only be connected via a medical isolation transformer. Do not connect any additional non-medical devices to a medical electrical system.
- ▶ Signal lines from devices with a functional connection, which are connected to different branches of the mains power supply, must be galvanically isolated at both ends.
- ▶ Connect the devices only to a power supply with a protective earth conductor.
- ▶ After installation of a medical electrical system, perform a test pursuant to IEC/EN/DIN 62353.

Risk to patient if the device power fails!

- ▶ Use the device only with an uninterruptible power supply. To guarantee a continuous power supply a medical emergency supply (UPS) is recommended.

Risk to patient due to incorrect application!

- ▶ Use the medical device only after instruction by the manufacturer or by an authorized person.
- ▶ The instructions for use enclosed with the individual EinsteinVision components and all products used (e.g., high frequency surgery) must be observed.
- ▶ Endoscopic procedures may only be carried out by specialists who have the necessary medical training, knowledge and experience.

Risk to patient and user due to premature wear!

- ▶ Handle and maintain the medical device properly.
- ▶ Use the medical device only in accordance with the intended purpose.

Possible restricted function if third party devices are used!

- ▶ Use the product with the recommended components and accessories.
- ▶ The full functionality of the product can only be guaranteed if the recommended accessory components are used.

Risk to the patient if the product is not used gently/a damaged product is used!

- ▶ Handle the product with appropriate care.
- ▶ Do not use the product if it has been subjected to strong mechanical stress or if it has been dropped, and send the product to the manufacturer or an authorized repair center for inspection.

Risk of damage to eyes from UV radiation!

- ▶ Do not use polarized 3D glasses as sun glasses.

Device failure due to incorrect storage and usage conditions!

- ▶ Store and operate the product only within the specified environmental conditions.

Risk to patient if the live image fails!

- Shock from defibrillation can lead to failure of the live-image. It can take up to 3 seconds for recovery of the live-image. Restart of the Camera Control Unit can possibly be necessary by switching off and on.
- ▶ Check live-image after defibrillation.

Improper operation caused by decreased electromagnetic immunity or increased electromagnetic emissions!

- ▶ Use product only with original or by the manufacturer specified transducers, cables and accessories.

Note

Additional equipment and/or peripheral devices which are connected to the interfaces of the product must show their relevant specifications (such as IEC/EN/DIN 60601-1).

Note

The medical electrical device may not be modified in any way.

Note

To guarantee optimal function of the EinsteinVision components, use in controlled environmental conditions is recommended (air-conditioned operating room).

Note

Always handle the medical device very carefully as it contains sensitive optic, mechanic and electronic components.

Note

Do not store the polarized 3D glasses in an environment with high temperatures, e.g., near a radiator.

Note

Ensure that all devices operated in the vicinity meet their relevant EMV requirements.

Note

The EinsteinVision components are used to visualize the inside of the body during minimally invasive procedures. Do not use the components for diagnostic purposes. This applies in particular with the use of image optimization algorithms.

Note

Images and videos recorded are for documentation purposes only. Do not use for diagnosis or findings, as the image quality is reduced in the event of data compression.

Note

Only operate the camera control unit when it is closed. Only operate the device in a horizontal position.

Note

Before putting into operation, check the compatibility of all components and accessories using the accessories list.

Note

All accessories and spare parts must only be procured from the manufacturer.

Note

Spare parts and fuses must only be replaced by authorized personnel.

Note

The connection of electrical devices to a multiple socket creates a system and can lead to a reduced level of safety.

Note

To disconnect the device completely from power, pull out the plug.

Note

Only operate the device with the original cables (included in scope of delivery) or with specified cables/components.

- ▶ To prevent damage caused by improper setup or operation, and in order not to compromise warranty and manufacturer liability:
 - Use the product only according to these instructions for use.
 - Follow the safety and maintenance instructions.
 - Only combine Aesculap products with each other.
 - Follow the application advisories acc. to standard, see Extracts from relevant standards.
- ▶ Ensure that the product and its accessories are operated and used only by persons with the requisite training, knowledge, or experience.
- ▶ Keep the instructions for use accessible for the user.
- ▶ Always adhere to applicable standards.

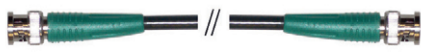
CAUTION

Federal law restricts this device to sale by, or on order of a physician!

3. Product description

3.1 Scope of supply

Designation	Art. no.
EV3.0 camera control unit	PV630
BNC cable (red, with ferrite bead) length 3.0 m	TA014803
BNC cable (green, with ferrite bead) length 3.0 m	TA014804
MIS Bus cable length 0.75 m	OP942
Instructions for use	TA014615
EMC brochure	TA014630
Open source license brochure	TPO003-121-01

Cables included	
BNC cable (red)	
BNC cable (green)	
MIS-Bus cable	

3.2 Components required for operation

The EV3.0 Camera Control Unit PV630 is intended for use with:

- PV631/PV632 EV3.0 camera head 10 mm 0°/30°

In order to guarantee full functionality, the joint use of the following components is recommended:

- OP950 LED light source
- PV646 3D monitor

For the 2D visualization the following components can also be connected:

- PV481 Full HD CMOS camera head with oscillating lens
- PV482 Full HD CMOS camera head with zoom lens
- PV485 Full HD 3CMOS camera head with zoom lens

If the recommended components are not used, the following errors/complications may occur, among others:

- ▶ Incorrect color display
- ▶ No/incorrect display of the camera menu
- ▶ No/incorrect 3D display
- ▶ Restricted operating function from the camera head
- ▶ Increased fogging on the outlet of the optical lens
- ▶ No switching between 2D and 3D display
- ▶ Damage to the sterile cover

Note

Furthermore, if the product is used with third party devices, unrestricted responsibility is passed to the system configurator or medical institution.

3.3 Intended use

The EV3.0 Camera Control Unit serves the 3D and 2D visualization of the interior of the body during minimally invasive surgical procedures. The application in combination with a camera head serves 3D and 2D visualization of the intracorporeal surgical area during endoscopic diagnostic and surgical procedures. The EV3.0 Camera Control Unit is only used for visualizing and not for diagnosis. Prefer additionally the intended use of the relevant camera heads and components

Indications

There are no medical indications associated with the use of the camera if separated from the endoscopic system.

Contraindications

The use of the EinsteinVision components and their accessories is contraindicated if endoscopic procedures are contraindicated for any reason. As is appropriate with any surgical procedure, consideration must be given to patient size and workspace volume when using the EinsteinVision components.

Depending on the patient's disease, there may exist contraindications that rely on the general condition of the patient or the specific disease pattern.

The decision to perform an endoscopic procedure rests with the responsible surgeon and should be made on the basis of an individual risk-benefit analysis.

3.4 Operating principle

CAUTION

Risk of intolerances such as vertigo, headache and nausea due to 3D visualization!

- ▶ Before first use, check tolerance. If intolerances occur during use, switch over to 2D imaging.

WARNING

Risk to patient if the image shown is misinterpreted!

- ▶ When using image optimization algorithms (e.g. Red Enhancement and Smoke Reduction) be aware that the image shown is not the original image.
- ▶ During use of these, never rely solely on the optimized image.

Note

The EV3.0 camera head can only be used for 3D visualization if the LED light source OP950 is used in addition to the EV3.0 camera control unit. For this, the camera control unit and light source must be correctly connected via a MIS-Bus cable.

This product is a camera control unit that is intended for use in medical endoscopy.

The camera control unit can be used with different camera heads. Together, the camera control unit and camera head comprise the camera. In combination with a suitable monitor, the camera can provide either two-dimensional or three-dimensional images. For three-dimensional visualization, the user requires polarized 3D glasses.

In addition to the general setting options, the camera offers different endoscopic application profiles. These profiles are preconfigured to offer an optimized image display; however, the profiles can also be modified and customized for your specific needs.

The current image display can also have effects added, for example to distinguish tissue types better from each other or to get a clearer picture in case of smoke production.

For quick access to frequently used functions, each profile can select favorites, which can also be configured individually in the same way as the profiles.

Current image settings and status indicators are shown on the screen in the form of symbols and messages.

The general setting options include language selection for the configuration menu and monitor displays, and the camera head button settings. The factory setting is English.

Main menu on screen

The main menu on the screen opens with a long press (≥ 2 seconds) of the menu button (3) on the camera control unit or connected camera head.

The main menu is used to operate and configure the device. The settings options are reached by selecting the relevant symbols.

The functions represented by the symbols are explained below.

	Information on the currently connected camera head. No further setting options
	Profiles (application profiles and configurable user profiles)
	Favorites
	Brightness
	Digital zoom
	Contrast
	Edge enhancement
	Image rotation 180° (only available in combination with EV3.0 camera head 30°)
	Switch to 2D/3D (only available in connection with the EV3.0 camera head)
	Access to the selectable image effects
	System settings, restoring the default settings

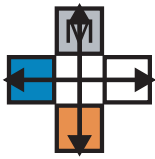
Aesculap®

EinsteinVision 3.0 camera control unit PV630

Navigation in the main menu

With a long press of the button (≥ 2 seconds), the main menu opens and is displayed on the screen.

Use the relevant buttons to navigate to the desired item in the menu or submenu.

Profiles

Profiles are preinstalled camera settings, optimized for the application.

The following four standard profiles are available when using 3D camera heads (PV631, PV632):

- LAP General Surgery
- LAP Uro
- LAP Gyn
- Cardio/Thoracic

The following 5 standard profiles are available when using 2D camera heads (PV48x):

- LAP General Surgery
- LAP Uro
- LAP Gyn
- Neuro
- Cardio/Thoracic

Image setting

The image display can be changed manually for the set profile. These changes are retained until the next time the type of camera head or profile is changed.

The following setting options are available:

- Brightness (-5 to +5 / increments of: 1)
- Digital zoom (1 to 1.8 / increments of: 0.2)
- Contrast (-5 to +5 / increments of: 1)
- Edge enhancement (-5 to +5 / increments of: 1)
- Image Rotation 180° (on/off)
- Toggle 2D/3D (on/off)

Effects

Effects can be added to the set profile, which can improve the displayed image under certain conditions of use.

The following effects are available:

- **Red Enhancement (RE):** Red Enhancement intensifies red color tones and enables a better visual distinction of different structures in environments with little differences in color, e.g. to accentuate the colors of vessels.
- **Smoke Reduction (SR):** Smoke Reduction produces a clearer picture in smoky environments and therefore a better view (e.g. during coagulation with high frequency surgical devices).

- **Picture out of Picture (PoP):** Picture out of Picture is not an image effect in itself but it splits the screen into two equal image areas for the same camera image. The reduced-size original image from the camera is displayed on the left and the reduced-size image with an effect is displayed on the right. With PoP, you can compare the original image with the image with an applied effect and choose between the different effects.
- **Status indicators on screen:** The currently activated image settings are displayed on the top right of the screen in the status area.

Note

The Red Enhancement and Smoke Reduction effects can be combined with each other.

The symbols displayed when the corresponding function is activated are explained below.

	AUTO light intensity control deactivated
	Zoom function activated
	Image rotation 180° activated (only available in connection with EV3.0 camera head 30°)
	2D mode activated (only available in connection with EV3.0 camera head)
	Picture out of Picture (PoP) activated
	Red Enhancement (RE) activated
	Smoke Reduction (SR) activated

Favorites

Depending on the profile and connected camera head, the device offers a pre-defined selection of functions as favorites for quick access.

The favorites in the standard profiles cannot be changed. If user profiles are created, the favorites can then also be individually defined.

Control buttons on the front of the device

The 4 control buttons on the front of the device in the shape of a cross have different functions depending on the operating mode.

The assignment of button functions is always the same on the front of the device and on the camera head.

Different functions can be assigned for a short press of the button (<2seconds) and a long press of the button (≥2 seconds).

Button assignment in live mode

The button assignment in live mode depends on whether a standard profile or user profile is set.

If a user profile is set, the buttons can be assigned individually, however the assignment for entering the menu cannot be changed.

In standard profile, the buttons are assigned as follows.

	Short press of button	Long press of button
	Image rotation 180° (only available in connection with EV3.0 camera head 30°)	Open menu
	Zoom	White balance
	Light source on	Light source on/off
	Single frame image	Video recording start/stop

Note

The camera control unit does not have an internal memory. Video and image recordings are only possible if an external documentation device is connected.

Button assignment when on-screen display is active

When on-screen display is active, the buttons can be used to navigate the menu, until the on-screen display function is left.

	Short press of button	Long press of button
	Up / control +	Leave the menu
	To the right / level deeper / save controller settings	Leave the menu
	Down / control -	Leave the menu
	To the left / level back / save controller settings	Leave the menu

Button assignment when on-screen keyboard is active

When the on-screen keyboard is active, the buttons can be used to navigate on the keyboard, until the on-screen keyboard function is exited.

	Short press of button	Long press of button
	Up	-
	Right	Confirm/enter
	Down	-
	Left	-

Aesculap®

EinsteinVision 3.0 camera control unit PV630

Display of faults

If the functionality of the device is restricted, faults are displayed in the top left area of the screen.

If one of the faults described below occurs, the main menu closes automatically and it cannot be accessed while the fault persists.

Therefore existing settings, for example an image rotation or activated effect, cannot be deactivated or undone while the fault persists.

Malfunction	Detection	Cause	Remedy
Overheating	On-screen message: "Overheating! Camera may switch off automatically."	The device features an internal temperature sensor.	Check ventilation slits on the bottom of the device. Set up the device so that sufficient ventilation is guaranteed. Check if light source is placed on top of the camera control unit. If the message appears repeatedly: contact B. Braun/Aesculap partner or Aesculap Technical Service.
Interrupted connection to light source	On-screen message: "Malfunction! Connect compatible light source OP950."	After connecting a 3D camera head, the device can only operate if the light source OP950 is properly connected and switched on.	Check MIS bus connection between the camera and light source. If the message appears despite correctly connected OP950: contact B. Braun/Aesculap partner or Aesculap Technical Service.
	On-screen message: "AUTO light quantity control deactivated".	While a 2D camera head is connected to the device, the light source OP950 is not properly connected or not active.	Reactivate AUTO light quantity control using the menu.
Failure of the 3D function	On-screen message: "Malfunction! One of the camera channels is out of service."	The 3D function works with two camera modules. If one camera module fails, the device switches to two-dimensional display.	Restart the device. If the message appears repeatedly: contact B. Braun/Aesculap partner or Aesculap Technical Service.
Malfunction of the camera control unit	The five control buttons on the front of the device flash quickly.	In case of a malfunction of the device, the connected monitor no longer displays an image.	Restart the device. If the message appears repeatedly: contact B. Braun/Aesculap partner or Aesculap Technical Service.

Status messages

Status messages are superimposed in the right upper area of the screen if there are changes to device status.

Status message on screen	Status/Status Change	Remedy
3D camera head 0° 3D camera head 30° 2D camera head	A compatible camera head is attached to the camera control unit and recognized.	–
Connect the camera head to the control unit.	No compatible camera head is connected to the control unit.	Check the plug connection between the camera head and the camera control unit. If the message appears repeatedly: contact B. Braun/Aesculap partner or Aesculap Technical Service.
Connect the compatible camera head to the control unit.	No compatible camera head is connected to the control unit.	Check compatibility between camera head and camera control unit. If the message appears repeatedly: contact B. Braun/Aesculap partner or Aesculap Technical Service.
Light source ON Light source OFF	The light source is switched on or off.	–
Check the optical fiber connection.	The optical fiber is not plugged at all or not correctly to the light source OP950.	Check the plug connection between the optical fiber connection and the light source. If the message appears repeatedly: contact B. Braun/Aesculap partner or Aesculap Technical Service.
Symbols of active image effects are not displayed.	The display of the status symbols in the user menu is deactivated before restarting the camera control unit. The status symbols are only displayed for a short time during activation/deactivation of the relevant function.	Reactivate status symbols in the menu.
White balancing successful White balancing failed	After activation of the white balance there is a display of whether this has been done successfully.	–
Tip heating defective	The tip heating of the camera head is out of order and the tip of the endoscope is not heated.	Restart the device. If the message appears repeatedly: contact B. Braun/Aesculap partner or Aesculap Technical Service.
Camera head buttons out of order due to magnetic field	An external magnetic field disrupts the function of the buttons on the camera head. Operation using the buttons on the camera control unit is still possible.	Restart the device. If the message appears repeatedly: contact B. Braun/Aesculap partner or Aesculap Technical Service.

Aesculap®

EinsteinVision 3.0 camera control unit PV630

4. Preparation and setup

Non-compliance with these rules will result in complete exclusion of liability on the part of the manufacturer.

- ▶ When setting up and operating the product, adhere to
 - The national installation and operator regulations
 - National regulations on fire and explosion protection

Note

The device is authorized for operation in hospitals.

Note

For the safety of patients and users it is essential that the mains power cord and, especially, the protective earth connection are intact. In many cases defective or missing protective earth connections are not registered immediately.

- ▶ Connect the device via the potential equalization terminal at the rear panel of the device to the potential equalization system of the room used for medical purposes.

WARNING

Risk to patient due to improperly conducted earth current due to lack of or defective earth!

- ▶ Do not touch the device and patient at the same time.

WARNING

Risk to patient due to failure of device!

- ▶ Keep an operationally ready replacement device on hand. If necessary change to conventional surgical methods.

WARNING

Danger of explosion from a device that is improperly set up!

- ▶ Make certain that the power plug is connected to the power supply outside of any areas where there is a risk of explosion.
- ▶ Do not use the product in potentially explosive areas or close to highly flammable or explosive gases (e.g. oxygen, anesthetic gases).

WARNING

Risk of fire or electric shock if moisture penetrates the device!

- ▶ Ensure that the device is set up protected from water drips and splashes.
- ▶ Do not immerse the housing in any fluid.
- ▶ Do not expose the device to rain or any moisture.
- ▶ Do not place objects on the product that contain any liquids.

WARNING

Interference of high frequency energy with the device!

- ▶ Do not use any mobile or portable equipment that emit high-frequency energy (e.g. mobile phones, cell phones, GSM phones) in the vicinity of the product.
- ▶ Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 100 to any part of the EV3.0 visualization system, including cables specified by the manufacturer. Otherwise, degradation of the performance of this equipment could result.

CAUTION

Risk to patient through poor/lack of view!

- ▶ Adjust the monitors and display elements so that they are highly visible to the user.

CAUTION

Risk of electric shock from incorrect installation of the device!

- ▶ First installation without power connection.

CAUTION

Risk to persons and risk of damage to equipment due to improper cable routing!

- ▶ Lay all cables and leads so that they do not present a tripping hazard.
- ▶ Do not place objects on the cables.

CAUTION

Risk to patient and user through condensation and short circuit!

- ▶ Before putting into operation make sure that all components used have enough time to adjust to the changed environmental conditions.

CAUTION

Malfunction of the device!

- ▶ For sufficient ventilation, set up the EV3.0 camera control unit so that the ventilation slots are always clear.
- ▶ Set up and operate the EV3.0 camera control unit horizontally.

Note

The equipotential plugs for all devices used must be connected with the equipotential tape IEC/EN/DIN 60601-1 or consistent with national standards).

Note

Unauthorized access to data or loss of data: Only use the EV3.0 camera control unit on secure networks.

Note

Ensure that the corresponding Interconnection Conditions are kept. Also, the relevant standards and the respective national deviations must be followed.

Personnel responsible for setting up or installing medical electrical equipment must possess the necessary qualifications and be familiar with the safety regulations in place at the installation location and the applicable official regulations.

4.1 Electromagnetic tolerance

Medical electrical equipment is subject to stringent electromagnetic compatibility (EMC) requirements.

Despite the device's high interference immunity and low emitted interference, compliance with EMC-related requirements is necessary with regard to installation, installation location, and ambient conditions.

It is therefore important that you observe and follow the related safety messages in this document and the EMC brochure included with the device.

4.2 Combination with medical electrical equipment

The product can be used with components from other manufacturers if they fulfill the requirements of IEC/EN/DIN 60601-1 for the safety of medical devices.

It is the operator's responsibility to check and make sure that the system is and remains fully operational.

When using devices from different manufacturers and when operating an endoscope and/or endoscopic accessories together with medical electrical equipment, it must be ensured that the applied part is properly isolated: Type CF, defibrillation-proof.

4.3 Setup

- ▶ Set up the device.
- ▶ When setting up the device, ensure that:
 - It is positioned on an even, non-slip surface that is able to support its weight
 - It is installed outside of the patient environment and the sterile area on a sufficiently stable support
 - It is protected from water drips and splashes

- It will not be subject to vibration during operation
 - The ventilation slots are not obstructed
- ▶ Connect the **equipotential bonding connection** on the rear of the device to the POAG connection at the installation site.

Stacking of units

- ▶ Do not exceed the maximum stack height of 450 mm.
- ▶ Place the units in a stable position.
- ▶ Aesculap devices must be stacked on top of each other facing identically.
- ▶ Never stack in an offset arrangement.

Connecting 3D monitors

Follow the instructions for use of the monitors used.

- ▶ Connect a 3D-DVI-connector to a 3D monitor or connect two 3G SDI connectors (R1/L1, R2/L2) to a 3D monitor.

Note

Either the DVI or the SDI cables are to be used for operation. The SDI cable comes standard. Both variants are represented in the following figures.

- ▶ When using the SDI variant: Place the ferrite beads on the SDI-cable close to the monitor.

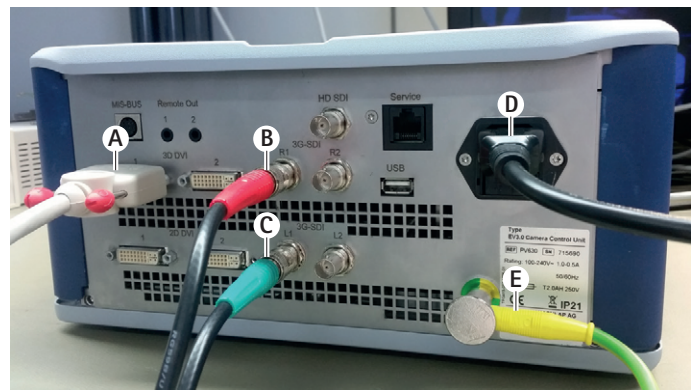
Fig. 1

Legend

- A DVI cable at 3D DVI output
- B SDI cable red at 3G SDI output (right channel)
- C SDI cable green at 3G SDI output (left channel)
- D Power cord
- E Equipotential bonding for connection of the camera control unit to the equipotential bonding strip

Aesculap®

EinsteinVision 3.0 camera control unit PV630

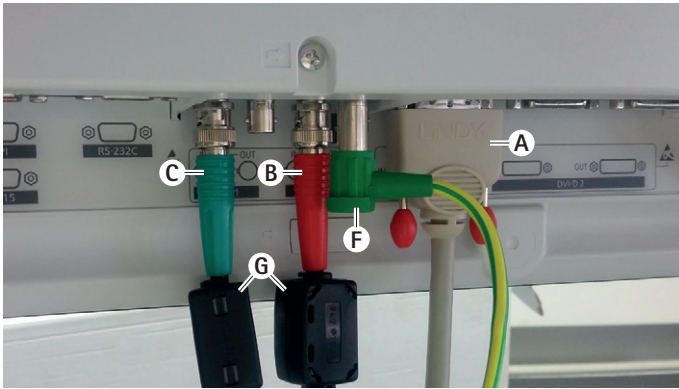
Fig. 2

Legend

- A DVI cable at DVI input (3D DVI)
- B SDI cable red at 3G SDI input (right channel)
- C SDI cable green at 3G SDI input (left channel)
- F Equipotential bonding for connection of the monitor to the equipotential bonding strip
- G Ferrite beads

Connecting 2D monitors

Follow the instructions for use of the monitors used.

- ▶ Connect a 2D- or 3D-DVI connection with a 2D monitor or a 3G SDI connection (R1/R2/L1/L2) with a 2D monitor.

Connecting a documentation system

In combination with an external documentation system (e.g., EDDY3D/EDDY), images and videos can be taken with the EV3.0 camera control unit.

Follow the instructions for use of the documentation system.

- ▶ Connect a remote connection (Remote 1 or Remote 2) of the camera control unit with the corresponding connection on the documentation system.

Connecting the light source OP950

The connected light source can be controlled using the buttons on the camera head.

- ▶ Connect MIS bus connection of the camera control unit and the light source.
- ▶ Ensure that the connectors engage in place when inserting them.

4.4 First use

WARNING

Risk of injury and/or product malfunction due to incorrect operation of the electromedical system!

- ▶ Adhere to the instructions for use of any medical device.

Requirement: Installation must have been completed.

- ▶ Connect the camera head.
- ▶ Switch on the device and all devices being used.
- ▶ If there is a focusing ring, focus the image.
- ▶ Set the white balance.
- ▶ Check that a good image is displayed on the monitor at working distance without vertical lines, color variations, or flickering. When using a 3D camera head, perform a function test with 3D polarization glasses.
- ▶ Ensure that the buttons on the camera head and camera control unit can be operated properly.
- ▶ If applicable, ensure that the light source can be switched on and off via the buttons on the camera head and camera control unit.

Voltage

The mains voltage must correspond to the voltage indicated on the type plate at the back of the unit.

5. Working with the camera control unit PV630

5.1 System set-up

WARNING

Risk of infection due to unsterile parts!

- ▶ Unsterile components must not enter the sterile area.
- ▶ Process products and accessories that are delivered unsterile before use and only use them with sterile accessories.

CAUTION

Risk to patient through addition of earth currents!

- ▶ If the product is used with medical electrical equipment and/or power-driven endoscope accessories, the leakage currents can multiply.

WARNING

Risk to patient and user through burns, sparks, or explosion!

- ▶ When using a high frequency surgical device during an endoscopic procedure, observe the safety instructions of the corresponding instructions for use.

Note

Before application of endoscopic high frequency surgery prepare the patient accordingly. Take measures to remove or avoid formation of combustible gases (e.g. gastro-intestinal tract/colonoscopy, urinary bladder/transurethral resection).

WARNING

Risk of injury due to improper use!

- ▶ Observe the intended use of the camera control unit.

CAUTION

Interference of electromagnetic emissions with the image quality (e.g. minor banding, minor color changes on the monitor image)!

- ▶ Check the image quality if used in combination with additional peripheral devices (e.g. monitor, documentation system).

Note

In combination with the EV3.0 camera control unit, the endoscope is classed as a type CF defibrillation-proof applied part.

Connecting the accessories

DANGER

Risk of injury due to unapproved configuration using additional components!

- ▶ For all applied components, ensure that their classification matches that of the application component (e.g. Typ CF defibrillation-proof) of the respective device.

Combinations of accessories that are not mentioned in the present instructions for use may only be employed if they are specifically intended for the respective application, and if they do not compromise the performance and safety characteristics of the products.

All devices connected to the interfaces must also be proven to fulfill the relevant IEC standards. IEC 60950 for data processing devices and IEC/EN/DIN 60601-1 for medical electrical devices).

All configurations must fulfill the fundamental standard IEC/EN/DIN 60601-1. The person connecting the devices with each other is responsible for the configuration and must ensure that the fundamental standard IEC/EN/DIN 60601-1 or relevant national standards.

- ▶ Please contact your B. Braun/Aesculap partner or the Aesculap Technical Service, see Technical Service.

Connecting the power supply

DANGER

Risk of death by electric shock!

- ▶ Connect the product only to a grounded power supply.

- ▶ Ensure that the mains voltage at the site of operation corresponds to that specified on the device nameplate.
- ▶ Connect the power cord to the mains connection socket on the rear of the device.
- ▶ Connect the power cord to the power supply.
- ▶ Lay the cable so that no-one can trip over it or get caught in it.

Visual inspection

Carry out a visual inspection prior to every procedure.

- ▶ Ensure that the device has been properly set up. Ensure that the equipotential bonding has been properly connected.
- ▶ Ensure that the housing is free from external damage.
- ▶ Ensure that the power cord and all other cables are not damaged.
- ▶ Ensure that the contacts in the connection sockets for the camera head are free from moisture and contamination.

5.2 Function checks

Risk of electric shock from defective cables or devices!

- ▶ Inspect the electrical equipment regularly.
- ▶ Replace loose, worn or defective cables/connections.

Risk to patient due to impaired or lack of 3D image impression!

- ▶ Check the fit of the polarized 3D glasses before use and if necessary use a suitable glasses strap or have them fitted by an optician (e.g. if loupe glasses are also to be used in combination with the 3D glasses).
- ▶ Do not store the polarized 3D glasses in an environment with high temperatures (e.g. near a radiator).
- ▶ Do not use worn, scratched, broken or damaged polarized 3D glasses.

Risk to patient and user from functional failure or electric shock!

- ▶ Use only dry components (e.g. plug connector for camera control unit, optical lens).

Risk to patient through improperly represented image!

- ▶ Prior to application and after every change of settings (e.g. image rotation or activation of algorithms) check the correct display of the live image. If necessary, perform white balance.

Note

Only perform the surgical procedure if all of the components being used are in a flawless condition.

Note

Before each use, after a power failure or any interruption, the functionality of all connected devices must be checked and all connections must be checked to see they are correct.

Note

The endoscope image must be focused, bright and clean at an appropriate working distance.

Note

Do not continue to use damaged products.

Carry out a function check prior to every procedure.

For 3D visualization, a direct alignment of the 3D monitor with the surgeon is required for a good display of the three-dimensional image.

- ▶ Connect the camera head to the camera control unit.
- ▶ Switch on the device and all devices being used.
- ▶ If there is a focusing ring, focus the image.
- ▶ Point the camera head to an object within the normal working distance.
- ▶ Set the white balance.
- ▶ Align the monitor to the surgeon.
- ▶ Put on the polarized 3D glasses and ensure that they fit well. If necessary, use a suitable glasses strap or have them fitted by an optician.
- ▶ Ensure that a good 3D live image is displayed on the monitor without vertical lines, color variations, or flickering.
- ▶ Ensure that the buttons on the camera head and camera control unit can be operated properly.
- ▶ If applicable, ensure that the light source can be switched on and off via the buttons on the camera head and camera control unit.

5.3 Configuration

System settings

The following settings can be applied under the menu item **System Settings**:

- Language: To select the language for the camera menu (default setting: English)
- 50 Hz/60 Hz: To select the image refresh rate
- Network setup: To create a network connection for the technical service
- Display status bar: To select whether symbols with information on the currently active image settings should be displayed in the status area of the screen.
- Acoustic signal Remote output (deactivated as standard): To select whether an acoustic signal should sound when an image is taken or when starting and ending a video recording.
- Restore default settings: To restore all device settings to the default settings

The following system information can also be displayed under this menu item:

- Software version of the camera control unit, connected camera head and connected light source
- Serial number of the camera control unit, connected camera head and connected light source

The submenu **Network Setup** can also be used to display an overview of the current network settings.

- ▶ Open menu and navigate to menu item **System Setup**.
- ▶ Navigate to the required item and apply the settings.

User profiles

In the **Profile** submenu up to three individual user profiles may be created for each type of camera head.

To create a user profile, an initial profile first needs to be selected, of which the image settings can then be changed and saved in the user profile.

In addition to the image settings, the following can also be set in the user profile:

- Activation/Deactivation of the AUTO light intensity control
- Button function assignment on the camera control unit and camera head
- Individual favorites menu

The user profile creation process ends with the assignment of a profile name for the created user profile using the on-screen keyboard. The profile name must not be more than 16 characters and it must be different from the standard profile names.

A created user profile can be changed or deleted when the relevant camera head type is connected.

- ▶ Open menu and navigate to menu item **Profile**.
- ▶ Navigate to **Create user profile** and select an initial profile.
- ▶ Apply the individual settings.
- ▶ Navigate to **Save and exit** and use the on-screen keyboard to enter a name for the user profile.
- ▶ With a long press of the white button (≥ 2 seconds), the user profile is saved under the entered name.
- ▶ To change, go to **Change user profile** and navigate to the entry that needs changing.
- ▶ To delete, go to **Delete user profile** and navigate to the user profile that needs deleting.

Favorites

For the standard profiles, the favorites are pre-defined and cannot be changed. For the individually configurable user profiles, the favorites can be selected individually and changed if necessary.

Up to 10 favorites can be defined if necessary.

- ▶ Open menu and navigate to menu item **Profile**.
- ▶ Navigate to **Create user profile** and then to **Configure favorites menu**.
- ▶ To change the settings, go to **Change user profile** and then **Configure favorites menu**.

Button function assignment for control buttons

For the standard profiles, the button functions are pre-defined and cannot be changed.

For the individually configurable user profiles, the button functions can be selected individually and changed if necessary.

If the button functions are changed, the control buttons on the front of the device and on the camera head automatically have the same functions.

- ▶ Open menu and navigate to menu item **Profile**.
- ▶ Navigate to **Create user profile** and then to **Button function assignment**.
- ▶ Select a function for a long press and short press for each button.
- ▶ To change the settings, go to **Change user profile** and then **Button function assignment**.

5.4 Safe operation

Excessive temperatures in combination with light sources

Risk of injury due to excessive temperature!

- ▶ Do not touch the light guide connector, distal end of the endoscope or the tip of the sterile cover during use or immediately after use.

Light sources, especially high-power light sources, emit large amounts of light energy and thermal energy.

The light guide connector, distal end of the endoscope and the tip of the sterile cover can become extremely hot.

Risks from use of light sources:

- Irreversible tissue damage or unwanted coagulation to the patient or user
- Burns or thermal damage to surgical equipment (e.g. surgical drapes, plastic materials, etc.)
- If the light source fails during use, this may endanger the patient.
- ▶ Always keep an operationally ready replacement light source on hand.

Safety precautions:

- ▶ Do not illuminate the inside of the patient for longer than necessary with the light source.
- ▶ Use the AUTO light intensity control or set the light source so that a bright, well-lit image is visible with the lowest possible light intensity.
- ▶ Do not allow the distal end of the endoscope or the light guide connector to come into contact with patient tissue or with combustible or heat-sensitive materials.

Aesculap®

EinsteinVision 3.0 camera control unit PV630

- ▶ Do not touch the distal end of the endoscope.
- ▶ Remove contamination from the distal end surface or light emission surface.

WARNING

Risk of injury and/or malfunction!

- ▶ Always carry out a function check prior to using the product.

WARNING

Risk to patient through fogging of the optics and limited view!

- ▶ To prevent or reduce fogging of the optical surfaces, it can be helpful to supply the insufflation gas via an other trocar than the camera trocar.

WARNING

Risk of injury from burns and unwanted deep penetration and risk of damage to product!

- ▶ Only switch on high frequency current when the corresponding applied part (electrode) can be seen through the endoscope and there is no contact between the two.

WARNING

Electric shock due to misdirected leakage currents!

- ▶ Before thoracic or cardiothoracic surgery, deactivate implanted defibrillators (ICD).
- ▶ Before any kind of defibrillation, remove the used camera head from the patient.

WARNING

Malfunction in connection with magnetic resonance!

- ▶ Do not use the product in a magnetic resonance environment.

CAUTION

Risk of injury due to development of heat at the tip of the endoscope and sterile cover!

- ▶ Make sure that the tip temperatures on the type of the endoscope and the sterile cover prescribed by IEC/EN/DIN 60601-2-18 are not exceeded.
- ▶ Do not put down the endoscope on the patient during use.
- ▶ In the patient's abdominal cavity, always leave a sufficient distance between the optical lens and the tissue surfaces and mucous membranes of the patient.
- ▶ Use the automatic light intensity control or set the light source so that a bright, well-lit image is visible with the lowest possible light intensity.
- ▶ Switch off the light source if the illumination is no longer needed or if the endoscope is out of the patient for a longer period of time.
- ▶ Do not touch the light guide connector, distal end of the endoscope or the tip of the sterile cover during use or immediately after use.

WARNING

Risk to patient through gas embolism!

- ▶ Avoid overinsufflation before high frequency surgery (e.g., with air or inert gas).

Note

If the 3D display fails, the operation can be continued with 2D display.

Note

The current assignment of the camera head buttons is displayed on the screen as **camera head information** after opening the camera menu.

Note

A thorough understanding of the principles and methods used for example in electro-surgical procedures is necessary to avoid risks of shock or burns to patients and users as well as damage to other equipment and instruments. Electrosurgical procedures may only be conducted by trained specialist staff.

Active profile after switching on.

If no camera head is connected when the camera control unit is switched on, the standard profile "LAP General Surgery" is activated on startup.

If a 2D or 3D camera head is connected before the unit is switched on, the last profile used for the relevant camera head is activated on startup.

Active profile after changing the camera head.

If the camera head is removed from the active camera control unit, the existing profile remains active if a camera head of the same type (2D or 3D) is connected, even if there is a switch between a 3D-0° and a 3D-30° camera head.

If another camera head type is attached, the standard profile "LAP General Surgery" is activated.

User profiles dependent on camera head

Three user profiles may be individually configured for both types of camera head (2D and 3D). The relevant user profiles can only be displayed, activated or changed if the corresponding camera head type is connected.

Switching On and Off

Switching on

- ▶ Press the ON switch.

Switching off

- ▶ Press the OFF switch.

Disconnecting the device from the mains power

- ▶ Pull out the power cord on the rear of the device.

Connecting the camera head

Note

The camera cable is fragile when bent, kinked, twisted, pulled, or pressed. These may damage the optical components (e.g. light guide fibers) and thus lead to inoperability.

For further details on the camera head, see the respective instructions for use.

- ▶ Insert the connection plug fully into the connection sockets on the camera control unit, until it engages.

Profile Selection

- ▶ Open menu and navigate to menu item **Profile**.
- ▶ Select the required profile.

Performing white balance

- ▶ Perform a white balance prior to every procedure.

Note

Note that ambient light that deviates significantly from the color temperature of the light source can affect the results of the white balance.

Note

When using a sterile cover, perform the white balance after mounting the sterile cover, otherwise the result of the white balance will be falsified.

- ▶ Connect the camera head being used for the procedure.
- ▶ Connect the light guide to the light source.
- ▶ When using a 3D camera head: use a sterile cover on the camera head.
- ▶ Switch on all system components.
- ▶ Point the endoscope to a white object at a distance of approx. 5 cm.

- ▶ Press the **White Balance** button on the front of the camera control unit or the relevant button on the camera head.

If the white balance is successful, the message "white balance successful" appears on the 3D monitor.

If the white balance is unsuccessful, the message "white balance unsuccessful" appears.

If the white balance is unsuccessful:

- ▶ Prevent overexposure.
 - Enlarge the distance to the white object.
 - Adjust the brightness of the light source.
- ▶ Repeat the white balance.

Image display

The image display can be changed manually for individually created profiles. These changes are retained until the next time the type of camera head or profile is changed.

The following changes are possible:

- Brightness
- Digital zoom
- Contrast
- Edge enhancement
- Image rotation 180° (only available in connection with EV3.0 camera head 30°)
- Switch to 2D/3D (only available in connection with the EV3.0 camera head)
- ▶ Open menu and navigate to the required menu item.
- ▶ Using the arrow buttons, apply the setting to the controller or switch a function on or off.

Effects

Effects can be added to the set profile.

- ▶ Open menu and navigate to menu item **Effects**.
- ▶ Select and activate the effect.
- ▶ Deactivate the effect in the same way.

Aesculap®

EinsteinVision 3.0 camera control unit PV630

6. Reprocessing procedure

6.1 General safety notes

Note

Adhere to national statutory regulations, national and international standards and directives, and local, clinical hygiene instructions for reprocessing.

Note

Successful reprocessing of this medical device can only be ensured through a validated reprocessing procedure. The operator/reprocessing technician is responsible for this.

Note

For the latest information on reprocessing and material compatibility see also the Aesculap extranet at www.aesculap-extra.net

WARNING

Risk of infection for patients and/or user due to:

- Residues of cleaning and disinfectants on the product.
- Insufficient or incorrect cleaning and disinfection of the product and accessories

CAUTION

Risk of defects, consequential damage or a shortened product service life!

- ▶ Follow and observe the manufacturer's requirements for reprocessing

6.2 General information

DANGER

Damage to the product due to improper reprocessing!

- ▶ Use only wipe disinfection for cleaning and disinfecting the camera control unit.
- ▶ Do not, under any circumstances, clean or disinfect the product in an ultrasonic cleaning bath.
- ▶ Use cleaning and disinfecting agents which are suitable and approved for the present product.
- ▶ Observe the manufacturer's cleaning and disinfecting instructions regarding concentration, temperature and exposure time.
- ▶ Only reconnect the camera control unit to power when all cleaned parts are completely dry.
- ▶ Do not sterilize the product under any circumstances.

Note

The housing of the EV3.0 camera control unit and the polarized 3D glasses are only suitable for wiping with a moist cloth (wipe disinfection). Do not immerse into liquid.

6.3 Cleaning/disinfection

Product-specific safety instructions for the reprocessing procedure

DANGER

Risk of electric shock and fire hazard!

- ▶ Unplug the device before cleaning.
- ▶ Do not use flammable or explosive cleaning or disinfecting solutions.
- ▶ Ensure that no fluids will penetrate the product.

Note

Do not sterilize the camera control unit or reprocess it in a cleaning/disinfecting machine.

Validated cleaning and disinfection procedure

Procedure	Special features	Reference
Wipe disinfection for electrical devices without sterilization	none	Chapter Wipe disinfection for electrical devices without sterilization
<ul style="list-style-type: none"> ■ EV3.0 camera control unit PV630 		

Wipe disinfection for electrical devices without sterilization

CAUTION

Risk of electric shock and damage to the device from moisture!

- ▶ Carry out the wipe disinfection with damp cloths only. Never use dripping wet cloths.
- ▶ Remove excess cleaning solution from the device immediately.
- ▶ Only reconnect the camera control unit to power when all cleaned parts are completely dry.

Phase	Step	T [°C/°F]	t [min]	Conc. [%]	Water quality	Chemicals
I	Wipe disinfection	RT	≥2	-	-	55 % propane-2-ol + 0,5 % quarternary ammonium chlorides*

* Super Sani-Cloth Disposable Wipes were used for validation.

Phase I

- ▶ Remove any visible residues with a disposable disinfectant wipe.
- ▶ Wipe all surfaces of the optically clean product with a fresh, disposable disinfectant wipe.
- ▶ Observe the specified application time (2 min minimum).

6.4 Inspection, maintenance and checks

- ▶ Inspect the product after each cleaning and disinfecting cycle to be sure it is:
 - Clean
 - Function
- ▶ Prior to each use, inspect the product for:
 - Irregular running noise
 - Excessive heat-up
 - Too strong vibrations
 - Damage
- ▶ Set aside the product if it is damaged.

6.5 Storage

- Damage to product due to improper storage!**
- ▶ Store the product protected from dust in a dry, well ventilated, and temperature-controlled room.
 - ▶ Store the product protected from direct sunlight, high temperatures, high air humidity or radiation.
 - ▶ Do not expose the device to direct UV light, radioactivity or strong electromagnetic radiation.
 - ▶ Store the product separately or use containers in which it can be fixed in place.
 - ▶ Always transport the product with care, even if it is in an equipment trolley.

7. Repairs

Send damaged products to the manufacturer or authorized repair center. Authorized repair centers can be inquired about from the manufacturer.

Risk of infection due to soiled or contaminated products!

- ▶ Clean, disinfect and sterilize the product and any accessories thoroughly before sending. Otherwise prepare the product in as far as is possible, and label it correspondingly.
- ▶ Prior to sending, remove the sterile cover from the camera head.
- ▶ Select suitable and safe packaging (ideally the original packaging).
- ▶ Package the product in such a way that the packaging does not become contaminated.

Note

The specialized repair company may refuse to repair soiled or contaminated products for safety reasons. The manufacturer reserves the right to return contaminated products to the sender.

8. Maintenance

This product is maintenance-free. It does not contain components which have to be maintained by the user or manufacturer. However the manufacturer recommends that a qualified person or a hospital technician inspects the product regularly in a preventive safety inspection. The product should be inspected at least once every 12 months. In case of preventative safety inspections, no special conditions and measures have to be observed.

When checking, perform the following tests pursuant to section 5 of IEC/EN/DIN 62353:

- 5.2 Visual inspection
- 5.3.2 Measuring of protective earth resistance (only applicable for devices of protection class I)
- 5.3.3 Leakage currents
- 5.3.4 Measurement of insulation resistance

The test results must be documented in full pursuant to section 6 of IEC/EN/DIN 62353 in a test report which can be taken from annex G of the standard.

After performing the preventive safety inspection, a functional test is to be carried out.

For relevant services contact your national B. Braun/Aesculap agency, see Technical Service.

9. Troubleshooting list

Problem	Cause	Remedy
No image on the monitor	No power supply	Set up the equipment properly, see Preparation and setup.
	Defective fuse on the camera control unit	Change the fuse according to the instructions in "Changing a fuse".
	Controller incorrectly connected to monitor	Connect the camera control unit properly, see Preparation and setup.
	Connection cable not connected or defective	Connect the camera head to the camera control unit. Check the plug for moisture. If the cable is defective, send it for repair.
	Monitor not tuned in to the correct input signal	Tune in the monitor to the correct input signal.
Poor color rendering	Monitor not configured correctly	Check monitor settings.
Colored stripes through the image	Video cable defective	Replace the video cable or send the product for repair.
Unnatural color rendering	White balance not set correctly	Set the white balance, see Preparation and setup.
Image is distorted	Settings not completed correctly	Set the device properly, see Preparation and setup.
	Monitor not set correctly	Readjust monitor settings according to Instructions For Use of monitor.
No operation	No power supply	Check the mains connection and connect the power if necessary. Check the fuse and replace if necessary.
Not enough light	Light guide cable not connected correctly	Check that the light guide cable is connected properly; correct if necessary.
	Light guide defective	Replace the light guide.

9.1 User maintenance

Fuse replacement

DANGER

Risk of death by electric shock!

- ▶ **Unplug the device before changing the fuses!**

Specified fuses: T2,00AH/250 V~

- ▶ Use a small screwdriver to release the clip on the fuse holder.
- ▶ Pull out the fuse holder.
- ▶ Replace both fuses.
- ▶ Reinsert the fuse holder in such a way that it audibly snaps into place.

Note

If the fuses burn out frequently, the device is faulty and should be repaired, see Technical Service.

Note

Follow the country-specific regulations of the fuse.

9.2 Repairs

For repairs (except changing the fuse), contact the manufacturer or authorized repair center. Authorized repair centers can be inquired about from the manufacturer.

For a fast processing of your service requests, send in the product indicating:

- Article number (REF)
- Serial number (SN)
- Detailed description of defects

Note

If a component in the visualization system is replaced (e.g. Service, Upgrade), the system start-up process must be performed again. It is also recommended to perform a safety inspection (STK) after any maintenance or replacement of a device.

10. Technical Service

WARNING

Risk of injury and/or malfunction!

- ▶ **Do not modify the product.**

CAUTION

Malfunction of the product due to damage in transport!

- ▶ **Select suitable and safe packaging (ideally the original packaging).**
- ▶ **Keep the original packaging for possible returns in case of service.**
- ▶ **Package the product in such a way that the packaging does not become contaminated.**

Note

Service and repair work must only be carried out by authorized trained personnel. For service and repairs, please contact your national B. Braun/Aesculap agency.

Modifications carried out on medical technical equipment may result in loss of warranty rights and applicable licenses.

Warranty

The manufacturer gives a 12 month warranty for the function of the product. This warranty is restricted to claims which are immediately made in writing within the specified warranty period starting from the date of invoice, if applicable with reference to repairs, and indicating the invoice number. Statutory warranty claims are not restricted by this warranty.

This warranty is only applicable to defects that cannot be attributed to normal wear and tear, misuse, incorrect handling, third party interference, lack of or incorrect reprocessing or Force Majeure.

Warranty claims shall not be accepted if the user himself or a non-authorized repair center carries out repairs or makes changes to the product. In the event a product has to undergo maintenance, the same applies for maintenance work which is not explicitly authorized.

Liability claims arising from improper use or combinations with other devices or accessories cannot be made.

Service addresses

Aesculap Technischer Service

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78532 Tuttlingen / Germany

Phone: +49 7461 95-1601

Fax: +49 7461 16-939

E-Mail: ats@aesculap.de

Or in the US:

Aesculap Inc.

Attn. Aesculap Technical Services

615 Lambert Pointe Drive

Hazelwood

MO, 63042 USA

Aesculap Repair Hotline

Phone: +1 800 214-3392

Fax: +1 (314) 895 -4420

Other service addresses can be obtained from the address indicated above.

11. Accessories/Spare parts

Art. no.	Designation
PV647	Stand for 3D monitor
PV646	3D monitor
PV621	Polarized 3D glasses (15 units)
PV622	3D anti-mist polarization glasses
PV623	Polarized 3D glasses clip
PV624	3D eye protection polarization glasses
A075305	Fuse T2,00AH/250 V~
OP942	MIS-Bus cable, 0.75 m
PV437	DVI cable, 3.0 m
TA014803	BNC cable, red, 3.0 m
TA014804	BNC cable, green, 3.0 m
GK535	Equipotentialization cable, 4.0 m
TA008205	Equipotentialization cable, 0.8 m
TE780	Power cord for Europe, black 1.5 m
TE730	Power cord for Europe, black 5.0 m
TE734	Power cord for United Kingdom, Ireland, black 5.0 m
TE735	Power cord for USA, Canada, Japan, gray 3.5 m
FS095	Power cord for Switzerland, orange
TE676	IEC extension cord for low-heat devices 1.0 m
TE736	IEC extension cord for low-heat devices 2.5 m

12. Technical data

Classification according to Directive 93/42/EWG

Art. no.	Designation	Class
PV630	EV3.0 camera control unit	I

Mains voltage ranges (Current intake) 100 - 240 V~ (1,0 A - 0,5 A)

Power consumption 120 W

Protection class (pursuant to IEC/EN/DIN 60601-1) I

Protection rating IP21

Device protection T2,00AH/250 V~

Applied part Type CF, defibrillation-proof

Frequency 50/60 Hz

Video signal
 2 x 3D via 3G SDI (1080p)
 2 x 3D via DVI-D (1080p)
 2 x 2D via DVI-D (1080p)
 1 x 2D via HD-SDI (1080i)

Weight 10,2 kg

Dimensions (W x H x D) 330 mm x 146 mm x 353 mm

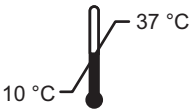
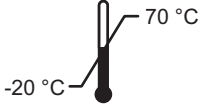
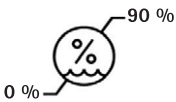
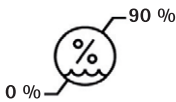
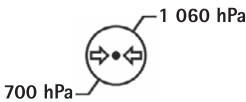
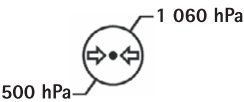
Operating mode Suitable for continuous operation

Conforming to standard IEC/EN/DIN 60601-1

EMC IEC/EN/DIN 60601-1-2

CISPR 11 Class A

12.1 Ambient conditions

	Operation	Storage and transport
Temperature		
Relative humidity		
Atmospheric pressure		

14. Distributor in the US/Contact in Canada for product information and complaints

Aesculap Inc.
 3773 Corporate Parkway
 Center Valley, PA, 18034,
 USA

13. Disposal

Note

The user institution is obliged to process the product before its disposal, see Reprocessing procedure.

 	<p>Adhere to national regulations when disposing of or recycling the product, its components and its packaging! The recycling pass can be downloaded from the Extranet as a PDF document under the respective article number. The recycling pass includes disassembling instructions for the product, as well as information for proper disposal of components harmful to the environment. Products carrying this symbol are subject to separate collection of electrical and electronic devices. Within the European Union, disposal is taken care of by the manufacturer as a free-of-charge service.</p>
--	---

- Detailed information concerning the disposal of the product is available through your national B. Braun/Aesculap agency, see Technical Service.

Note

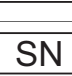
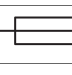
Incorrect disposal can lead to environmental damage.

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Aesculap AG | Am Aesculap-Platz | 78532 Tuttlingen | Germany