
# Aesculap Surgical Instruments


- GB** Instructions for use  
bipoJet
- USA**
- D** Gebrauchsanweisung  
bipoJet
- F** Mode d'emploi  
bipoJet
- E** Instrucciones de manejo  
bipoJet
- I** Istruzioni per l'uso  
bipoJet
- P** Instruções de utilização  
bipoJet
- NL** Gebruiksaanwijzing  
bipoJet
- RUS** Инструкция по применению  
bipoJet
- CN** 使用说明书  
bipoJet

1

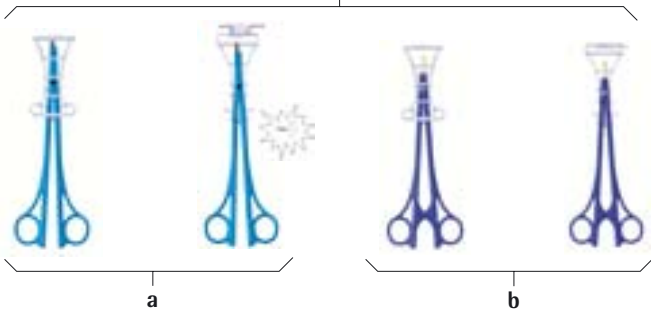2
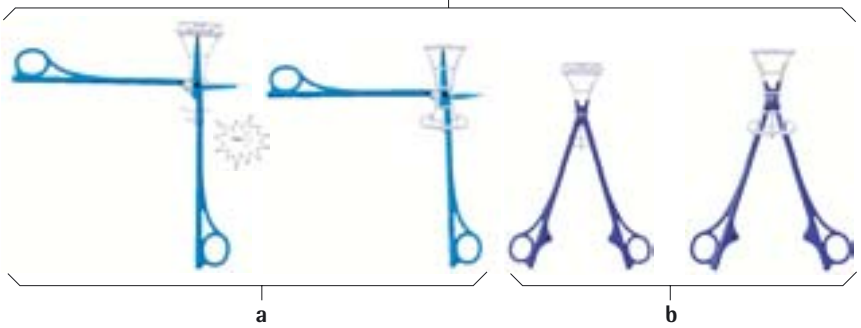
## bipoJet

## Legend

- 1 Bipolar scissors **a**/clamp **b** with ceramics protector in closed condition (for sterilization and transport)
- 2 Bipolar scissors **a**/clamp **b** with ceramics protector in open condition (for transport in the tray after operation and for mechanical cleaning)

## Symbols on product

|                                                                                   |                                     |
|-----------------------------------------------------------------------------------|-------------------------------------|
|  | Attention, see instructions for use |
|-----------------------------------------------------------------------------------|-------------------------------------|

## Intended use

The bipolar instruments are used for clamping, cutting, preparing and coagulating of tissues in surgical operations.

## Safe handling and preparation

## CAUTION

**Federal law restricts this device to sale by or on order of a physician!**

- Read, follow and keep the instructions for use.
- Use the product only in accordance with its intended use, see Intended use.
- Clean the new product either manually or mechanically prior to the initial sterilization.
- Store any new or unused products in a dry, clean and safe place.
- Prior to each use, inspect the product for: loose, bent, broken, cracked, worn, or fractured components.
- Do not use the product if it is damaged or defective. Set aside the product if it is damaged.
- Replace any damaged components immediately with original spare parts.

**WARNING**

**Danger of injury from burns and/or explosion from flammable gases!**  
Application of the HF device according to its intended use can involve spark formation.

- Follow the safety notes in the instructions for use of the HF device.

**WARNING**

**Risk of burns due to tissue outside the coagulation site touching non-insulated areas of the bipolar instruments!**

- Be careful that only insulated areas of the bipolar instruments come into contact with tissue.

#### **WARNING**

**Risk of injury to the patient/ medical personnel due to incorrect handling!**

- **Operate the instrument with authentic Aesculap HF cables only.**
- **Do not connect the instrument to the monopolar output.**
- **Do not use instruments with damaged ceramic parts.**

#### **CAUTION**

**Damage to the working ends and ceramics caused by incorrect handling!**

- **Apply caution when operating the instrument.**
  - **Do not apply excessive force.**
  - **Do not try cutting or removing clamps or clips.**
  - **Protect the instrument against knocks and impacts.**
  - **Transport and process the instrument only in the ceramics protector.**
  - **Prior to each application, inspect the ceramic insulation of the instrument for missing parts and/or cracks.**
- Adjust the HF power output to the intended operation. Take into account clinical experience or reference values.
- Select the lowest possible HF power output.
- Keep the product's contact surfaces clean during the operation. Wipe off encrusted tissue residues or body fluids, using a moistened swab.

The plug end of the product is fitted with the following connector: Pin 1.6 mm

Refer to our brochures to find the compatible cable.

The insulation of the product is rated for a maximum recurrent peak voltage of  $600 V_p \approx 1,200 V_{pp}$  in coagulation mode.

The insulation of the accessories (e.g. HF cable, instruments) must be of adequate rating for the maximum peak output voltage, see IEC 60601-2-18.

To avoid HF burns:

- Set the HF device for a maximum peak output voltage smaller than or equal to the peak voltage rating specified for the product.
- The product's working tip must be in the user's field of vision whenever HF power is activated.
- Prior to activating the HF device, make certain that the working end of the product is not touching any electrically conductive accessories.
- Prior to each use, visually check the product for: damage to, or surface changes on, the insulation.
- Never put down the product on or next to the patient.
- Follow the instructions for use of the HF device.

## bipoJet

## Function checks

- Connect the HF cable to the instrument.
- Connect the plug of the HF cable to the bipolar output of the HF device.
- Push the tips of the slightly opened instrument onto a cloth soaked in a saline solution.
- Activate the HF device, making sure that the instrument is only touched where it is insulated. The onset of steam or sparks means that the instrument is ready for operation.
- If no steam or sparks become visible, check the connections between the HF cable and the HF device.
- Slowly increase the current supplied.

## Safe operation

**Risk of injury due to unintentional hemorrhage caused by a hasty cut or the HF current too high or too low!**

- Cut slowly.
- Adjust the HF current to a level appropriate for the application.

## Validated processing procedure

*Note*

*Observe all relevant national regulations and standards with regard to processing.*

*Note*

*For patients with Creutzfeldt-Jakob disease (CJD), suspected CJD or possible variants of CJD, observe the relevant national regulations concerning the reprocessing of the products.*

*Note*

*Mechanical processing should be preferred over manual cleaning because of the better and more reliable cleaning results of mechanical processing.*

*Note*

*Successful processing of this medical product can only be ensured through a validated processing procedure. The user/processor is responsible for the validation.*

*Due to process tolerances, manufacturer's specifications can only serve as guide values for the processing procedures applied by the individual user.*

*Note*

*Up-to-date information on processing can be found on the Aesculap Extranet at [www.aesculap-extra.net](http://www.aesculap-extra.net)*

### General notes

Encrusted or fixated residues from surgery can make the cleaning process more difficult or ineffective, and can cause corrosion of "stainless" steels. To avoid this, the time interval between application and processing should not exceed 6 h, and neither fixating pre-cleaning temperatures > 45 °C nor any fixating disinfecting agents (active ingredient: aldehyde, alcohol) should be used.

Excessive doses of neutralizers or basic detergents can cause chemical deterioration and/or fading of the laser engraving on stainless steel.

Residues containing chlorine or chlorides, as contained e.g. in surgical residues, tinctures, medicines, saline solutions, service water used for cleaning, and in cleaning/disinfecting agents, cause corrosion damage (pitting, stress corrosion) and thereby destruction of stainless steel products. To remove such residues, the products must be rinsed sufficiently with fully desalinated water and dried thoroughly.

Only use process chemicals recommended by the chemicals' manufacturers as effective cleaning/disinfecting agents compatible with the materials used in the respective product. Strictly observe all application instructions, e.g. regarding temperatures, concentrations, exposure times, etc. Failure to do so can result in the following problems:

- Optical changes to the material, e.g. fading or discoloration of titanium or aluminum. For aluminum, the application/process solution only needs to be of pH > 8 to cause visible surface changes or
- Material damage, e.g. corrosion, cracks, fracturing, premature aging or swelling.
- Do not use process chemicals that cause stress cracking of plastics such as PPSU, or attack and cause brittleness e.g. of silicone softeners.

- Clean the product immediately after use.  
Encrusted debris on HF instruments is softened most effectively and conservatively by 5 minutes immersion treatment in a 3 % H<sub>2</sub>O<sub>2</sub> solution. Encrustations can be removed by hand, with a soft brush. After that treatment, proceed through the normal steps of processing.
- Further detailed advice on hygienically safe and material-/value-preserving reprocessing can be found at [www.a-k-i.org](http://www.a-k-i.org), Publications, Red Brochure – Proper maintenance of instruments.
- Use suitable cleaning/disinfecting agents if the product is put away in wet condition. Prior to mechanical cleaning and disinfecting, rinse the product thoroughly with running water.

### Preparations at the place of use

- Open jointed products.
- Remove visible residues as completely as possible, using a lint-free, single-use cleaning tissue.

### Preparation prior to cleaning

- Have the product dry in a disposal container and ready for immediate cleaning and disinfecting within 30 min after use.

## bipoJet

## Cleaning/Disinfecting

## CAUTION

Damage to the product due to inappropriate cleaning/disinfecting agents and/or excessive temperatures!

- Use cleaning and disinfecting agents approved for the product (e.g. aluminum, Peek, PTFE, surgical steel) according to manufacturer's instructions.
- Observe specifications regarding concentration, temperature and exposure time.

## Manual cleaning/disinfecting

- Clean hinged or jointed products in open and closed positions.
- Check visible surface for residues after manual cleaning/disinfecting.
- Repeat the cleaning process if necessary.

## Manual cleaning with immersion disinfecting and cleaning with brush

| Stage | Step                     | T<br>[°C/°F] | t<br>[min] | Conc.<br>[%] | Water<br>quality | Chemical                                                 |
|-------|--------------------------|--------------|------------|--------------|------------------|----------------------------------------------------------|
| I     | Disinfecting<br>cleaning | RT<br>(cold) | 15         | 2            | D-W              | BBraun Stabimed; aldehyde phenol<br>and QAV free; pH = 9 |
| II    | Intermediate rinse       | RT<br>(cold) | 1          | -            | D-W              | -                                                        |
| III   | Final rinse              | RT<br>(cold) | 0.5        | -            | FD-W             | -                                                        |
| IV    | Drying                   | RT           | -          | -            | -                | -                                                        |

D-W: Drinking water

FD-W: Fully desalinated water (demineralized)

RT: Room Temperature

### Stage I

- Fully immerse the product in the disinfecting solution. Make certain that all accessible surfaces are moistened.
- Clean the product under running tap water, if necessary with a suitable cleaning brush (e.g. toothbrush), until all visible residues have been removed from the surface.
- Brush through all surfaces that are not accessible to visual inspection, e.g. in products with hidden crevices, lumens or complex geometries, for at least 5 min or until no more residues can be removed. Mobilize non-rigid components, e.g. set screws, links, etc. during cleaning.
- After that, thoroughly rinse (at least 5 times) with the cleaning solution these difficult to access areas of the product, using a single-use syringe (20 ml).
- To avoid the risk of corrosion, do not use a metal brush or other abrasives, which would damage the surfaces, for cleaning the product.

### Stage II

- Rinse the product completely (all accessible surfaces) under running water.

### Stage III

- Carry out a full rinse of the product (all accessible surfaces) under running water.
- Allow water to drip off for a sufficient length of time.

### Stage IV

- Dry the product completely with compressed air.

## bipoJet

**Mechanical cleaning/disinfecting**

- Put the ceramics protector on the product.
- Place the product on a tray that is suitable for cleaning (avoid rinsing blind spots).
- Store products with hinges or joints on the tray in such a way that the joints are kept open.

**Mechanical cleaning/disinfecting with manual pre-cleaning****Manual pre-cleaning with brush**

| Stage     | Step                             | T<br>[°C/°F] | t<br>[min] | Conc.<br>[%] | Water<br>quality | Chemical                                                 |
|-----------|----------------------------------|--------------|------------|--------------|------------------|----------------------------------------------------------|
| <b>I</b>  | <b>Disinfecting<br/>cleaning</b> | RT<br>(cold) | 15         | 2            | D-W              | BBraun Stabimed; aldehyde phenol<br>and QAV free; pH = 9 |
| <b>II</b> | <b>Intermediate rinse</b>        | RT<br>(cold) | 1          | -            | D-W              | -                                                        |

D-W: Drinking water

RT: Room Temperature

**Stage I**

- Fully immerse the product in the disinfecting solution. Make certain that all accessible surfaces are moistened.
- Clean the product with a suitable cleaning brush (e.g. toothbrush) until all visible residues have been removed from the surface.
- Brush through all surfaces that are not accessible to visual inspection, e.g. in products with hidden crevices, lumens or complex geometries, for at least 5 min or until no more residues can be removed. Mobilize non-rigid components, e.g. set screws, links, etc. during cleaning.
- After that, thoroughly rinse (at least 5 times) with the cleaning solution these difficult to access areas of the product, using a single-use syringe (20 ml).
- To avoid the risk of corrosion, do not use a metal brush or other abrasives that would damage the surfaces for cleaning the product.

**Stage II**

- Rinse the product completely (all accessible surfaces) under running water.

### Mechanical alkaline cleaning and thermal disinfecting

Machine type: Single-chamber cleaning/disinfecting machine without ultrasound

- Place the product on a tray that is suitable for cleaning (avoid rinsing blind spots).
- Store products with hinges or joints on the tray in such a way that the joints are kept open.

| Stage | Step                    | T<br>[°C/°F] | t<br>[min] | Water<br>quality | Chemical                                                                                       |
|-------|-------------------------|--------------|------------|------------------|------------------------------------------------------------------------------------------------|
| I     | Prerinse                | < 25/77      | 3          | D-W              | -                                                                                              |
| II    | Cleaning                | 55/131       | 10         | FD-W             | - concentrate, alkaline:<br>pH = 10.9<br>< 5 % anionic tensides<br>- 1% solution:<br>pH = 10.5 |
| III   | Neutralization          | 20/68        | 2          | FD-W             | - concentrate, acid:<br>pH = 2.6<br>Basis: Citric acid<br>- 0.1% solution:<br>pH = 3.0         |
| IV    | Intermediate rinse      | 70/158       | 1          | FD-W             | -                                                                                              |
| V     | Thermal<br>disinfecting | 94/201       | 10         | FD-W             | -                                                                                              |
| VI    | Drying                  | 90/194       | 40         | -                | -                                                                                              |

D-W: Drinking water

FD-W: Fully desalinated water (demineralized)

## bipoJet

**Control, care and inspection**

- Allow the product to cool down to room temperature.
- Lightly lubricate moving parts (e.g. joints and links) with a maintenance oil suitable for the respective sterilization process (e.g. Aesculap STERILIT® spray JG600 or maintenance oil JG598).
- Inspect the product after each cleaning and disinfecting cycle to be sure it is clean, functioning properly, not damaged, has intact insulation and does not have any loose, bent, broken, cracked, worn, or fractured components.
- Check for compatibility with associated products.
- Set aside the product if it is damaged.

**Packaging**

- Appropriately protect products with fine working tips.
- Sort the product into its appropriate storage device or put it on a suitable tray.  
Make certain that any cutting edges are protected.
- Package trays appropriately for the sterilization process (e.g. in Aesculap sterile containers).
- Pack the product in such a way that the packaging will prevent recontamination of the product in the period between reprocessing and reuse.

**Sterilization method and parameters***Note*

*To prevent breakage due to stress crack corrosion, sterilize the instruments with the locking mechanism open, or locked on the first ratchet tooth.*

- Make certain that all external and internal surfaces will be exposed to the sterilizing agent (e.g. by opening all valves and faucets).
- Validated sterilization process
  - Steam sterilization through fractionated vacuum process
  - Steam sterilization acc. to EN 285/ANSI/AAMI/ISO 11134–1993, ANSI/AAMI ST 46–1993, validated acc. to EN ISO 17665 or EN 554/ISO 13683
  - Sterilization through fractionated vacuum process at 134 °C/holding time 5 min
- When sterilizing several products at the same time in one steam sterilizer: Make certain that the maximum allowable load capacity of the steam sterilizer, as specified by the manufacturer, is not exceeded.

### Sterilization for the US market

- Sterilization of the device may be accomplished by steam.
- Aesculap does not recommend the device be sterilized by "Flash" or chemical sterilization.
- Surgical instruments may also be placed within an Aesculap rigid sterilization container (sterile container) for processing under generally accepted hospital in-use conditions.

The recommended sterilization parameters are as follows:

| Sterilization method | Temp.          | Minimum exposure time |                               |
|----------------------|----------------|-----------------------|-------------------------------|
|                      |                | Wrapped               | In a sterile container system |
| Pre-vacuum           | 270—<br>275 °F | 4 min                 | 4 min                         |

### WARNING for the US market

**If this device is/was used in a patient with, or suspected of having Creutzfeldt-Jakob Disease (CJD), the device cannot be reused and must be destroyed due to the inability to reprocess or sterilize to eliminate the risk of cross-contamination!**

### Storage

- Store processed products in germ-proof packaging in a dust-protected, dry, dark and temperature-controlled room.

### Technical Service

For service, maintenance and repairs, please contact your national B. Braun/Aesculap agency.

Modifications carried out on medical technical equipment may result in loss of guarantee/warranty rights and forfeiture of applicable licenses.

### Service addresses

Aesculap Technischer Service  
Am Aesculap-Platz  
78532 Tuttlingen / Germany  
Phone: +49 7461 95-1602  
Fax: +49 7461 16-5621  
E-mail: [ats@aesculap.de](mailto:ats@aesculap.de)

Or in the US:  
Aesculap Inc.  
Attn. Aesculap Technical Services  
615 Lambert Pointe Drive  
Hazelwood, MO 63042  
Aesculap Repair Hotline  
Phone: +1 800 214-3392  
Fax: +1 314 895-4420  
Other service addresses can be obtained from the address indicated above.

### Distributor in the US/Contact in Canada for product information and complaints

Aesculap Inc.  
3773 Corporate Parkway  
Center Valley, PA 18034  
USA

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**Aesculap AG**