Aesculap Endoscopic Technology

(GB) Instructions for use/Technical description

(USA) EinsteinVision 3.0 Camera Heads PV631/PV632 and sterile covers for PV641SU/PV642SU

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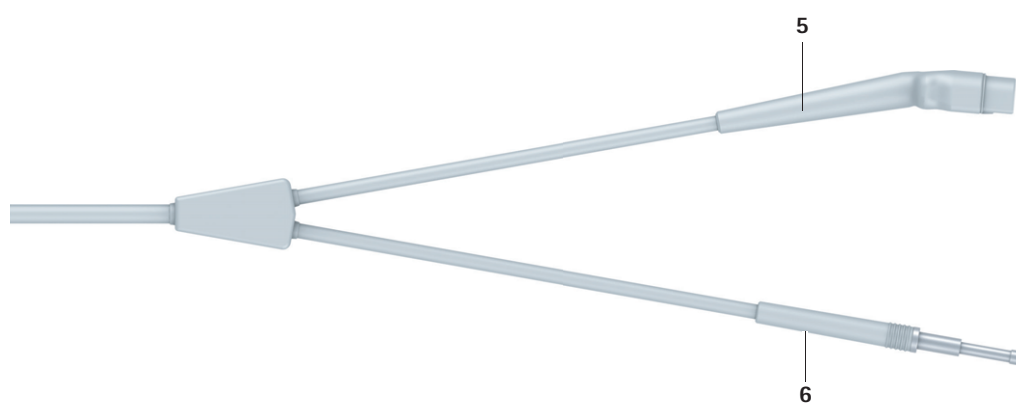
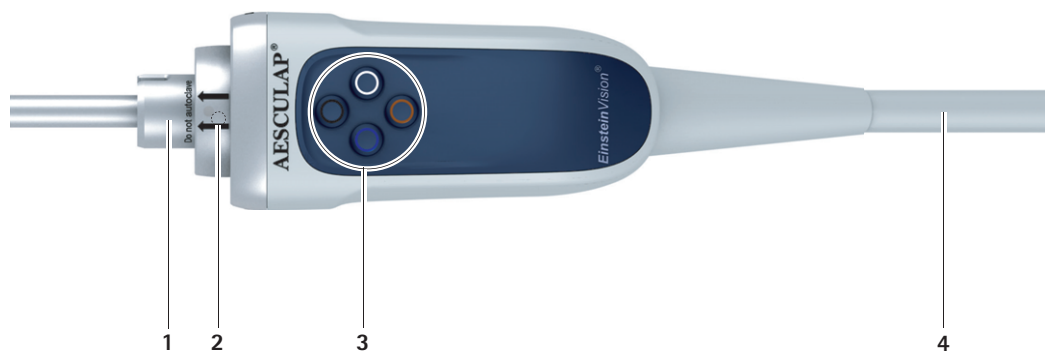
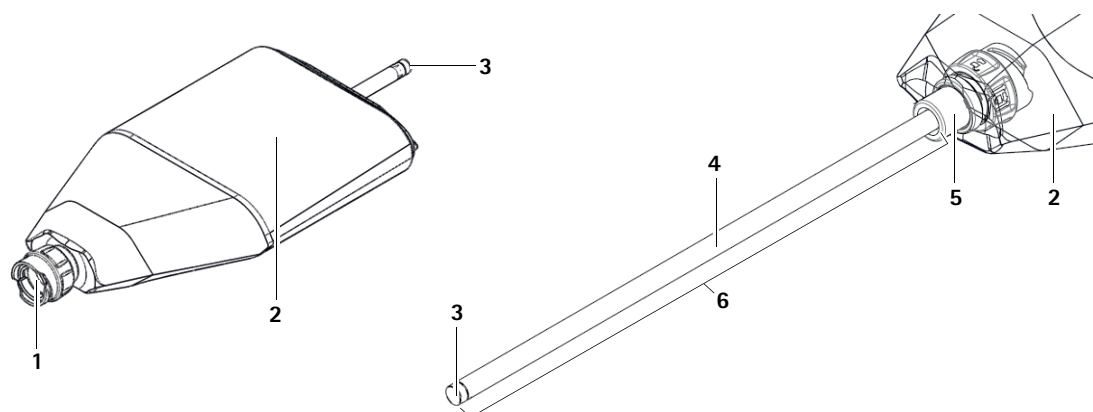
B

Aesculap®

EinsteinVision 3.0 Camera Heads PV631/PV632 and sterile covers for PV641SU/PV642SU

Legend

A

- 1 Endoscope shaft
- 2 Coupling for sterile cover
- 3 Camera head buttons
- 4 Connection cable to camera control unit/light source
- 5 Connector for camera control unit
- 6 Connector for light source

B

- 1 Coupling for connection to camera head
- 2 Sterile sleeve
- 3 Tip segment
- 4 Shaft (315 mm working length, 10 mm outer diameter)
- 5 Adapter for fastening onto endoscope holder
- 6 Insertion area of sterile cover

Symbols on product and packages

	Signal word: WARNING Indicates a hazard. If not avoided, the hazard can result in death or serious injury.
	Signal word: CAUTION Indicates a potential hazard. If not avoided, this hazard may result in injury and/or damage to the product.
	Follow the instructions for use
	Caution (IEC/EN/DIN 60601-1 3rd edition)/ Attention, read accompanying paperwork (IEC/EN/DIN 60601-1 2nd edition)
	Defibrillation protected application part type CF pursuant to IEC/EN/DIN 60601-1
	Permissible storage temperature
	Permissible relative air humidity during storage
	Permissible atmospheric pressure during storage

	Not permitted for use in a magnetic resonance environment
	Caution, fragile
	Sterilization using irradiation
	Not for reuse as per the manufacturer-defined intended use
	Do not re-sterilize
	Do not use if the packaging is damaged
	Use by
	Batch number
	Article number
	Serial number
	Manufacturer
	Date of manufacture
	Federal Law (US) restricts this device to sale to or on behalf of a physician
	Marking of electric and electronic devices according to directive 2002/96/EC (WEEE). see Disposal
	Non sterile
	Do not autoclave

Contents

1.	Applicable to	3
2.	Safe handling	4
3.	Product description	6
3.1	Scope of supply	6
3.2	Components required for operation	6
3.3	Intended use	6
3.4	Operating principle	7
4.	Preparation	8
4.1	Combination with medical electrical equipment	9
4.2	Inspections	9
5.	Working with the camera heads	10
5.1	General notes on use	10
5.2	System set-up	10
5.3	Function checks	12
5.4	Safe operation	13
6.	Reprocessing procedure	16
6.1	General safety instructions	16
6.2	General information	16
6.3	Dismantling prior to carrying out the reprocessing procedure	17
6.4	Preparation before cleaning	17
6.5	Cleaning/disinfection	17
6.6	Wipe disinfection for electrical devices without sterilization	18
6.7	Manual cleaning/disinfecting	19
6.8	Mechanical cleaning/disinfecting	21
6.9	Material-compatible cleaning and disinfecting agents	23
6.10	Inspection, maintenance and checks	23
6.11	Packaging	23
6.12	Storage	23
7.	Service	24
8.	Maintenance	24
9.	Troubleshooting list	25
9.1	Repairs	25
10.	Technical Service	26
11.	Accessories/Spare parts	26
12.	Technical data	26
12.1	Camera head	26
12.2	Sterile cover	27
12.3	Ambient conditions	27
13.	Disposal	27
14.	Distributor in the US/Contact in Canada for product information and complaints	27

1. Applicable to

These instructions for use apply to the following devices:

Designation	Art. no.
EV3.0 camera head 10 mm 0°	PV631
EV3.0 camera head 10 mm 30°	PV632
Sterile cover for Camera Head 0°	PV641SU
Sterile cover for camera head 30°	PV642SU

These instructions for use are an integral component of the device and contain all the information required by users and operators for safe and proper use.

Target group

These instructions for use are intended for physicians, medical assistants, medical technicians and employees of sterile services who are entrusted with the installation, operation, maintenance and reprocessing of the device.

Using and storing this document

These instructions for use must be stored in a defined location so that it may be accessed at all times by the target group.

In the event of the sale of this device or its relocation, this document must be handed over to the new owner.

Supplementary documents

For the use and assignment of the camera head buttons and the related options and possible functions, observe the instructions for use of the camera control unit concerned.

The instructions for use of the sterile cover and instructions for use of all other devices used must be observed to ensure the safe use of the device.

► For item-specific instructions for use and information on material compatibility, see also the Aesculap Extranet at <https://extranet.bbraun.com>

EinsteinVision 3.0 Camera Heads PV631/PV632 and sterile covers for PV641SU/PV642SU

2. Safe handling

CAUTION

Federal law restricts this device to sale by, or on order of a physician!

WARNING

Risk to patient from a defective sterile barrier!

- ▶ Use sterile cover for EV3.0 Camera Head (PV641SU and PV642SU) only with EV3.0 Camera Head (PV631 and PV632).
- ▶ If the sterile cover for EV3.0 Camera Head is inserted in the patient without a trocar, this must be done completely without friction. No lateral pressure must be applied to the shaft in such a case.
- ▶ The sterile cover is delivered sterile. Cleaning, disinfection and/or sterilization of the sterile cover are prohibited.
- ▶ To prevent overheating, the sterile cover for EV3.0 Camera Head (PV641SU and PV642SU) must only be used with a EV3.0 Camera Head (PV631 and PV632), if the LED light source OP950 is used in addition to the EV3.0 Camera Control Unit.

WARNING

Risk to user and patient from non-observance of instructions, warnings and precautions!

- ▶ Use the product only according to these instructions for use.
- ▶ Check the product is in good working order before every use.
- ▶ Carry out function checks before using the equipment.
- ▶ Do not use the product if irregularities are observed during inspections and checks.

WARNING

Risk to patient due to incorrect application!

- ▶ Use the product only after instruction by the manufacturer or by an authorized person.
- ▶ The instructions for use enclosed with the individual EinsteinVision components (e.g., high frequency surgery) and all products used must be observed.
- ▶ Endoscopic procedures may only be carried out by specialists who have the necessary medical training, knowledge and experience.

WARNING

Risk to patient and user due to premature wear!

- ▶ Handle and maintain the product properly.
- ▶ Use the product only in accordance with the intended purpose.

WARNING

Risk of injury due to electric shock!

- ▶ When installing a medical electrical system, there is a risk of fire, short-circuit or electric shock. Installation must be performed by qualified personnel only.
- ▶ For the combination of electrical devices, follow annex 1 of IEC/EN/DIN 60601-1. Non-medical devices that comply with the applicable IEC safety standards can be connected only via a medical isolating transformer. Do not connect any additional non-medical devices to a medical electrical system.
- ▶ Signal lines from devices with a functional connection, which are connected to different branches of the mains power supply, must be galvanically isolated at both ends.
- ▶ Connect the devices only to a power supply with a protective earth conductor.
- ▶ After installation of a medical electrical system, perform a test pursuant to IEC/EN/DIN 62353.

WARNING

Risk to the patient if the product is not used gently / a damaged product is used!

- ▶ Handle the product with appropriate care.
- ▶ Do not use the product if it has been subjected to strong mechanical stress or if it has been dropped, and send the product to the manufacturer or an authorized repair center for inspection.

WARNING

Risk of damage to eyes from UV radiation!

- ▶ Never use polarized 3D glasses as sunglasses.

WARNING

Risk to patient if the device power fails!

- ▶ Use the product only with an uninterruptible power supply. A medical emergency power supply (UPS) is recommended here.

WARNING

Restricted function if third party devices are used!

- ▶ Use the product only with the recommended components and accessories.
- ▶ Full functionality is only guaranteed if the recommended components and accessories are used.

CAUTION

Device failure due to incorrect storage and usage conditions!

- ▶ Store and operate the product only within the specified environmental conditions.

WARNING

Risk to patient from a defective sterile barrier!

- ▶ To prevent overheating, the sterile cover for EV3.0 Camera Head (PV641SU and PV642SU) must only be used with a EV3.0 Camera Head (PV631 and PV632), if the LED light source OP950 is used in addition to the EV3.0 Camera Control Unit.

WARNING

Risk to patient due to failure of a device!

- ▶ Keep an operationally ready replacement device and sterile cover on hand. If necessary change to conventional surgical methods.

Note

Only use the sterile cover with a trocar system designed for 10 mm endoscopes.

Note

Additional equipment and/or peripheral devices which are connected to the interfaces of the product must show their relevant specifications (such as IEC/EN/DIN 60601-1).

Note

The product may not be modified in any way.

Note

To guarantee optimal function of the EinsteinVision components, use in controlled environmental conditions is recommended (e.g., air-conditioned operating room).

Note

Always handle the product very carefully as it contains sensitive optic, mechanic and electronic components. Do not strike or drop the EV3.0 camera head.

Note

Do not store the polarized 3D glasses in an environment with high temperatures (e.g., near a radiator).

Note

All devices operated in the vicinity must meet the relevant EMC requirements.

Note

The EinsteinVision components are used to visualize the inside of the body during minimally invasive procedures. Do not use the components for diagnostic purposes. This applies in particular with the use of image optimization algorithms.

Note

Images and videos recorded are for documentation purposes only. Do not use for diagnosis or findings, as the image quality is reduced in the event of data compression.

Note

Before putting into operation, check the compatibility of all components using the accessories list.

Note

All accessories and spare parts must only be procured from the manufacturer.

Note

A thorough understanding of the principles and methods used for example in electro-surgical procedures is necessary to avoid risks of shock or burns to patients and users as well as damage to other equipment and instruments. Electro-surgical procedures must only be performed by trained personnel.

- ▶ Clean the new product (manually or mechanically) after removal of the transport packaging.
- ▶ Prior to use, check that the product is in good working order.
- ▶ Observe "Notes on Electromagnetic Compatibility (EMC)", see TA014630.
- ▶ To prevent damage caused by improper setup or operation, and in order not to compromise warranty and manufacturer liability:
 - Use the product only according to these instructions for use.
 - Follow the safety and maintenance instructions.
 - Only combine Aesculap products with each other.
 - Follow the application advisories acc. to standard, see Extracts from relevant standards.

EinsteinVision 3.0 Camera Heads PV631/PV632 and sterile covers for PV641SU/PV642SU

- ▶ Ensure that the product and its accessories are operated and used only by persons with the requisite training, knowledge, or experience.
- ▶ Keep the instructions for use accessible for the user.
- ▶ Always adhere to applicable standards.

The products PV641SU and PV642SU are gamma-sterilized and are supplied in sterile packaging.

- ▶ The product must not be reused.
- ▶ Do not use products from open or damaged sterile packaging.
- ▶ Before using the sterile cover, inspect it for damage such as rough surfaces, sharp edges or protrusions to avoid injuries to the patient.
- ▶ Cleaning, disinfection and/or (re-)sterilization of the sterile cover are prohibited.
- ▶ Prior to each use, inspect the product for: Loose, bent, broken, cracked, worn, or fractured components.
- ▶ Do not use the product if it is damaged or defective. Set aside the product if it is damaged.
- ▶ Do not use the product after its use-by date.

These components are optimally aligned with one another and offer the best possible quality and unrestricted functionality.

CAUTION

If the recommended components are not used, the following errors/complications may occur, among others:

- Incorrect color display
- No/incorrect display of the camera menu
- No/incorrect 3D display
- Restricted operating function from the camera head
- Increased fogging on the outlet of the optical lens
- Damage to the sterile cover
- No switching between 2D and 3D display

Note

If the product is used with third party devices, unrestricted responsibility is passed to the system configurator or medical institution.

3. Product description

3.1 Scope of supply

Designation	Art. no.
EV3.0 Camera Head 10 mm 0°/30°	PV631/PV632
Instructions for use	TA014618

3.2 Components required for operation

Camera Heads

The EV3.0 Camera Heads

- EV3.0 Camera Head 10 mm 0°/30° (PV631/PV632)

are used for 3D visualization during endoscopic surgery.

The integrated light guide serves to transmit the light between light source and the EV3.0 Camera Head. The EV3.0 Camera Head 10 mm 0°/30° (PV631/PV632) has to be used in combination with the Sterile Cover for EV3.0 Camera Head 0°/30°.

The EV3.0 Camera Heads are intended for use with:

- PV630 EV3.0 Camera Control Unit
- PV641SU/PV642SU Sterile Cover for EV3.0 Camera Head 0°/30° (for PV631/PV632)
- OP950 LED Light Source

Sterile Cover

The disposable sterile cover for EV3.0 Camera Head 0°/30° is part of the sterile concept for the EV3.0 Camera Head 10 mm 0°/30°. The sterile cover will be shifted over the scope and is the sterile barrier between the unsterile scope and patient. It guarantees the safe use of the endoscope under sterile conditions. The sterile cover is available with a 0° (PV641SU) or 30° (PV642SU) tip.

3.3 Intended use

The EV3.0 Camera Head 10 mm 0°/30° and the other EinsteinVision components are used for 3D visualization during:

- Laparoscopic standard endoscopy*
- Minimally invasive thoracic and cardiothoracic surgical procedures*

*Currently existing clinical data with regard to the effectiveness and safety of an application of EinsteinVision components to children are limited. Any application to this patient population should therefore be carried out at the discretion of the surgeon and after thorough consideration of the individual risk-benefit situation. The decision should in particular be based on the anatomical facts/proportions as well as the size and place of the surgical access needed.

In order to assure a safe application, standards/guidelines released by official organizations such as scientific societies of endoscopic surgery have to be considered.

Observe the additional intended uses of the relevant components.

Indications

The EV3.0 Camera Head 10 mm 0°/30° in combination with the sterile cover for EV3.0 Camera Head 0°/30° and the other EinsteinVision components are indicated for visualization during:

- urological, gynecological and general laparoscopy
- video-assisted minimally invasive mitral valve surgery
- video-assisted minimally invasive tricuspid valve surgery
- video-assisted minimally invasive ablation (radiofrequency cryoablation) of atrial fibrillation
- video-assisted minimally invasive ASD (Atrial Septal Defect) and VSD (Ventricular Septal Defect) closure
- video-assisted thoracic surgery lobectomy

Contraindications

The use of the EinsteinVision components and its accessories is contraindicated if endoscopic procedures are contraindicated for any reason. As is appropriate with any surgical procedure, consideration must be given to patient size and workspace volume when using the EinsteinVision components.

Depending on the patient's disease, there may exist contraindications that rely on the general condition of the patient or the specific disease pattern. The decision to perform an endoscopic procedure rests with the responsible surgeon and should be made on the basis of an individual risk-benefit analysis.

3.4 Operating principle

Note

The EV3.0 Camera Head can only be used for 3D visualization if the LED light source OP950 is used in addition to the EV3.0 Camera Control Unit. For this, the Camera Control Unit and light source must be correctly connected via a MIS-Bus cable.

CAUTION

Risk of intolerances such as vertigo, headache and nausea due to three-dimensional visualization!

- Before first use, check tolerance. If intolerances occur during use, switch over to two-dimensional imaging.

The camera heads within the scope of these instructions for use can be used with the EinsteinVision 3.0 Camera Control Unit and may only be connected to this unit.

Together, the camera control unit and camera head comprise the camera. In combination with a suitable monitor, the camera can provide either two-dimensional or three-dimensional images. For three-dimensional visualization, the user requires polarized glasses.

The camera head features four configurable camera head buttons that can be used to control image display, record images and video, or to navigate and change the settings in the configuration menu of the camera control unit.

The EV3.0 Camera Heads (0°/30°) must be used with the applicable sterile cover (0°/30°).

The sterile covers within the scope of application of these instructions for use can be used with the EV3.0 Camera Head 10 mm 0° (PV631) and the EV3.0 Camera Head 10 mm 30° (PV632).

The sterile covers are used during surgical procedures to create a sterile barrier between the camera head being used and the patient. The sterile cover is securely locked on the camera head. After being correctly mounted, the sterile cover fully covers the camera head and camera cable.

Camera head buttons

The camera head buttons have different functions depending on the operating mode.

The assignment of the buttons is always synchronous with the assignment of the buttons on the front of the camera control unit.

Different functions can be assigned for a short press of the button (<2 sec.) and a long press of the button (≥2 sec.).

The procedure for changing the assigned functions is described in the instructions for use of the EV3.0 Camera Control Unit PV630.

Button assignment in live mode

The button assignment in live mode depends on whether a standard profile or user profile is set. If a user profile is set, the buttons can be assigned individually, however the assignment for entering the camera control unit menu cannot be changed. In standard profile, the buttons are assigned as follows.

	Short press of button	Long press of button
	Image rotation 180°	Open menu
	Zoom	White balance
	Light source on	Light source on/off
	Single image	Video recording start/stop

Aesculap®

EinsteinVision 3.0 Camera Heads PV631/PV632 and sterile covers for PV641SU/PV642SU

Button assignment when on-screen display is active

When on-screen display is active, the buttons are used to navigate the menu until the menu is left.

	Short press of button	Long press of button
	up/control +	Leave the menu
	right/next level across/save control setting	
	down/control -	
	left/back a level/save control setting	

Button assignment when on-screen keyboard is active

When the on-screen keyboard is active, the buttons are used to navigate the keyboard until the keyboard is left.

	Short press of button	Long press of button
	up	-
	right	Confirm/enter
	down	-
	left	-

Usable trocars

The endoscope with corresponding sterile cover is intended for use with a trocar system designed for 10 mm endoscopes.

Due to possible deviations of the inner diameter of the trocars, not all 10 mm trocars may be compatible. This must be checked before the surgical procedure. The sterile cover must be able to be inserted into the trocar without using force.

4. Preparation

If the following instructions are not observed, the manufacturer assumes no responsibility for possible consequences.

Before installation and use, ensure that the:

- electrical installations comply with the relevant technical regulations,
- relevant regulations concerning fire and explosion protection are observed.

Note

The device is authorized for operation in hospitals.

Note

For the safety of the user, it is essential that the main power cord and, especially, the protective earth connection are intact. In many cases defective or missing protective earth connections are not registered immediately.

WARNING

Risk to patient from misdirected leakage current due to incorrect or defective earthing!

- Do not touch the product and patient at the same time.

WARNING

Risk to patient due to failure of a device!

- Keep an operationally ready replacement device and sterile cover on hand and if necessary change to conventional surgical methods.

WARNING

Danger of explosion from a device that is improperly set up!

- Make certain that the power plug is connected to the power supply outside of any areas where there is a risk of explosion.
- Do not use the product in potentially explosive areas or close to highly flammable or explosive gases (e.g., oxygen, anesthetic gases).

CAUTION

Risk to patient and user due to injury and infection!

- Do not use products with damaged fiber optics, damaged glass surfaces or stubborn deposits which cannot be removed by cleaning.

WARNING

Risk to patient due to use of damaged products!

- Do not use products if they have sharp edges or other hazardous surface damage.

WARNING

Risk of injury to the eyes!

- When checking the light guide, do not plug the light guide connector into the light source.

WARNING

Risk to patient and user due to condensation and short-circuit!

- ▶ Before putting into operation make sure that all components used have enough time to adjust to the changed environmental conditions.

CAUTION

Risk to persons and risk of damage to equipment due to improper cable routing!

- ▶ Lay all cables and leads so that they do not present a tripping hazard.
- ▶ Do not place objects on the cables.

CAUTION

Risk to patient due to poor/lack of vision!

- ▶ Adjust the monitors and display elements so that they are highly visible to the user.

CAUTION

Interference of high frequency energy with the device!

- ▶ Do not use any mobile or portable equipment that emit high-frequency energy (e.g., mobile phones, cell phones, GSM phones) in the vicinity of the product.

CAUTION

Reduced illumination of the operating field due to damaged light guide!

- ▶ Wind combination cable with integrated light guide only with a bend radius of more than 7 cm.
- ▶ Do not manipulate combination cable with sharp objects.
- ▶ Do not apply mechanical force to the combination cable.
- ▶ Do not modify the light guide.
- ▶ Do not damage the fiber surfaces of the light guide ends.

Note

The product including the combination cable is fragile when bent, kinked, twisted, pulled, or pressed. These may damage the optical components (e.g., light guide fibers) and thus lead to inoperability.

Note

Handle the camera head, combination cable and sterile cover carefully. Do not exert any mechanical forces and avoid mechanical impacts. Please contact Service if the light output decreases.

Note

The equipotential plugs for all devices used must be connected with the equipotential tape IEC/EN/DIN 60601-1-1 or consistent with national standards.

Note

Ensure that the corresponding interconnection conditions, standards and respective national deviations are observed.

4.1 Combination with medical electrical equipment

The product can be used with components from other manufacturers if they fulfill the requirements of IEC/EN/DIN 60601-1 for the safety of medical devices. It is the operator's responsibility to check and make sure that the system is and remains fully operational. When using devices from various manufacturers and when using an endoscope and/or endoscopic accessory with electrical medical devices, the electrical shielding of the part must be indicated: Type CF defibrillation protected.

4.2 Inspections

Carry out the following inspection steps before reprocessing and immediately before use of the device:

Inspection of the glass surfaces

- ▶ Before each use, inspect the endoscope and sterile cover for damage such as rough surfaces, sharp edges or protrusions to avoid injuries to the patient.
- ▶ Perform a visual inspection of the glass surfaces. The surfaces must be clean and smooth.
- ▶ If there are impairments/damage, see Troubleshooting list.

Inspection of the light guide

- ▶ Before each use, inspect the light guide for damage and ascertain whether there is effective light transmission.
- ▶ Expose the distal end of the endoscope to the ambient light or move it to and fro in front of a bright ceiling lamp or similar (no cold light source) and hold the light guide connector close to the eye.

The transmission properties of the light guide fibers change with time. Black spots on the end surface indicate broken glass fibers or wear of the light guide fibers. If only isolated light guide fibers are affected, this is safe.

If the light output has noticeably reduced due to a too high breakage rate (see Figure 2 "Damaged light guide fibers"), it is recommended to send the endoscope for repair.

Aesculap®

EinsteinVision 3.0 Camera Heads PV631/PV632 and sterile covers for PV641SU/PV642SU

Never continue to use endoscopes with damaged light guides or combination cables.

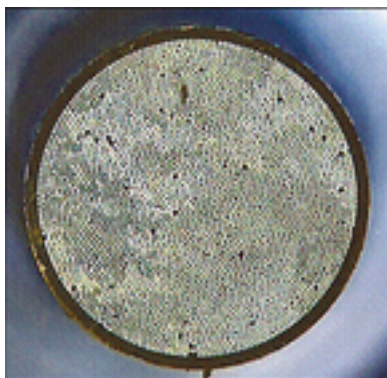
Fig. 1 Intact light guide fibers

Fig. 2 Damaged light guide fibers

First use

WARNING

Risk of injury and/or product malfunction due to incorrect operation of the electromedical system!

- Follow the instructions for use of all products required.

5. Working with the camera heads

5.1 General notes on use

The camera cable is fragile when bent, kinked, pulled, pressed, compressed and twisted. Observe the maximum bend radius of 7 cm and handle the cable with care. Do not make any changes to the cable and do not load with sharp-edged objects.

To avoid or reduce fogging of the optical surfaces, it can be helpful to supply the insufflation gas via another trocar than the camera trocar.

If the optical lens or tip of the sterile cover fogs over or gets soiled during the procedure, clean the lens if necessary with anti-fogging solution or sterile 0.9 % sodium chloride solution.

Depending on the condition and output of the light source and the light guide, combinations can occur which lead to the excessive development of heat at the proximal end of the light guide and damage the light cable, thus significantly reducing the transmission. In this case contact the manufacturer.

During cardiotoracic ablation, make sure that the ablation probe does not come into direct contact with the sterile cover.

5.2 System set-up

WARNING

Risk of infection due to unsterile parts!

- Ensure that unsterile components do not enter the sterile area.
- Process products and accessories that are delivered unsterile before use and only use them with sterile accessories.

WARNING

Risk of infection due to unsterile, soiled or contaminated products!

- Disinfect the camera head and camera cable after each use.
- Only use the camera head in the sterile area if it is in sterile condition. To do so, use the camera head with the correctly applied sterile cover.

WARNING

Risk of injury due to broken glass!

- Only mount the appropriate sterile cover on the camera head (0° on 0°, 30° on 30°).

WARNING

Risk to patient due to lack of vision of operating site!

- When mounting the sterile cover ensure that this noticeably engages on the camera head coupling.

WARNING

Risk of infection due to damaged sterile barrier!

- Do not use the sterile cover and EV3.0 Camera Head together with external tip heating.

WARNING

Risk to patient due to addition of leakage currents!

- If the product is used with medical electrical equipment and/or power-driven endoscope accessories, the leakage currents can multiply.

WARNING

Risk to patient and user from burns, sparks or explosion!

- When using a high frequency surgical device during an endoscopic procedure, observe the safety instructions of the corresponding instructions for use.

WARNING

Risk of injury due to improper use!

- Observe the intended use of the 3D endoscope and the sterile cover.
- Do not use the 3D endoscope and sterile cover as a lever.

CAUTION

Internal fogging due to humidity between endoscope/camera head and sterile cover!

- Mount the sterile cover on a dry endoscope/camera head only.
- If necessary, dry the endoscope with a sterile, lint-free cloth.

CAUTION

Interference of electromagnetic emissions with the image quality (e.g., minor banding, minor color changes on the monitor image)!

- Check the image quality if used in combination with additional peripheral devices (e.g., monitor, documentation system).

WARNING

Risk of infection due to unsterile cover!

- Check the integrity of the packaging of the sterile cover before opening.
- Do not use products from open or damaged sterile packaging.
- Do not use the product after the expiry date.

Note

In combination with the EV3.0 Camera Control Unit, the endoscope is classified as a type CF defibrillation-proof applied part.

Note

Before application of endoscopic high frequency surgery (HF surgery) prepare the patient accordingly.

Take measures to remove or avoid formation of combustible gases (e.g., gastro-intestinal tract/colonoscopy, urinary bladder/transurethral resection).

Connecting the accessories

DANGER

Risk of injury due to unapproved configuration using additional components!

- For all applied components, ensure that their classification matches that of the application component (e.g., Type BF or Type CF) of the respective device.

Combinations of accessories that are not mentioned in the present instructions for use may only be employed if they are specifically intended for the respective application, and if they do not compromise the performance and safety characteristics of the products.

All devices connected to the interfaces must also be proven to fulfill the relevant IEC standards (such as IEC 60950 for data processing devices and IEC/EN/DIN 60601-1 for medical electrical devices).

All configurations must fulfill the fundamental standard IEC/EN/DIN 60601-1. The person connecting the devices with each other is responsible for the configuration and must ensure that the fundamental standard IEC/EN/DIN 60601-1 or relevant national standards.

- Please address your B. Braun/Aesculap partner or Aesculap Technical Service with any inquiries in this respect; for a contact address, see Technical Service.

Note

The EV3.0 Camera Head can only be operated with the EV3.0 Camera Control Unit.

EinsteinVision 3.0 Camera Heads PV631/PV632 and sterile covers for PV641SU/PV642SU

5.3 Function checks

**Risk to patient due to incorrectly displayed image!**

- ▶ Before mounting the sterile cover, ensure that the endoscope tip and flat glass of the sterile cover are clean.
- ▶ Remove any contamination on the optical lens before the surgery. Clean the endoscope with a cloth soaked in alcohol (70 % ethanol) or with a neutral cleaning agent.
- ▶ Remove any contamination immediately during the surgery. To do this, clean the optical surfaces, e.g., with a damp sterile cloth and/or anti-fogging solution.
- ▶ Prior to application and after changing the settings (e.g., image rotation or activation of algorithms) check the correct display of the live image. If necessary, perform white balance.

**Risk of infection and injury due to defective sterile cover!**

- ▶ Before mounting the sterile cover, ensure that the endoscope tip and flat glass of the sterile cover are intact.
- ▶ After mounting the sterile cover, check it visually to ensure it is intact.
- ▶ At the end of the surgery, check the tip of the sterile cover to ensure it is intact.

**Risk of damage to eyes if light is beamed directly into eyes!**

- ▶ Do not look directly into the light exit at the distal end of the endoscope.
- ▶ When checking the light guide, do not plug the light guide connector into the light source.

**Risk to patient due to impaired or lack of 3D image impression!**

- ▶ Check the fit of the polarized 3D glasses before use and if necessary use a suitable glasses strap or have them fitted by an optician (e.g., if loupe glasses are also to be used in combination with the polarized 3D glasses).
- ▶ Do not store the polarized 3D glasses in an environment with high temperatures, e.g., near a radiator.
- ▶ Do not use worn, scratched, broken or damaged polarized 3D glasses.

**Risk to patient and user from functional failure or electric shock!**

- ▶ Use only dry components (e.g., plug connector for camera control unit, optical lens).

Note

Only perform the surgical procedure if all of the components being used are in a flawless condition.

Note

Before each use, after a power failure or any interruption, the functionality of all connected devices must be checked and all connections must be checked to see they are correct.

Note

The endoscope image must be focused, bright and clean at an appropriate working distance.

Note

In case of any noticeable mechanical damage (e.g., damage to the silicone tube), the camera head must not be used further with the combination cable.

Note

Do not continue to use damaged products.

Note

Differences in temperature and air humidity can cause the optical lens to fog up on insertion into the patient.

5.4 Safe operation

WARNING

Risk of burns due to development of heat at the tip of the endoscope and sterile cover!

- ▶ Make sure that the tip temperatures on the type of the endoscope and the sterile cover prescribed by IEC/EN/DIN 60601-2-18 are not exceeded.
- ▶ Do not put down the endoscope on the patient during use.
- ▶ In the patient's abdominal cavity, always leave a sufficient distance between the optical lens and the tissue surfaces and mucous membranes of the patient.
- ▶ Use automatic light intensity control or set the light source so that a bright, well-lit image is visible with the lowest possible light intensity.
- ▶ Switch off the light source if the illumination is no longer needed or if the endoscope is out of the patient for a longer period of time.
- ▶ Do not touch the light guide connector, distal end of the endoscope or the tip of the sterile cover during use or immediately after use.

WARNING

WARNING

WARNING

Risk of injury and/or malfunction!

- ▶ Always carry out a function check prior to using the product.

Risk to patient due to fogging of the optical lens and restricted view!

- ▶ To avoid or reduce fogging of the optical surfaces, it can be helpful to supply the insufflation gas via an other trocar than the camera trocar.

Risk of infection due to damaged sterile barrier!

- ▶ Before mounting the sterile cover, ensure that the light source is deactivated.
- ▶ Do not exert increased forces (e.g., rotation forces) to any part of the sterile cover (shaft and coupling).
- ▶ Before use, ensure that the sterile cover and the trocar system (10 mm) available in the OR are compatible.
- ▶ Only insert the endoscope and sterile cover in the trocar if it is possible to do so without friction and without exerting force.
- ▶ Do not insert an endoscope and sterile cover into any damaged trocar.
- ▶ Do not expose the sterile cover to excessive temperatures.
- ▶ Avoid direct contact between high frequency electrodes and the sterile cover.
- ▶ The camera head and sterile cover must be gripped in the holding arm only at the point intended for this.
- ▶ Do not allow the ablation probe and other applied parts to come into contact with the endoscope shaft and sterile cover during use.
- ▶ Do not use the sterile cover and EV3.0 Camera Head together with external tip heating.
- ▶ Only mount the appropriate sterile cover on the camera head (0° on 0°, 30° on 30°).
- ▶ Before mounting the sterile cover, ensure that the endoscope tip and flat glass of the sterile cover are intact.

Excessive temperatures in combination with light sources

Light sources, especially high-power light sources, emit large amounts of light energy and thermal energy. The light guide connector, distal end of the endoscope and the tip of the sterile cover can become extremely hot as a result of this.

Risks from use of light sources:

- Irreversible tissue damage or unwanted coagulation to the patient or user
- Burns or thermal damage to surgical equipment (e.g., surgical drapes, plastic materials)
- If the light source fails during use, this may endanger the patient. Therefore keep an operationally ready replacement light source on hand.

Safety precautions

- ▶ Do not illuminate the inside of the patient for longer than necessary with the light source.
- ▶ Use the automatic light intensity control or set the light source so that a bright, well-lit image is visible with the lowest possible light intensity.
- ▶ Do not allow the distal end of the endoscope or the light guide connector to come into contact with patient tissue or with combustible or heat-sensitive materials.
- ▶ Do not touch the distal end of the endoscope.
- ▶ Avoid contamination of the distal end surface or light emission surface.

Aesculap®

EinsteinVision 3.0 Camera Heads PV631/PV632 and sterile covers for PV641SU/PV642SU

WARNING

Risk of injury from burns and unwanted deep penetration and risk of damage to product!

- Only switch on high frequency current when the corresponding applied part (electrode) can be seen through the endoscope and there is no contact between the two.

WARNING

Risk of infection due to contamination of the sterile area!

- Fasten the combination cable covered by the sterile sleeve in the operating field so that it does not come into contact with unsterile surfaces.
- Fasten the combination cable so that the rear end of the combination cable, which is not covered by the sterile sleeve, does not enter the sterile area.

WARNING

Risk of injury due to misdirected leakage currents!

- Before thoracic or cardiothoracic surgery, deactivate implanted defibrillators (ICD).
- Before any kind of defibrillation, remove the camera head from the patient.

WARNING

Risk of fire and injury and risk to patient due to overheating at the endoscope light exit!

- Never leave the EV3.0 Camera Head unattended when light is being transmitted from the light source.
- During operation, do not place the EV3.0 camera head on the patient, in the operating field, or on temperature-sensitive or flammable objects.

WARNING

Risk to patient from gas embolism!

- Avoid overinsufflation e.g., with air or inert gas before high frequency surgery.

WARNING

Risk to patient from damage to the optical lens or sterile cover!

- Withdraw the endoscope and sterile sheath slowly from the used trocar.
- If a trocar with a manual opening valve is being used, only withdraw the endoscope and sterile sheath from the trocar when the valve is open.

CAUTION

Impairment of vision due to reflections!

- Before mounting the sterile cover, ensure that the light source is deactivated.

CAUTION

Malfunction in connection with magnetic resonance!

- Do not use the product in a magnetic resonance environment.

WARNING

Risk of infection due to unsterile cover!

- Replace the sterile cover if it falls on the ground or if it is touched by unsterile objects/persons.

Note

If the 3D display fails, the operation can be continued with 2D display.

Note

The current assignment of the camera head buttons is displayed on the screen as "Camera head information" after opening the camera menu.

Connection to camera control unit

- Insert the connection plug fully into the connection socket of the camera control unit until it engages in place.

Connection to the light source

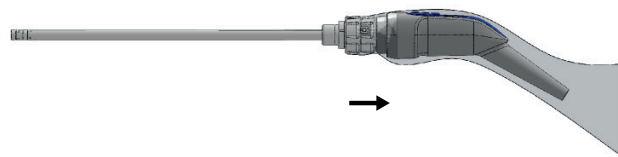
- Insert the connection plug fully into the light guide socket on the light source until it engages in place.

Mounting the sterile cover

- ▶ Mounting is carried out by one sterile and one unsterile person.
- ▶ The unsterile person opens the sterile cover packaging.
- ▶ The sterile person removes the sterile cover from the packaging and holds it while it is mounted on the shaft.
- ▶ The unsterile person does not hold the sterile cover from the front end of the shaft but only by the coupling.
- ▶ The unsterile person inserts the endoscope into the sterile cover shaft. Then they lock the locking mechanism by turning the coupling (70° counterclockwise).

- ▶ The unsterile person pulls the plastic sleeve of the sterile cover over the camera head and camera cable. When doing so, make certain that the plastic sleeve is only held at the every rear by the cardboard cover.

- ▶ Pull the sterile sleeve of the sterile cover over the camera cable as far as possible so that it is completely covered in the sterile area.
- ▶ Perform a visual inspection to see that the flat glass is in a flawless condition. When doing so, do not look directly into the light emitted from the exit of the light fiber.
- ▶ Then test the image.
- ▶ Only perform the surgical procedure if all of the EinsteinVision components being used are in a flawless condition. Do not hit or exert mechanical forces on the endoscope, camera head or sterile cover.

Removing the sterile cover

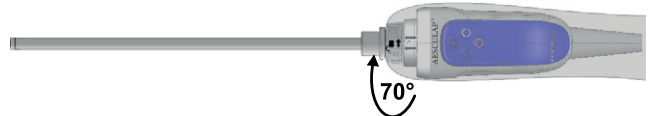
WARNING


Risk of infection/injury to patients and users!

- ▶ Use the sterile cover for EV3.0 Camera Head only once.
- ▶ Do not reprocess the sterile cover.

- ▶ Remove the sterile cover by unlocking the locking mechanism at the coupling. To do this, open the locking mechanism by turning the coupling (70° clockwise).

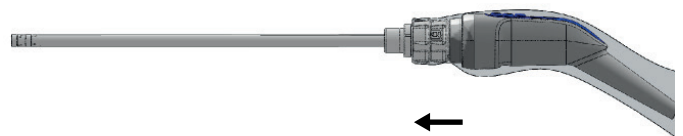- ▶ Then remove the sterile cover from the combination cable and camera head and dispose of according to the regulations.

EinsteinVision 3.0 Camera Heads PV631/PV632 and sterile covers for PV641SU/PV642SU

6. Reprocessing procedure

6.1 General safety instructions

Note

Adhere to national statutory regulations, national and international standards and directives, and local, clinical hygiene instructions for reprocessing.

Note

For patients with Creutzfeldt-Jakob disease (CJD), suspected CJD, or possible variants of CJD, observe the relevant national regulations concerning the reprocessing of products.

Note

The EV3.0 Camera Heads PV631/PV632 are not autoclavable and cannot be sterilized. Therefore, the EV3.0 Camera Heads PV631/PV632 and the sterile cover PV641SU/PV642SU must not be used in patients suspected of having Creutzfeldt-Jakob disease or who have already been diagnosed with Creutzfeldt-Jakob disease.

Note

When using the EV3.0 Camera Heads PV631/PV632 the sterile barrier is assured by the EV3.0 sterile cover. The sterile cover must be used for each application and must only be used once. It is recommended to perform wipe disinfection of the camera head and endoscope after use with the sterile cover (after removing the used sterile cover).

If severe contamination occurs due to the removal of the sterile cover or due to further handling in the OR, manual cleaning with immersion disinfection or mechanical alkaline cleaning and thermal disinfection of the EV3.0 Camera Heads PV631/PV632 can be performed.

Note

Successful reprocessing of this medical device can only be ensured through a validated reprocessing procedure. The operator/reprocessing technician is responsible for this.

Note

For up-to-date information about reprocessing and material compatibility, see also the Aesculap Extranet at www.aesculap-extra.net

WARNING

Risk of cross contamination during reprocessing!

If the device is potentially contaminated with the CJD agent, it is not suitable for reprocessing.

- ▶ Do not reprocess devices that are potentially contaminated.
- ▶ Devices that are potentially contaminated must be disposed of as clinical waste.

WARNING

Risk of infection for patient and/or user!

- ▶ Ensure that there is no residue from cleaning and disinfecting agents on the product.
- ▶ Avoid insufficient or incorrect cleaning and disinfection of the product and accessory.

CAUTION

Risk of damage to the device due to incorrect reprocessing!

- ▶ Do not autoclave the product.
- ▶ EV3.0 Camera Heads PV631/PV632: Do not sterilize.

CAUTION

Risk of defects, consequential damage or a shortened product service life!

- ▶ Follow and observe the manufacturer's requirements for reprocessing.
- ▶ Heated devices are very fragile. Handle camera head with care and do not subject it to any impact.
- ▶ Quick temperature changes can damage the product. Do not rapidly chill the product. Let it cool down without any additional cooling.

6.2 General information

Dried or affixed surgical residues can make cleaning more difficult or ineffective and lead to corrosion. Therefore the time interval between application and reprocessing should not exceed 6 h; also, neither fixating pre-cleaning temperatures >45 °C nor fixating disinfecting agents (active ingredient: aldehydes/alcohols) should be used.

Excessive neutralizing agents or basic cleaners may result in a chemical attack and/or fading and the laser marking becoming unreadable either visually or by machine in case of non stainless steel.

Residues containing chlorine or chlorides, e.g., in surgical residues, medicines, saline solutions, and in the service water used for cleaning, disinfection, and sterilization, will cause corrosion damage (pitting, stress corrosion) and result in damage to the products in case of non stainless steel. These must be removed by rinsing thoroughly with demineralized water and then drying.

Perform additional drying, if necessary.

Only process chemicals that have been tested and approved (e.g., VAH or FDA approval or CE mark) and which are compatible with the product's materials according to the chemical manufacturers' recommendations may be used for reprocessing the product. All the chemical manufacturer's application specifications must be strictly observed. Failure to do so can result in the following problems:

- Optical changes in materials, e.g., fading or discoloration of titanium or aluminum. For aluminum, the application/process solution only needs to be of pH >8 to cause visible surface changes.
- Material damage such as corrosion, cracks, fracturing, premature aging, or swelling.
- ▶ Do not use metal cleaning brushes or other abrasives that would damage the product surfaces and could cause corrosion.
- ▶ Further detailed advice on hygienically safe and material-/value-preserving reprocessing can be found at www.a-k-i.org, link to Publications, Red Brochure – Proper maintenance of instruments.

WARNING

Damage to the product due to improper reprocessing!

- ▶ Do not, under any circumstances, clean or disinfect the product in an ultrasonic cleaning bath.
- ▶ Use cleaning and disinfecting agents which are suitable and approved for the present product.
- ▶ Observe the manufacturer's cleaning and disinfecting instructions regarding concentration, temperature and exposure time.
- ▶ Only reconnect the system to power when all cleaned parts are completely dry.

WARNING

Danger for patient and the user due to high temperature after reprocessing!

- ▶ Ensure that the product has cooled down sufficiently before usage.

6.3 Dismantling prior to carrying out the reprocessing procedure

- ▶ Uncouple and completely remove the sterile cover from the camera head: Remove the sterile cover by unlocking the locking mechanism at the coupling. Then remove it and dispose of it according to regulations. The sterile cover is only intended for single use.
- ▶ Unplug the camera plug from the camera control unit.
- ▶ Unplug the light guide connector from the light source.

6.4 Preparation before cleaning

- ▶ Remove any visible surgical residues to the extent possible with a damp, lint-free cloth (damp with tap water or fully desalinated water). For this, use a soft cloth to not scratch the lens.
- ▶ Perform wipe disinfection.
- ▶ If the product is contaminated or potentially contaminated, transport the dry product in a sealed waste container for cleaning and disinfection within 6 hours.

6.5 Cleaning/disinfection

Product-specific safety instructions for the reprocessing procedure

DANGER

Risk of electric shock and fire hazard!

- ▶ Unplug the device before cleaning.
- ▶ Do not use flammable or explosive cleaning or disinfecting solutions.

CAUTION

Damage to the product due to inappropriate cleaning/disinfecting agents and/or excessive temperatures!

- ▶ Use cleaning and disinfecting agents according to the manufacturer's instructions. The cleaning and disinfecting agent must be approved for plastics and high-grade steel, and it must not attack softeners (e.g., silicone).
- ▶ Observe specifications regarding concentration, temperature and exposure time.
- ▶ Do not exceed the maximum permitted cleaning temperature of 60 °C.

Note

There may be significant discoloration of the EV3.0 Camera Head when using Cidex OPA. However, this discoloration has no influence on the functionality or safety of the product.

Aesculap®

EinsteinVision 3.0 Camera Heads PV631/PV632 and sterile covers for PV641SU/PV642SU

Validated cleaning and disinfection procedure

Validated procedure	Special features	Reference
Wipe disinfection for electrical devices without sterilization ■ EV3.0 Camera Heads PV631/PV632	none	Chapter Wipe disinfection for electrical devices without sterilization
Manual cleaning with immersion disinfection ■ EV3.0 Camera Heads PV631/PV632	■ Drying phase: Use a lint-free cloth or medical compressed air.	Chapter Manual cleaning/disinfecting and subsection: Chapter Manual cleaning with immersion disinfection
Valid for the US market: Manual cleaning ■ EV3.0 Camera Heads PV631/PV632	■ Valid for the US market	Chapter Manual cleaning/disinfecting and subsection: Chapter Valid for the US market: Manual cleaning
Mechanical alkaline cleaning and thermal disinfecting ■ EV3.0 Camera Heads PV631/PV632	■ Place the product in a tray that is suitable for cleaning (avoiding rinsing blind spots).	Chapter Mechanical cleaning/disinfecting and subsection: Chapter Mechanical alkaline cleaning and thermal disinfection

6.6 Wipe disinfection for electrical devices without sterilization

Phase	Step	T [°C/°F]	t [min]	Conc. [%]	Water quality	Chemicals
I	Wipe disinfection	RT	≥2	-	-	55 % propane-2-ol + 0,5 % quarternary ammonium chlorides*

* Super Sani-Cloth Disposable Wipes were used for validation.

Phase I

- ▶ Remove any visible residues with a disposable disinfectant wipe.
- ▶ Wipe all surfaces of the optically clean product with a fresh, disposable disinfectant wipe.
- ▶ Observe the specified application time (2 min minimum).

6.7 Manual cleaning/disinfecting

Note

Easily scratched surfaces. Risk of corrosion. Do not use metal cleaning brushes, metal objects or abrasive cleaners.

- After manual cleaning/disinfection, check visible surfaces visually for residues.
- Repeat the cleaning/disinfection process if necessary.

Manual cleaning with immersion disinfection

Phase	Step	T [°C/°F]	t [min]	Conc. [%]	Water quality	Chemicals
I	Cleaning	35-45/ 95-113	2-5	0.8	D-W	Enzymatic detergent Cidezyme/Enzol (Johnson & Johnson)
II	Intermediate rinse	RT (cold)	3 x 1	-	D-W	-
III	Disinfection	20-25/ 68-77	12	0.55	-	Ortho-phthalaldehyde solution Cidex OPA (Johnson & Johnson)
IV	Final rinse	RT (cold)	3 x 2	-	FD-W	-
V	Drying	RT	-	-	-	-

D-W: Drinking water

FD-W: Fully desalinated water (demineralized, low microbiological contamination: drinking water quality at least)

RT: Room temperature

Phase I

- Fully immerse the product in the cleaning/disinfecting solution for at least 2 to 5 minutes. Ensure that all accessible surfaces are moistened.
- After soaking in the cleaning solution, clean all external surfaces while the items are fully emerged in the cleaning solution – using a clean, lint-free, soft, non-sterile cloth or a soft bristled brush until all visible soiling is removed. Brush for at least 1 minute or until no more residues can be removed.

Phase II

- Fully immerse all parts in a tap water bath and thoroughly rinse all accessible surfaces twice for a minimum of one minute each.
- Use fresh water for each rinse.
- Allow water to drip off for a sufficient length of time.

Phase III

- Immerse all parts completely in disinfectant solution for at least 12 minutes. All accessible surfaces must remain immersed in the disinfectant solution bath throughout the entire disinfecting time
- Remove all adherent air bubbles from the component surfaces.

Phase IV

- Rinse/flush the product three times thoroughly (all accessible surfaces) for 2 minutes.
- Use fresh water for each rinse.
- Allow water to drip off for a sufficient length of time.

Phase V

- Thoroughly dry all parts with a clean, lint-free wipe or a lint-free surgical towel.
- Dry all accessible surfaces, especially the channels with medical-quality filtered compressed air ($p_{\max} = 0,5$ bar).

Aesculap®

EinsteinVision 3.0 Camera Heads PV631/PV632 and sterile covers for PV641SU/PV642SU

Valid for the US market: Manual cleaning

Phase	Step	T [°C/°F]	t [min]	Conc. [%]	Water quality	Chemicals
I	Cleaning	35-45/ 95-113	2-5	0.8	D-W	Enzymatic detergent Cidezyme/Enzol (Johnson & Johnson)
II	Rinse	RT (cold)	≥1	-	D-W	-
III	Final rinse	RT (cold)	3 x 2	-	FD-W	-
IV	Drying	RT	-	-	-	-

D-W: Drinking water

FD-W: Fully desalinated water (demineralized, low microbiological contamination: drinking water quality at least)

RT: Room temperature

Phase I

- Fully immerse the product in the cleaning/disinfecting solution for at least 2-5 minutes. Ensure that all accessible surfaces are moistened.
- After soaking in the cleaning solution, clean all external surfaces while the items are fully emerged in the cleaning solution – using a clean, lint-free, soft, non-sterile cloth or a soft bristled brush until all visible soiling is removed. Brush for at least 1 minute or until no more residues can be removed.

Phase II

- Fully immerse the instrument in a tap water bath for a minimum of 1 minute.

Phase III

- Thoroughly rinse the device three times for at least 2 minutes.
- Use fresh water for each rinse.
- Allow water to drip off for a sufficient period of time.

Phase IV

- Thoroughly dry all parts with a clean, lint-free wipe or a lint-free surgical towel.
- Dry all accessible surfaces, especially the channels with medical-quality filtered compressed air ($p_{\max} = 0,5$ bar).

6.8 Mechanical cleaning/disinfecting

Note

The washer/disinfector must be of tested and approved effectiveness (e.g., FDA approval or CE mark according to DIN EN ISO 15883).

Note

The washer/disinfector used for processing must be serviced and checked at regular intervals.

Note

Place the EV3.0 Camera Head in the tray in such a way to avoid rinsing blind spots.

Note

Only use the reprocessing tray for reprocessing. Do not place the camera head in the tray for transportation.

Inserting the camera head into the reprocessing tray

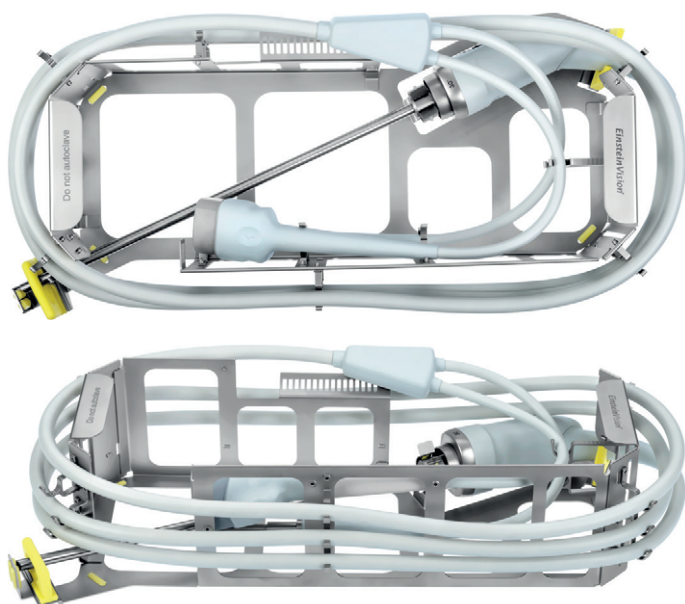
Fig. 3

- Put the camera head and endoscope sideways into the inlay of the reprocessing tray as pictured on the base of the inlay.
- Insert the cable into the hangers/holders/carriers and wrap it around the inlay with slight tension starting from the bottom to the top.
- Lay the Y-piece down on the specially designed tilted shelf and lead the thinner cables to the opposite holders on the inside of the inlay. The cables may not have direct contact to each other and to the camera head.
- Avoid rinsing blind spots.

Fig. 4

- Place the camera connector in the upper and the light guide connector in the lower holder. Make sure that the light guide connector is not touching the endoscope and that both plugs are properly secured.

Fig. 5

- Ensure that the tip of the endoscope is not touching the forward silicone protector.

Aesculap®

EinsteinVision 3.0 Camera Heads PV631/PV632 and sterile covers for PV641SU/PV642SU

Mechanical alkaline cleaning and thermal disinfection

Machine type: single-chamber cleaning/disinfection device without ultrasound

Phase	Step	T [°C/°F]	t [min]	Water quality	Chemistry/Note
I	Pre-rinse	Cold	2	D-W	-
II	Cleaning	55/131	10	D-W	Alkaline detergent Neodisher MediClean forte (Dr. Weigert) (5mL/L) pH >10
III	Intermediate rinse	-	1	D-W	In accordance with the standard cycle program for the cleaning and disinfecting machine
IV	Intermediate rinse	-	1	D-W	In accordance with the standard cycle program for the cleaning and disinfecting machine
V	Thermal disinfection	90/194	5	FD-W	-
VI	Drying	≤100/ 212	-	-	In accordance with the standard cycle program for the cleaning and disinfecting machine

D-W: Drinking water

FD-W: Fully desalinated water (demineralized, low microbiological contamination: drinking water quality at least)

- Check visible surfaces for residues after mechanical cleaning/disinfecting.
- Repeat the cleaning/disinfection process if necessary.

6.9 Material-compatible cleaning and disinfecting agents

Wipe disinfection

Meliseptol HBV wipes 50 % propan-1-ol (B. Braun)

Chlorine disinfection (Wetwipe)

Manual immersion disinfection

Stabimed fresh (B. Braun)

Anioxyde 1000 (Anios) – only for immersion disinfection of endoscope shaft

6.10 Inspection, maintenance and checks

- ▶ Inspect the product after each cleaning and disinfecting cycle to be sure it is: clean, functional, and undamaged.
- ▶ Check the product for damage or excessive heating.
- ▶ Set aside the product if it is damaged.

6.11 Packaging

- ▶ Always observe the instructions for use of the respective packaging.
- ▶ Ensure that the packaging will prevent a recontamination of the product.

6.12 Storage

CAUTION

Damage to product due to improper storage!

- ▶ Store the product protected from dust in a dry, well ventilated, and temperature-controlled room.
- ▶ Store the product protected from direct sunlight, high temperatures, high air humidity or radiation.
- ▶ Protect the product from direct UV light, radio-activity or strong electromagnetic radiation.
- ▶ Store the product separately or use containers in which it can be fixed.
- ▶ Always transport the product with care.

CAUTION

Damage to the optical lens due to inappropriate transport packaging!

- ▶ Use the tray for EV3.0 Camera Head for endoscopes only and not for instruments.
- ▶ Use the tray for EV3.0 Camera Head only within the clinic.
- ▶ Only place a disinfected EV3.0 Camera Head in the tray.

- ▶ Store sterile packaged single-use products in germ-proof packaging in a dust-protected, dry, dark and temperature-controlled room.

Camera head holder

The camera head holder PV636 is used for storing the camera head between uses.

- ▶ Carefully insert the camera head without sterile cover into the holder until the coupling for the camera head is fully positioned in the holder adapter.
- ▶ Ensure the shaft is not tilted and carefully insert the endoscope tip into the holder. Strong impacts on the endoscope tip can cause damage.
- ▶ Hang the camera cable loosely on the holder in big loops (bend radius of combination cable >7 cm).

WARNING

Risk of infection due to contaminated products!

- ▶ Only insert the camera head in the camera head holder if the camera head is dry and clean.

CAUTION

Damage to product due to improper handling!

- ▶ Do not strike the EV3.0 Camera Head against the camera head holder.
- ▶ Only store the EV3.0 Camera Head in the camera head holder with deactivated light source.
- ▶ When transporting the EV3.0 Camera Head over thresholds or uneven ground remove the EV3.0 Camera Head out of the camera head holder beforehand.

Aesculap®

EinsteinVision 3.0 Camera Heads PV631/PV632 and sterile covers for PV641SU/PV642SU

7. Service

Send damaged products to the manufacturer or authorized repair center. Authorized repair centers can be inquired at the manufacturer.

WARNING

Risk of infection due to soiled or contaminated products!

- ▶ Clean, disinfect and sterilize the product and/or any accessories thoroughly before sending.
- ▶ Remove the sterile cover from the camera head before sending.
- ▶ Select suitable and safe packaging (ideally the original packaging).
- ▶ Package the product in such a way that the packaging does not become contaminated.

Note

If for urgent reasons it is not possible to thoroughly clean/disinfect and sterilize the product, reprocess the product as far as possible and label it accordingly.

The specialized repair company may refuse to repair soiled or contaminated products for safety reasons.

The manufacturer reserves the right to return contaminated products to the sender.

8. Maintenance

This product is maintenance-free. It does not contain components which have to be maintained by the user or manufacturer.

However the manufacturer recommends that a qualified person or a hospital technician inspects the product regularly in a preventive safety inspection. The product should be inspected at least once every 12 months.

In case of preventative safety inspections, no special conditions and measures have to be observed.

When checking, perform the following tests pursuant to section 5 of IEC/EN/DIN 62353:

- 5.2 Visual inspection
- 5.3.2 Measuring of protective earth resistance (only applicable for devices of protection class I)
- 5.3.3 Leakage currents
- 5.3.4 Measurement of insulation resistance

The test results must be documented in full pursuant to section 6 of IEC/EN/DIN 62353 in a test report which can be taken from annex G of the standard.

After performing the preventive safety inspection, a functional test is to be carried out.

For relevant services contact your national B. Braun/Aesculap agency, see Technical Service.

9. Troubleshooting list

Problem	Possible cause	Remedy
Image cloudy	Glass surfaces soiled	Clean the glass surfaces according to the Cleaning and Disinfection section.
	Stubborn residue on the glass surfaces	Remove residues according to the Cleaning and Disinfection section; check water quality.
	Leaky, defective lens system	Send in endoscope system/camera head for repair.
Image too dark, too little illumination	Glass surfaces soiled	Clean the glass surfaces according to the Cleaning and Disinfection section.
	Light guide not attached to light source correctly	Check that the light guide is connected properly.
	Defective fiber optics	Check the glass surfaces.
	Light guide defective	Replace camera head.
Illumination has yellow tint	Soiled fiber optics	Clean the glass surfaces, see Cleaning/disinfection.
	Light guide dirty, defective	Check light guide (e. g. illuminate onto a white surface).
Corrosion, staining, discoloration	Inadequate cleaning (e.g., remaining protein residue)	Clean again, rub thoroughly if necessary.
	Excessive chloride concentration	Check water quality.
	Heavy metal ions and/or silicates, increased content of iron, copper, manganese in water	Check water quality; if applicable, only use deionized (fully desalinated) water.
	Too high concentration of mineral substances (e.g., calcium) or organic substances	
	Contaminated cleaning/disinfecting solutions, too frequent use	The cleaning and disinfecting solutions should be replaced regularly.
	Extraneous rust (e.g., resulting from previously damaged or non-corrosion-resistant instruments being reprocessed at the same time)	Check supply systems; in case of joint reprocessing, check material compatibility, for existing damage and avoid mutual contact.
	Contact corrosion	Avoid contact with other components.

9.1 Repairs

For repairs, contact the manufacturer or authorized repair center. Authorized repair centers can be inquired about from the manufacturer.

For a fast processing of your service requests, send in the product indicating:

- Article number (REF)
- Serial number (SN) / Batch number (LOT)
- Detailed description of defects

Note

If a component of the visualization system is replaced (e. g., service, upgrade), the system start-up process must be performed again. It is also recommended to perform a safety inspection after any maintenance or replacement of a device.

EinsteinVision 3.0 Camera Heads PV631/PV632 and sterile covers for PV641SU/PV642SU

10. Technical Service

WARNING

Risk of injury and/or malfunction!

- ▶ Do not modify the product.

CAUTION

Malfunction of the product due to damage in transport!

- ▶ Select suitable and safe packaging (ideally the original packaging).
- ▶ Keep the original packaging for possible returns in case of service.
- ▶ Package the product in such a way that the packaging does not become contaminated.

Note

Service and repair work must only be carried out by authorized trained personnel. For service and repairs, please contact your national B. Braun/Aesculap agency.

Modifications carried out on medical technical equipment may result in loss of warranty rights and applicable licenses.

Service addresses

Aesculap Technischer Service

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78532 Tuttlingen / Germany

Phone: +49 7461 95-1601

Fax: +49 7461 14-939

E-Mail: ats@aesculap.de

Or in the US

Aesculap Inc.

Attn. Aesculap Technical Services

615 Lambert Pointe Drive

Hazelwood

MO, 63042 USA

Aesculap Repair Hotline

Phone: +1 800 214-3392

Fax: +1 (314) 895 -4420

Other service addresses can be obtained from the address indicated above.

11. Accessories/Spare parts

Art. no.	Designation
PV641SU	Sterile Cover for camera head 0°
PV642SU	Sterile cover for Camera Head 30°
PV636	EV3.0 Camera Head holder
JF444R	Tray for 3D camera head
JF454R	Lid for tray JF444R

12. Technical data

Classification acc. to Directive 93/42/EEC

Art. no.	Designation	Class
PV631	EV3.0 Camera Head 10 mm 0°	Ila
PV632	EV3.0 Camera Head 10 mm 30°	Ila
PV641SU	Sterile cover for Camera Head 0°	Ila
PV642SU	Sterile cover for Camera Head 30°	Ila

12.1 Camera head

Protection rating	IPX7
Viewing direction	0°/30°
Diameter	PV631/PV632: 10 mm (in combination with sterile cover PV641SU/642SU)
Working length	PV631/PV632: 315 mm (in combination with sterile cover PV641SU/PV642SU)
Sensor format	Native full HD 1/3"
Scanning system	Progressive scan
Brightness	Automatic
Weight (excl. cable)	PV631/PV632: 624 g
Camera cable length	4 118 mm
Light guide length	4 118 mm
Dimensions (L x W x H)	PV631: 491 mm x 56.5 mm x 47.6 mm PV632: 492 mm x 56.5 mm x 47.6 mm






Applied part	Type CF, defibrillation-proof (in combination with camera control unit PV630)
Conforming to standard	IEC/EN/DIN 60601-1
EMC	IEC/EN/DIN 60601-1-2
CISPR 11	Class A

12.2 Sterile cover



Viewing direction	0°/30°
External diameter	10 mm
Working length	315 mm
Total length (without sterile sleeve) 0°/30°	366 mm/367 mm
Weight	90 g
Sterile sleeve dimensions (L x D)	3 000 mm x 255 mm

12.3 Ambient conditions

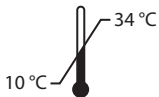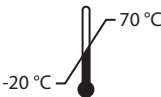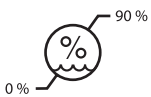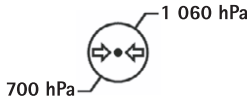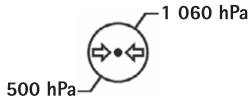Camera Head

	Operation	Storage and transport
Temperature		
Relative humidity		
Atmospheric pressure		

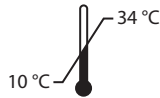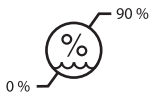Sterile cover

	Operation
Temperature	
Relative humidity	
Atmospheric pressure	

13. Disposal

Note

The user institution is obliged to reprocess the product before its disposal, see *Reprocessing procedure*.

	Adhere to national regulations when disposing of or recycling the product, its components and its packaging! The recycling pass can be downloaded from the Extranet as a PDF document under the respective article number. (The recycling pass includes disassembling instructions for the product, as well as information for proper disposal of components harmful to the environment.) Products carrying this symbol are subject to separate collection of electrical and electronic devices. Within the European Union, disposal is taken care of by the manufacturer as a free-of-charge service.

- Detailed information concerning the disposal of the product is available through your national B. Braun/Aesculap agency, see Technical Service.

Note

Incorrect disposal can lead to environmental damage.

14. Distributor in the US/Contact in Canada for product information and complaints

Aesculap Inc.
3773 Corporate Parkway
Center Valley, PA, 18034,
USA

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