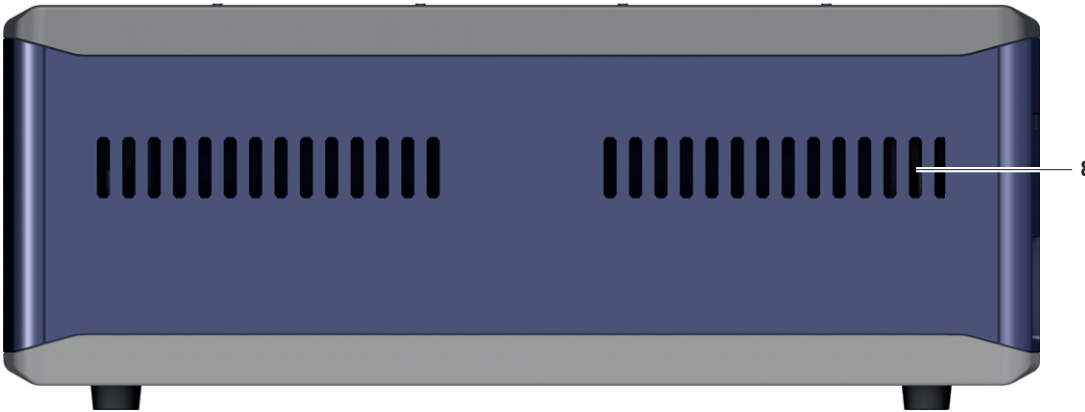
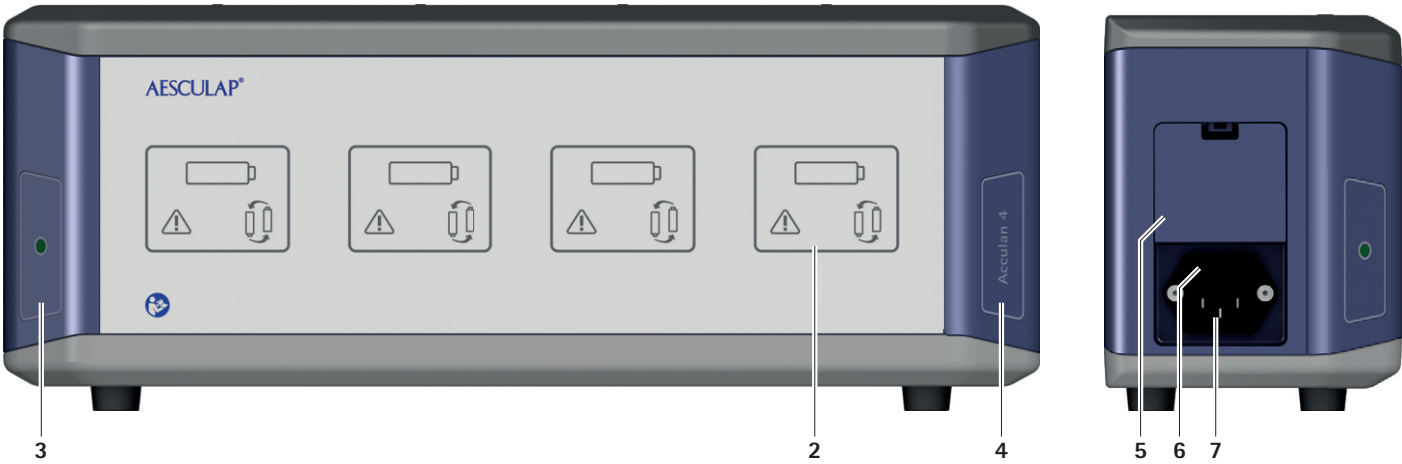
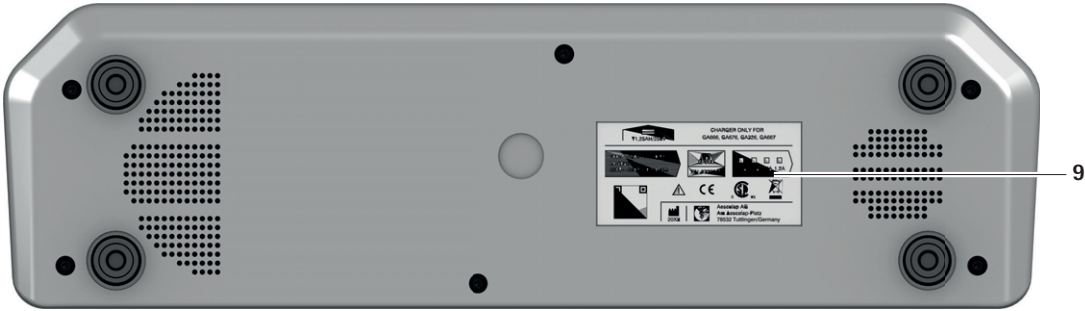
Aesculap® Acculan 4

- (GB)** Instructions for use/Technical description
- (USA)** Acculan 4 Charger GA320
- (D)** Gebrauchsanweisung/Technische Beschreibung
Acculan 4 Ladegerät GA320
- (F)** Mode d'emploi/Description technique
Chargeur Acculan 4 GA320
- (E)** Instrucciones de manejo/Descripción técnica
Cargador Acculan 4 GA320
- (I)** Istruzioni per l'uso/Descrizione tecnica
Caricabatterie Acculan 4 GA320
- (P)** Instruções de utilização/Descrição técnica
Carregador Acculan 4 GA320
- (NL)** Gebruiksaanwijzing/Technische beschrijving
Acculan 4 laadapparaat GA320

- (S)** Bruksanvisning/Teknisk beskrivning
Acculan 4-laddare GA320
- (RUS)** Инструкция по применению/Техническое описание
Зарядное устройство Acculan 4 GA320
- (CZ)** Návod k použití/Technický popis
Nabíječka Acculan 4 GA320
- (PL)** Instrukcja użytkowania/Opis techniczny
Ładowarka Acculan 4 GA320
- (SK)** Návod na použitie/Technický opis
Nabíjačka Acculan 4 GA320
- (TR)** Kullanım Kılavuzu/Teknik açıklama
Acculan 4 şarj cihazı GA320
- (KR)** 사용 설명서 / 기술 설명
Acculan 4 충전기 GA320

Aesculap® Acculan 4

Acculan 4 Charger GA320

Legend

- 1 Charging bays
- 2 Indicator fields for charging bays
- 3 "Power On" indicator (lit up green)
- 4 Product name Acculan 4
- 5 Cover/interlock connections for service. Use is reserved solely for the manufacturer or authorized service.
- 6 Fuse holder
- 7 Device socket
- 8 Ventilation slots
- 9 Type plate

Symbols on product and packages

	Caution Refers to the necessity for the user to review the instructions for use for important safety-related information, such as warning notices and precautionary measures, which cannot be attached to the medical product itself for a variety of reasons.
	Follow the instructions for use
	"Power On" indicator (lit up green)
	Indicator/symbol of the charging bays Charging progress indicator (illuminated in green)
	Indicator/symbol of charging bays "Action request" symbol (illuminated in orange)
	Indicator/symbol of charging bays "Battery change recommended" symbol (illuminated in orange)
	Fuse
	Alternating current
	Manufacturer's serial number
	Manufacturer and date of manufacture (year)

Labeling of electric and electronic devices according to Directive 2012/19/EU (WEEE), see chapter on disposal

Date of manufacture

Manufacturer's batch designation

Manufacturer's article number

Temperature limit
The temperature limit values to which the medical product can be safely subjected are designated.

Humidity, limit
Designates the humidity range to which the medical product can be safely subjected

Air pressure, limit
Designates the atmospheric pressure to which the medical product can be safely subjected

Contents

1.	Applicable to	3
2.	Intended use	3
2.1	Task/Function	3
2.2	Application Environment	3
2.3	Set-up location	3
3.	Safe handling	3
4.	Product description	4
4.1	Scope of supply	4
4.2	Components required for operation	4
4.3	Operating principle	4
5.	Preparation	5
6.	Working with the device	6
6.1	System set-up	6
6.2	Safe operation	6
7.	Validated reprocessing procedure	7
7.1	General safety notes	7
7.2	General information	7
7.3	Preparation before cleaning	7
7.4	Cleaning/disinfection	8
7.5	Wipe disinfection for electrical devices without sterilization ..	8
7.6	Inspection, maintenance and checks	9
8.	Maintenance	9
9.	Troubleshooting list	9
10.	Replacing fuses	10
11.	Technical Service	10
12.	Accessories/Spare parts	10
13.	Technical data	10
13.1	Performance data, information about standards	10
13.2	Ambient conditions	11
14.	Disposal	11
15.	Distributor in the US/Contact in Canada for product information and complaints	11

1. Applicable to

- For item-specific instructions for use and information on material compatibility, see also the Aesculap Extranet at <https://extranet.bbraun.com>

2. Intended use

2.1 Task/Function

The charger is used to charge and monitor the Aesculap Acculan 4 Li-Ionen battery GA346. The Aesculap Acculan 3Ti NiMH battery GA676 and GA666 can also be charged and monitored with this charger.

2.2 Application Environment

The charger is approved for operation in confined spaces, outside of the patient environment, in the non-sterile zone and outside of the explosion risk zone (such as areas with pure oxygen or anesthesia gases).

2.3 Set-up location

The charger must be placed on a table or shelf that ensures stability.
The charger may not be subjected to direct sunlight or moisture.

3. Safe handling

CAUTION

Federal law restricts this device to sale by, or on order of a physician!

DANGER

Risk of death by electric shock!

- Do not open the product.
- Connect the product only to a grounded power supply.

WARNING

Risk of injury and damage to property due to improper handling of the product!

- Observe the instructions for use of the Acculan 4 battery and of the Acculan3 Ti rechargeable battery.
- Follow the instructions for use of all products used.

WARNING

Risk of injury and material damage if this product is not used as intended!

- Use the product only in accordance with the intended purpose.

- Inspect the new product after removing its transport packaging and prior to first use to ensure it is in good working order.
- See "Information on electromagnetic compatibility (EMC) for Acculan 4 charger" TA022461; see Aesculap Extranet at <https://extranet.bbraun.com>

Aesculap® Acculan 4

Acculan 4 Charger GA320

Note

The Acculan 4 Charger GA320 fulfills the requirements of CISPR 11, Class A.

- ▶ To prevent damage caused by improper setup or operation, and to not compromise the manufacturer warranty and liability:
 - Use the product only according to these instructions for use.
 - Follow the safety and maintenance instructions.
 - Only combine Aesculap products with each other.
 - Follow the application advisories acc. to standard, see Extracts from relevant standards.
- ▶ Ensure that the product and its accessories are operated and used only by persons with the requisite training, knowledge, or experience.
- ▶ Keep the instructions for use accessible for the user.
- ▶ Always adhere to applicable standards.
- ▶ Ensure that the electrical installation of the room is consistent with the requirements of IEC/DIN EN.
- ▶ Do not pull any of the cables by the actual cable but only by the connector when disconnecting the device from power.
- ▶ Do not operate the product in explosion-hazard areas.
- ▶ Use the product in non-sterile zones only.
- ▶ Do not use the product if it is damaged or defective.
- ▶ Set aside the product if it is damaged.
- ▶ Follow the instructions for use of the corresponding Acculan battery.

4. Product description

4.1 Scope of supply

Art. no.	Designation
GA320	Acculan 4 Charger
TA014535	Instructions for use, Acculan 4 Charger GA320
C63401	Instructions for use, Validated Reprocessing Procedure, German
C63402	Instructions for use, Validated Reprocessing Procedure, English

4.2 Components required for operation

- Acculan 4 Charger
- Power cable, see accessories

4.3 Operating principle

The Acculan 4 Charger is designed for a mains voltage range of 100 V~ to 120 V~ (0.8 A to 0.7 A) or 220 V~ to 240 V~ (0.5 A to 0.4 A). The mains voltage is converted into a protective extra-low voltage. To make the charger ready for operation, connect it to mains power using the mains power cord. No other action is required to activate or deactivate the console. The charger performs a self-test, which takes approx. 3 seconds. During this test, all indicator elements at the front panel of the charger are illuminated.

Charging bays

The Acculan 4 Charger has four separate, independent slots that are compatible with the Acculan 4 battery GA346, and with the Acculan battery GA676 and GA666. Each charging bay has an indicator field on the charger front panel assigned to it.

Automatic start of charging process

Charging starts automatically as soon as a battery is inserted in a charging bay. The charging time in all four charging bays depends on the charge state and capacity of the batteries.

Li-Ion battery GA346

Charging principle

The batteries are charged very gently. The charge state of the battery is continuously monitored throughout charging. Monitoring of the charging curve ensures that the battery is fully charged without overcharging. Additionally, the battery temperature and charging time are monitored as well.

The Li-Ionen batteries GA346 undergo an automatic battery check. This check assesses the battery condition.

If the "Battery change recommended" symbol lights up, replacement of the battery is recommended. This battery can still be used. However, it must be taken into account that an intraoperative battery change may be necessary due to the reduced performance, or a safety function may activate.

Battery temperature

Each charging bay is equipped with a fan, which works independently of the battery temperature. The charging of the battery is aborted at a battery temperature of $>45^{\circ}\text{C}$. Excessive battery temperature is indicated by the "Action request" symbol illuminated in indicator field. When the battery temperature drops, the "Action request" symbol in indicator field **2** still remains illuminated. The battery must be removed from the respective charging bay and plugged in again, to restart charging.

Charging time

Charging is stopped as soon as the maximum charging time is reached. The permanent monitoring of the battery also allows detecting faults in the battery block. When a battery fault or a charging fault is detected, this is indicated by the "Action request" symbol illuminated in indicator field **2** and charging is aborted.

NiMH battery GA676 and GA666

Charging principle

The batteries are charged very gently with constant current pulses. The charge state of the battery is continuously monitored throughout charging. Monitoring of the charging curve ensures that the battery is fully charged without overcharging. Additionally, the battery temperature and charging time are monitored as well.

The NiMH GA676 and GA666 batteries undergo an automatic battery check. This check assesses the battery condition.

If the "Battery change recommended" symbol lights up, replacement of the battery is recommended. This battery can still be used. However, it must be taken into account that an intraoperative battery change may be necessary due to the reduced performance, or a safety function may activate.

Battery temperature

Each charging bay is equipped with a fan, which works independently of the battery temperature. The charging of the battery is aborted at a battery temperature of $>47^{\circ}\text{C}$. Excessive battery temperature is indicated by the "Action request" symbol illuminated in indicator field **2**. When the battery temperature drops, the "Action request" symbol in indicator field **2** still remains illuminated. The battery must be removed from the respective charging bay and plugged in again, to restart charging.

Charging time

Charging is stopped as soon as the maximum charging time is reached. The permanent monitoring of the battery also allows detecting faults in the battery block. When a battery fault or a charging fault is detected, this is indicated by the "Action request" symbol illuminated in indicator field **2** and charging is aborted. After a successfully completed charging process, the charger switches to trickle charging.

Note

The Acculan 4 Charger has a function for maintaining the maximum charge when a battery is left in the charger after complete recharging.

5. Preparation

Non-compliance with the following instructions will preclude all responsibility and liability in this respect on the part of Aesculap.

- ▶ When setting up and operating the product, adhere to:
 - The national regulations for installation and operation
 - Instructions for use according to IEC/DIN EN regulations.
 - national regulations on fire and explosion protection.

Note

The Acculan 4 Charger GA320 may not be subjected to direct sunlight or moisture.

Note

After being installed and put into operation, the Acculan 4 Charger GA320 must not be transported or moved to a different set-up location.

- ▶ Make certain that the ventilation slots in the bottom of the housing and in the rear panel of the charger are not covered.
- ▶ Make certain the device is set up on a sufficiently stable support (table, shelf).
- ▶ Make sure that the device socket is accessible to the user.
- ▶ Make sure that the charging bays are accessible to the user.

Aesculap® Acculan 4

Acculan 4 Charger GA320

6. Working with the device

6.1 System set-up

Connecting to the power supply

Risk of death by electric shock!

- ▶ **Connect the product only to a grounded power supply.**

- ▶ Carry out a visual inspection before connecting the product to the mains power:
 - Check the mains power cord for any damage.
 - Check the charger for any damage (e.g. bent, short-circuited contacts of the charging bays).
- ▶ Do not use the product if it is damaged, defective or dirty.
- ▶ Set aside the product if it is damaged.
- ▶ Clean product if dirty, see Validated reprocessing procedure.
- ▶ Plug the mains power cord into the device connector 7.
- ▶ Plug in the mains plug at the mains power socket.

After being connected to the mains power, the charger performs a self-test. During this time the "Power On" 3 indicator and the indicator fields for charging bays 2 must be illuminated on the front panel of the device.

- ▶ Carry out a visual inspection of the "Power On" indicator 3 and the indicator fields for charging bays 1.

The indicator fields for the charging bays 2 go out and only the "Power On" 3 indicator remains lit.

The charger is ready.

Note

In case of any deviation from the sequence described above, see Troubleshooting list

Putting out of operation

Note

The Acculan 4 Charger GA320 does not have a power ON/OFF switch. The all-pole disconnection of the product from the main power supply is only guaranteed when the power cord/mains connection line is unplugged.

- ▶ Unplug the power cord from the mains power socket.
- ▶ Unplug power cord from the device connector 7.

6.2 Safe operation

WARNING

Risk of injury and material damage due to inappropriate use of batteries!

- ▶ Only insert batteries that are indicated in the Intended Use.
- ▶ Insert batteries in their correct position.

WARNING

Risk of injury and material damage due to damaged batteries!

- ▶ Check batteries for any damage before charging.
- ▶ Do not use or recharge damaged batteries.

CAUTION

Total battery failure can occur by leaving the batteries in the charger while disconnected from the power supply!

- ▶ Always leave charger connected to the power supply. Indicator lamp for readiness must illuminate.
- ▶ Always remove batteries from the charger when it is disconnected from the power supply.

CAUTION

Loss of capacity/performance due to prolonged storage of discharged batteries!

- ▶ After a successfully completed charging process, leave the battery in the operational charger (trickle charging).

- ▶ As excessive battery temperatures will prolong the charging time, allow batteries to cool down prior to charging.
- ▶ To remove the battery, if necessary use a cloth as the batteries become heated during the charging process.

Charge battery GA346, GA666 and GA676

- Insert the battery into the charging bay with the contacts pointing downwards, so that the marking on the battery **A** corresponds with the marking on the charging bay **B**, see Fig. 1.

The "Charging progress" indicator is illuminated in green.

The charging process starts automatically.

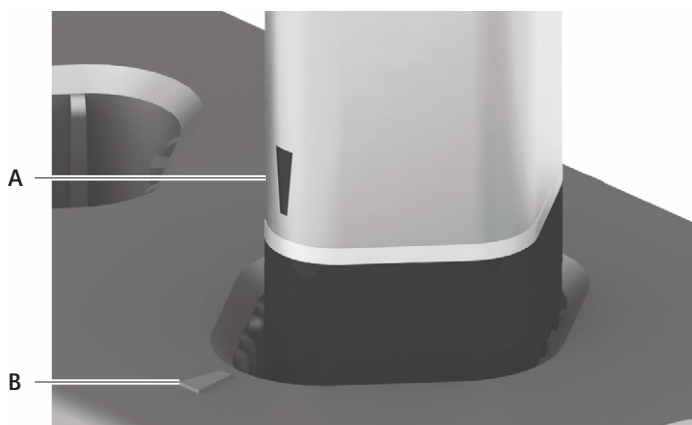
Fig. 1 Insert battery into charging bay

Indicator fields for charging bays

"Charging progress" indicator

The "Charging progress" indicator is illuminated in green as soon as a battery is inserted into the corresponding charging bay. It indicates the charging status of the battery. At full illumination, the charging process is complete and the battery can be removed from the charging bay.

"Action request" indicator

The "Action request" indicator is illuminated in orange when a charging process cannot be properly carried out. The "Charging progress" indicator goes out simultaneously.

Possible causes: bad contacts or excessive battery temperature during charging. see Troubleshooting list

"Battery change recommended" indicator

The "Battery change recommended" indicator is illuminated in orange in addition to the "Charging progress" indicator when the charging process is complete. The corresponding battery has failed the automatic battery check. This battery can nevertheless still be used. However, it must be taken into account that an intraoperative battery change may be necessary due to the reduced performance.

7. Validated reprocessing procedure

7.1 General safety notes

Note

Adhere to national statutory regulations, national and international standards and directives, and local, clinical hygiene instructions for reprocessing.

Note

For patients with Creutzfeldt-Jakob disease (CJD), suspected CJD, or possible variants of CJD, observe the relevant national regulations concerning the reprocessing of products.

Note

It should be noted that successful reprocessing of this medical device can only be guaranteed following prior validation of the reprocessing method. The operator/reprocessing technician is responsible for this.

The recommended chemistry was used for validation.

7.2 General information

Only process chemicals that have been tested and approved (e.g. VAH or FDA approval or CE mark) and which are compatible with the product's materials according to the chemical manufacturers' recommendations may be used for processing the product. All the chemical manufacturer's application specifications must be strictly observed. Failure to do so can result in the following problems:

- Visual material changes
- Material damage (such as corrosion, cracks, breaks, premature aging or swelling).

7.3 Preparation before cleaning

- The battery must be removed from the respective charging bay.
- Unplug the power cord from the mains power socket.
- Unplug power cord from the device connector 7.

Aesculap® Acculan 4

Acculan 4 Charger GA320

7.4 Cleaning/disinfection

Product-specific safety instructions for the reprocessing procedure

DANGER

Risk of electric shock and fire hazard!

- ▶ Unplug the power cord before cleaning.
- ▶ Do not use flammable or explosive cleaning or disinfecting solutions.
- ▶ Ensure that no fluids will penetrate the product.

CAUTION

Damage to, or destruction of the product caused by mechanical cleaning/disinfection!

- ▶ Only clean and disinfect the product manually.
- ▶ Do not sterilize the product under any circumstances.

CAUTION

Risk of damage to the product due to inappropriate cleaning/disinfecting agents!

- ▶ Only use cleaning/disinfecting agents approved for surface cleaning. Follow the manufacturer's instructions for the respective cleaning/disinfecting agent.
- ▶ Do not introduce any fluid into the charging bays.

7.5 Wipe disinfection for electrical devices without sterilization

Phase	Step	T [°C/°F]	t [min]	Conc. [%]	Water quality	Chemicals
I	Wipe disinfection	RT	≥1	-	-	Meliseptol Wipes sensitive

RT: Room temperature

Phase I

- ▶ Remove any visible residues with a disposable disinfectant wipe.
- ▶ Wipe all surfaces of the optically clean product with a fresh, disposable disinfectant wipe.
- ▶ Observe the specified application time (1 min minimum).
- ▶ Clean the contacts in the charging bay with isopropanol or ethyl alcohol and a cotton swab. Do not use corrosive chemicals.

7.6 Inspection, maintenance and checks

- Inspect the product after each cleaning and disinfecting cycle to be sure it is clean, functional and undamaged.
- Set aside the product if it is damaged.

8. Maintenance

For technical service, please contact your national B. Braun/Aesculap agency, see Technical Service.

9. Troubleshooting list

- Have defective products repaired by Aesculap Technical Service, see Technical Service.

Issue	Detection	Cause	Remedy
Charger not working	"Power On" indicator not illuminated	Mains power cord not plugged in	Plug mains power cord into the device connector 7 on the charger and into the socket of the mains power supply.
		Power cord faulty	Replace power cord.
		Fuses blown	Replace fuses, see Replacing fuses.
		Charger defect	Replace product.
Defective indicator	During the switch-on check, not all indicator elements are illuminated ("Power On" indicator and indicator fields for charging bays)	Charger defect	Replace product.
Battery does not charge	Battery inserted, charging progress indicator not illuminated	Charger contacts dirty	Clean the contacts, see Validated reprocessing procedure.
		Battery defective	Repeat the charging process in another charging bay. If the problem occurs again, replace product.
		Charger defect	The charging process does not start up either with a different battery. Replace product.
	Battery inserted, "Action request" indicator illuminated	Excessive battery temperature detected during charging	Remove the battery from the charging bay, let it cool down and restart charging. If the problem occurs again, replace battery.
		Battery defective	Repeat the charging process in another charging bay. If the problem occurs again, replace battery.
		Charger defect	The charging process does not start up either with a different battery. Replace product.

Aesculap® Acculan 4

Acculan 4 Charger GA320

10. Replacing fuses

DANGER

Risk of death by electric shock!

- Pull the mains plug before replacing the fuse sets.

Mandatory fuse set: 2x – T 1.25 AH, 250 V

- Use a small screwdriver to release the clip on the fuse holder.
- Pull out the fuse holder.
- Replace both fuse sets.
- Reinstall the fuse holder so that it audibly clicks into place.

Note

If the fuses burn out frequently, the device is faulty and should be repaired, see *Technical Service*.

11. Technical Service

DANGER

Danger to life of patients and users if the product malfunctions and/or protective measures fail or are not used!

- Do not modify the product.

Modifications carried out on medical technical equipment may result in loss of guarantee/warranty rights and forfeiture of applicable licenses.

- For service and repairs, please contact your national B. Braun/Aesculapagency.

Service addresses

Aesculap Technischer Service

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78532 Tuttlingen / Germany

Phone: +49 (7461) 95 -1601

Fax: +49 (7461) 14 -939

E-Mail: ats@aesculap.de

Or in the US:

Aesculap Inc.

Attn. Aesculap Technical Services

615 Lambert Pointe Drive

Hazelwood

MO, 63042 USA

Aesculap Repair Hotline

Phone: +1 (800) 214 -3392

Fax: +1 (314) 895 -4420

Other service addresses can be obtained from the address indicated above.

12. Accessories/Spare parts

Art. no.	Approvals	Color	Length [m]
TE780	Europe	Black	1.5
TE730	Europe	Black	5.0
TE734	Great Britain	Black	5.0
TE735	USA, Canada, Japan	Gray	3.5

Art. no.	Designation
TA020112	Fuse link T 1.25 AH, 250 V
TA014534	Instructions for use of Acculan 4 Charger GA320 (A4 for ring binder)

13. Technical data

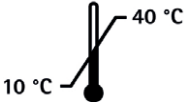
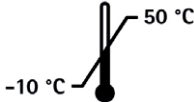
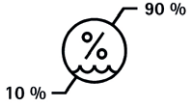
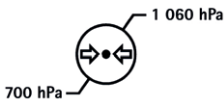
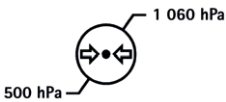
13.1 Performance data, information about standards

Acculan 4 Charger GA320

Classification according to Directive 93/42/EWG	I
Protection class according to IEC/DIN EN 60601-1	I
Housing protection rating according to IEC/DIN EN 60529	IP21
For charging Aesculap battery types	GA346 and GA666, GA676
Mains voltage ranges (Current intake)	100 V~ to 120 V~ (0.8 A to 0.7 A) / 220 V~ to 240 V~ (0.5 A to 0.4 A)
Frequency	50 Hz to 60 Hz
Charging voltage/Output voltage	max. 15 V
Charging current/Output current	max. 1.2 A ± 0.3 A
Operational mode	Continuous operation
Device protection	T 1.25 AH, 250 V Form factor: 5 mm x 20 mm
Weight without battery	approx. 2.3 kg

Dimensions (L x W x H) without battery	410 mm x 120 mm x 160 mm ±5 mm
Dimensions (L x W x H) with battery	410 mm x 252 mm x 160 mm ±5 mm
Conforming to standard	IEC/DIN EN 60601-1
EMC	IEC/DIN EN 60601-1-2
CISPR11	Class A

13.2 Ambient conditions

	Operation	Storage and transport
Temperature		
Relative humidity		
Atmospheric pressure		

Note

The atmospheric pressure of 700 hPa corresponds to a maximum operating altitude of 3 000 m.

14. Disposal

Note

The user institution is obliged to process the product before its disposal, see Validated reprocessing procedure.

Adhere to national regulations when disposing of or recycling the product, its components and its packaging!

The recycling pass can be downloaded from the Extranet as a PDF document under the respective article number. (The recycling pass includes disassembling instructions for the product, as well as information for proper disposal of components harmful to the environment.)

Products carrying this symbol are subject to separate collection of electrical and electronic devices. Within the European Union, disposal is taken care of by the manufacturer as a free-of-charge service.

- Detailed information concerning the disposal of the product is available through your national B. Braun/Aesculap agency, see Technical Service.

15. Distributor in the US/Contact in Canada for product information and complaints

Aesculap Inc.
3773 Corporate Parkway
Center Valley, PA, 18034,
USA

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