

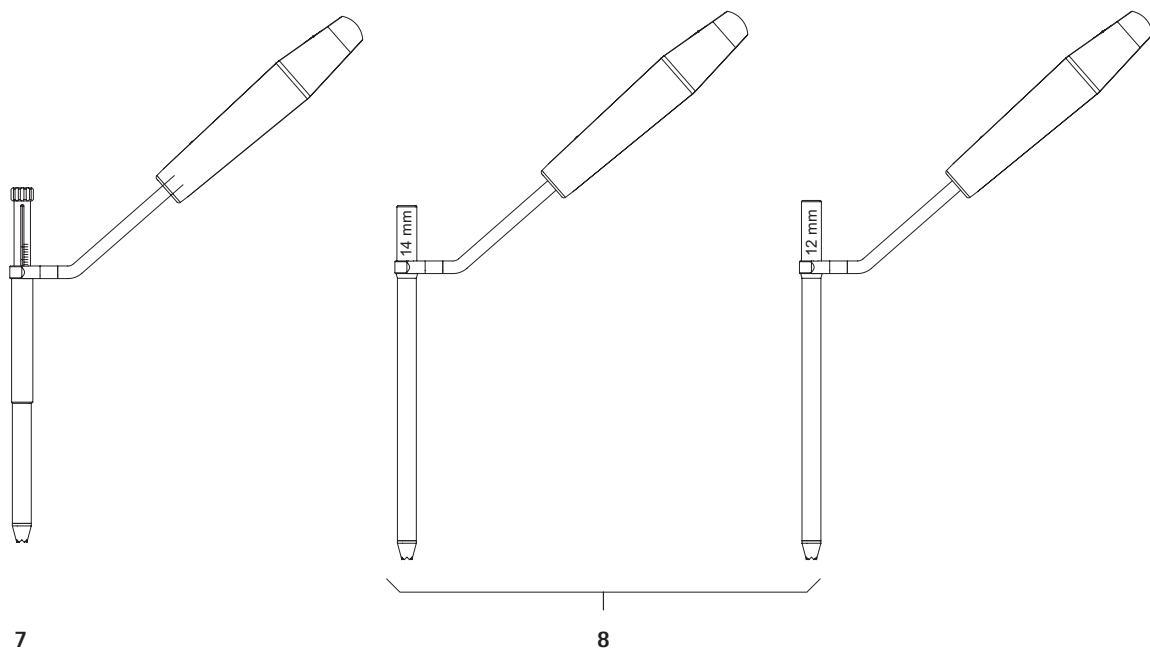
Aesculap® S⁴® Cervical

Aesculap Spine

- (GB)** Instructions for use/Technical description
- (USA)** S⁴® Cervical Instruments
- (D)** Gebrauchsanweisung/Technische Beschreibung
S⁴® Cervical Instrumente
- (F)** Mode d'emploi/Description technique
Instruments S⁴® Cervical
- (E)** Instrucciones de manejo/Descripción técnica
Instrumentos S⁴® Cervical
- (I)** Istruzioni per l'uso/Descrizione tecnica
Strumenti S⁴® Cervical
- (P)** Instruções de utilização/Descrição técnica
Instrumentos S⁴® Cervical
- (NL)** Gebruiksaanwijzing/Technische beschrijving
S⁴® Cervical Instruments

- (S)** Bruksanvisning/Teknisk beskrivning
S⁴® Cervical Instrument
- (RUS)** Инструкция по применению/Техническое описание
Инструменты S⁴® Cervical
- (CZ)** Návod k použití/Technický popis
Nástroje S⁴® Cervical
- (PL)** Instrukcja użytkowania/Opis techniczny
Instrumenty S⁴® Cervical
- (SK)** Technical alterations reserved
Nástroje S⁴® Cervical
- (TR)** Kullanım Kılavuzu/Teknik açıklama
S⁴® Cervical Aletler

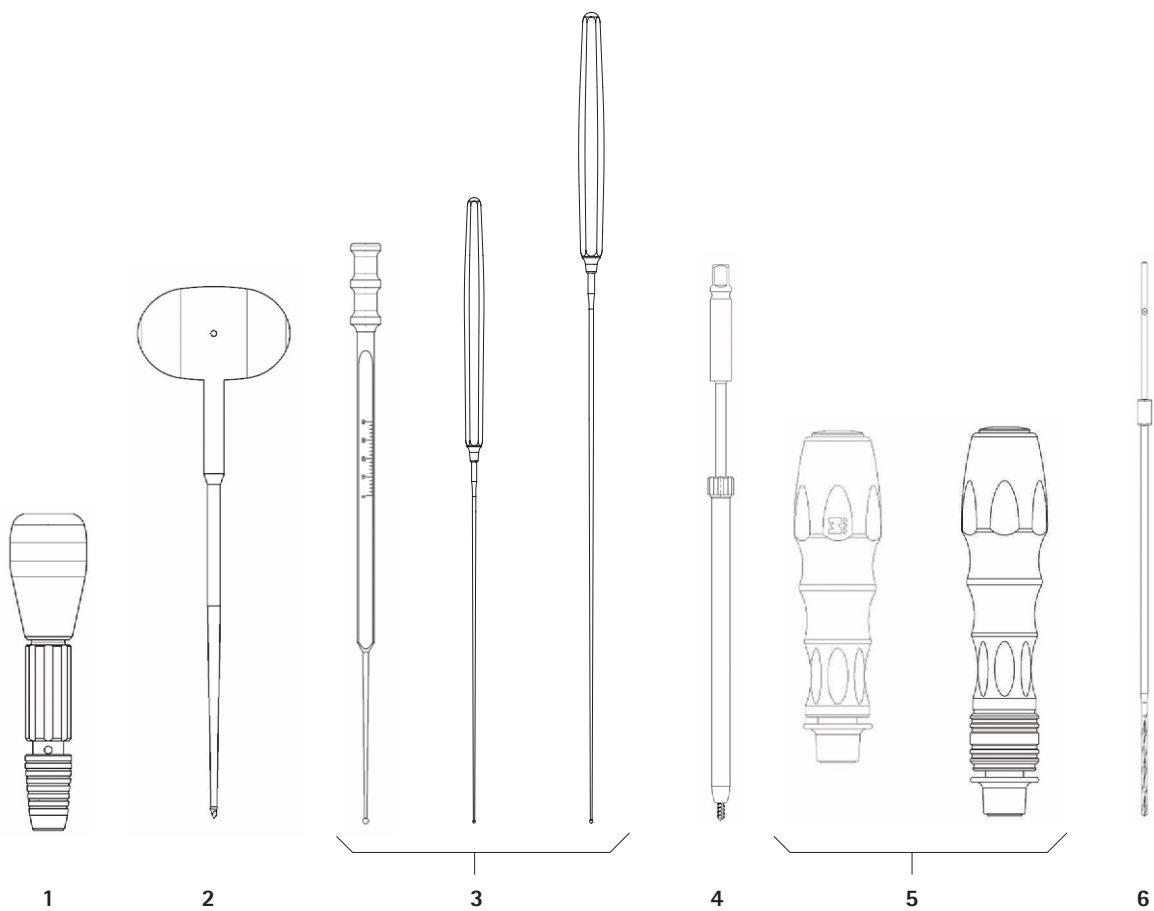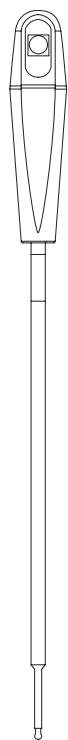
9


10

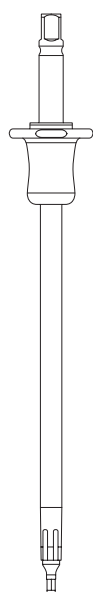
11

12

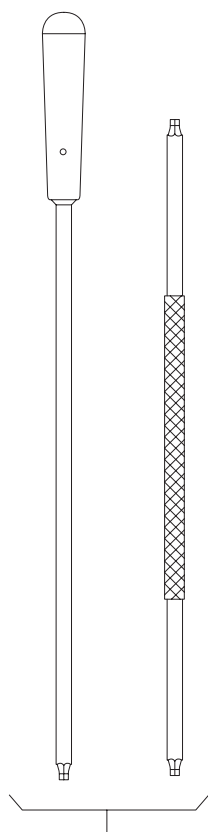
13

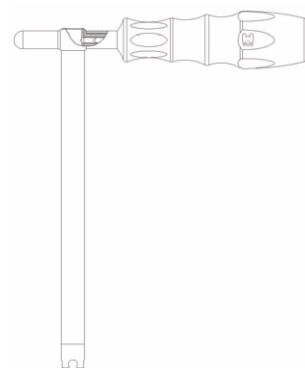
14


15

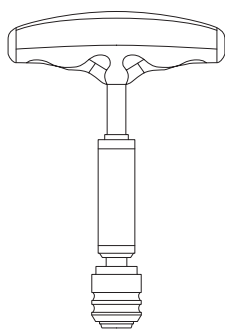16

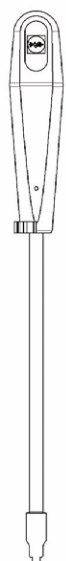
17

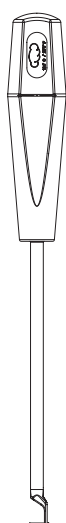
18


19

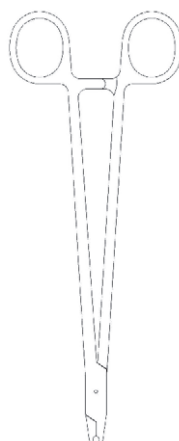
20

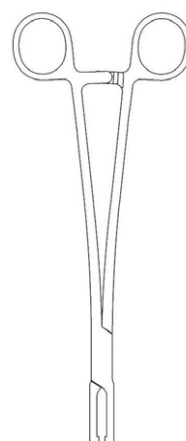
21

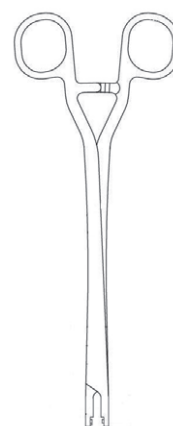
22

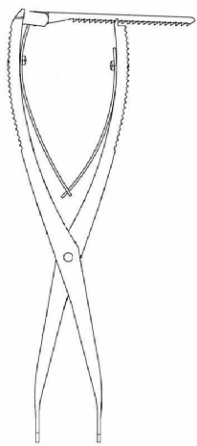
23

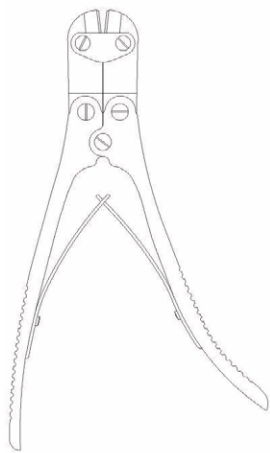
24

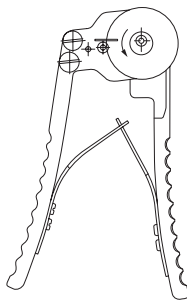
25

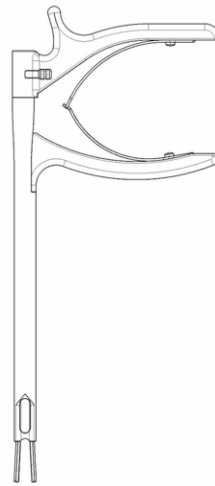
26

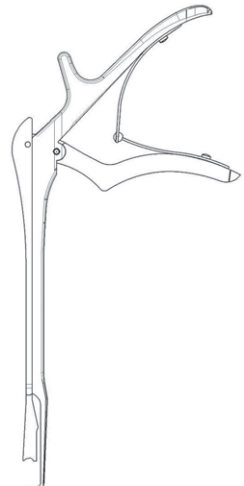27

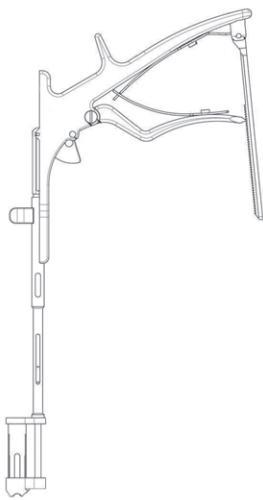
28

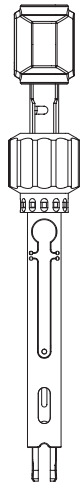
29

30

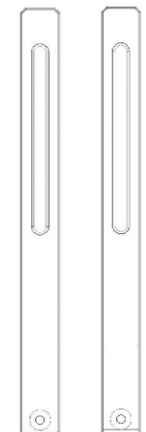
31

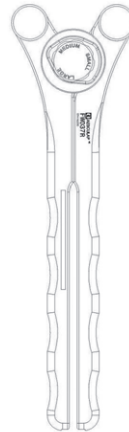
32

33

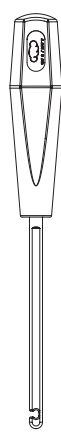34

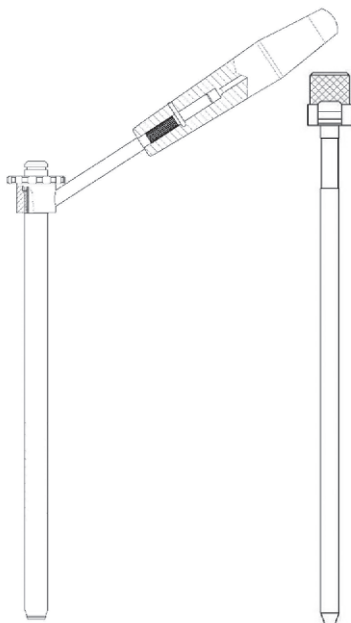35

36


37


38


39

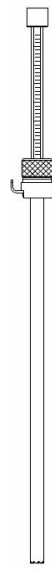40

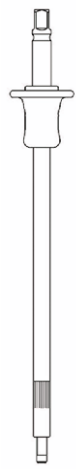41

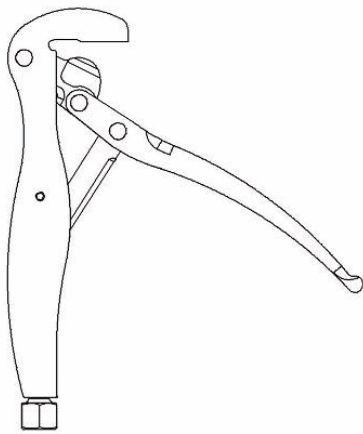42

43

44




45

46

47

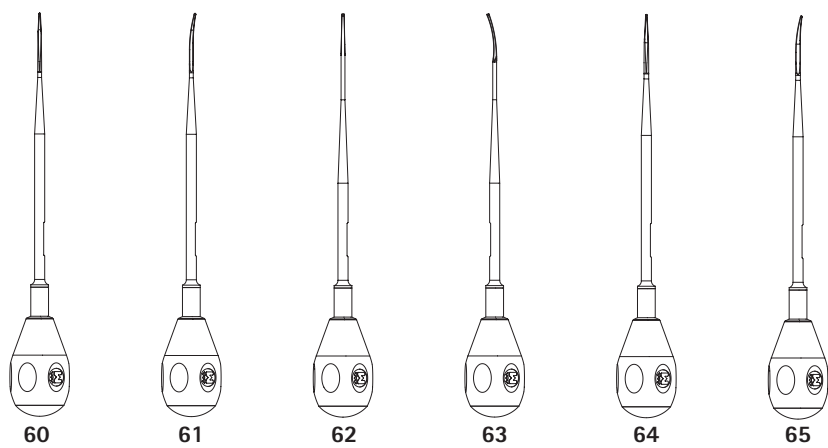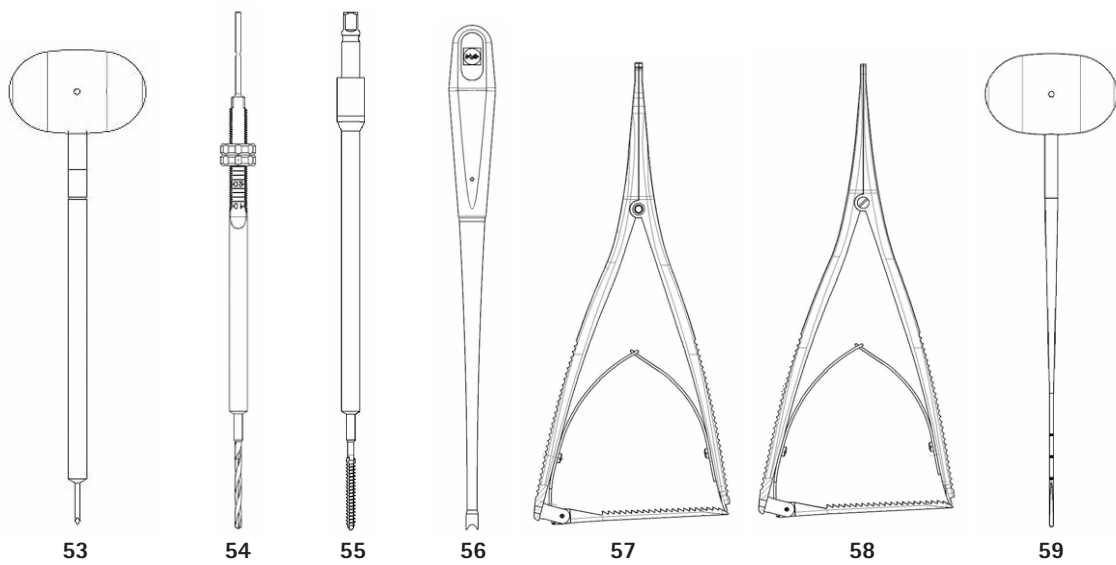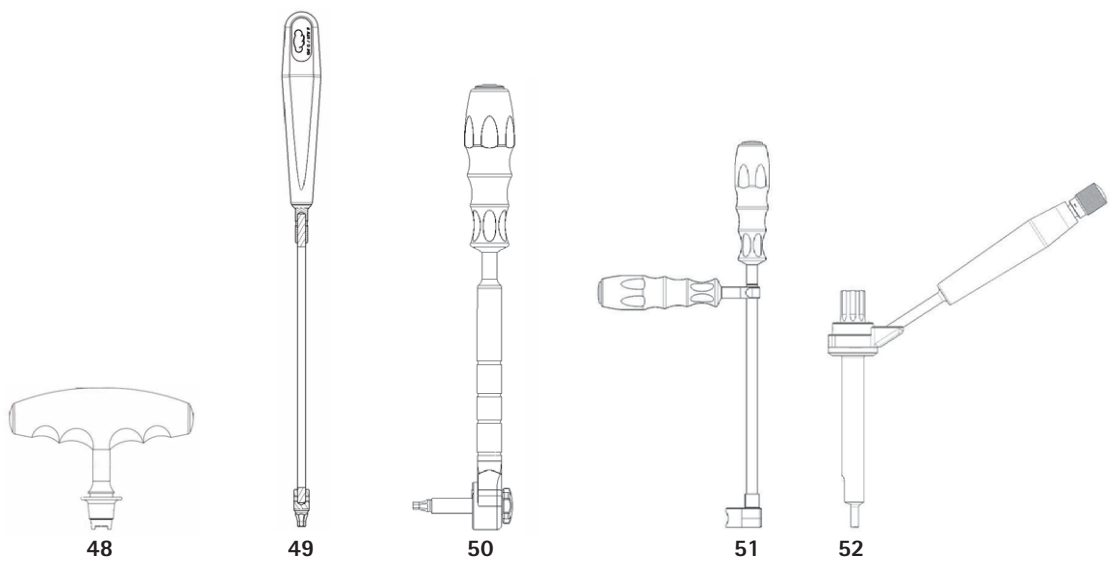Aesculap® S⁴® Cervical

S⁴® Cervical Instruments

Legend

- 1 Twist drill handle FJ839R
- 2 Bone awl FW041R
- 3 Depth gauge FW042R, sounder FW044R & long version sounder FW671R
- 4 Tap FW046R & FW047R
- 5 Ratchet handle FW165R (non ratchet handle FW067R)
- 6 Drill FW051SU for 3.5 mm screws (FW052SU for 4.0 mm screws)
- 7 Variable drill guide FW053R
- 8 Fixed drill guide FW049R (14 mm) & FW669R (12 mm)
- 9 Ball end screwdriver FJ968R (long version for favored angle screws FJ988R)
- 10 Set screw removal screwdriver FW064R
- 11 Self-holding polyaxial screwdriver FW070R and FW131R (long versions for favored angle screws FW069R and FW132R)
- 12 Threaded screwdriver FW128R
- 13 Set screw starters FW058R, FW059R, FW133R and FW134R
- 14 Counter torque FW062R
- 15 Torque indicating screwdriver FW061R
- 16 Torque limiting T-handle FW129R
- 17 Screw body manipulator FW065R
- 18 Polyaxial screw revision instrument FW135R
- 19 Lamina preparator FW071R
- 20 Rod holding forceps FW076R
- 21 Straight hook holding forceps FW422R
- 22 Curved hook holding forceps FW528R
- 23 Compression forceps FW427R
- 24 Rod cutter FW082R
- 25 Straight cut rod cutter FW672R
- 26 Rod persuader FW084R
- 27 Rod persuader FW077R
- 28 Rod persuader ME084R
- 29 Linear rod persuader FW673R
- 30 Rod templates FW078R, FW080R and FW081R
- 31 Rod bending plates FW036R
- 32 Rod bender FW037R
- 33 In-Situ rod benders FW073R and FW074R
- 34 In-Situ rod bender FW136R
- 35 Favored angle guide sleeve FW066R
- 36 Apfelbaum C1-C2 obturator FJ983R
- 37 Apfelbaum trocar FJ984R
- 38 Favored angle screw tap FW089R
- 39 Favored angle screw drill FW088SU
- 40 Apfelbaum inner sleeve guide FJ985R
- 41 Occipital screwdriver FW213R
- 42 Occipital plate bending pliers FW090R
- 43 Drill bit for occipital plate screws 4.5 mm FW091SU (5.5 mm FW092SU)
- 44 Tap for occipital plate screws 4.5 mm FW093R (5.5 mm FW094R)
- 45 Occipital drill guide double ended for 4.5 mm screws FW095R (for 5.5 mm screws FW096R), Occipital tap guide double ended for 4.5 mm screws FW097R (for 5.5 mm screws FW098R)
- 46 Occipital screw remover FW099R
- 47 Occipital screw removal shaft FW101R
- 48 Occipital T-handle for screw removal FW116R
- 49 Occipital flexible set screw starter FW109R
- 50 Occipital torque wrench FW103R
- 51 Occipital counter torque handle FW104R
- 52 Guide tube for smooth shank screws FW054R
- 53 Smooth shank bone awl FW085R
- 54 Smooth shank drill with depth gauge FW086SU
- 55 Smooth shank screw tap FW087R
- 56 Rod pusher FW063R
- 57 Distraction forceps, curved FW428R
- 58 Distraction forceps, straight FW523R
- 59 Pedicle probe FW045R
- 60 Bone probe straight FW674R
- 61 Bone probe curved FW675R
- 62 Bone probe blunt straight FW676R
- 63 Bone probe blunt curved FW677R
- 64 Bone probe thoracic straight FW678R
- 65 Bone probe thoracic curved FW679R

Symbols on product and packages

	Sterilization using irradiation
	Not for reuse in intended applications as defined by the manufacturer
	Use by
	Caution, general warning symbol Caution, see documentation supplied with the product
	Date of manufacture

Contents

1.	Applicable to	3
2.	Intended use	3
3.	Safe handling and preparation	3
4.	Safe operation	4
4.1	Preparing the holes for S ⁴ Cervical screws	4
4.2	Drilling the holes for S ⁴ Cervical screws	4
4.3	Tapping (optional)	4
4.4	Positioning and temporary fixation of an S ⁴ Cervical screw	4
4.5	Adapting the S ⁴ Cervical rods	4
4.6	Rod persuader (26, 27, 28, 29)	4
4.7	Insertion of set screw	4
4.8	Removing the S ⁴ Cervical implants	4
4.9	Pseudo rod revision instrument	4
5.	Disassembling	5
5.1	Twist drill handle (1)	5
5.2	S ⁴ Cervical tap (4)	5
5.3	S ⁴ Cervical self holding polyaxial screwdriver (11)	5
5.4	Rod persuader (26)	5
5.5	S ⁴ Cervical TH Threaded screwdriver (12)	5
5.6	Linear rod persuader (5)	5
6.	Assembling	6
6.1	Twist drill handle (1)	6
6.2	S ⁴ Cervical tap (4)	6
6.3	S ⁴ Cervical self holding polyaxial screwdriver (11)	6
6.4	Rod persuader (26)	6
6.5	S ⁴ Cervical TH Threaded screwdriver (12)	6
6.6	Linear rod persuader (5)	6
6.7	Straight cut rod cutter (25)	6
7.	Validated reprocessing procedure	7
7.1	General safety instructions	7
7.2	Single-use products (FW051SU, FW052SU, FW086SU, FW088SU, FW091SU, FW092SU)	7
7.3	General information	7
7.4	Disassembling the product before carrying out the reprocessing procedure	7
7.5	Preparations at the place of use	7
7.6	Preparation before cleaning	7
7.7	Cleaning/disinfection	7
7.8	Manual cleaning/disinfection	7
7.9	Mechanical cleaning/disinfecting	7
7.10	Mechanical cleaning/disinfection with manual pre-cleaning	7
7.11	Inspection, maintenance and checks	7
7.12	Packaging	7
7.13	Steam sterilization	7
7.14	Sterilization for the US market	7
7.15	Storage	7

8.	Technical Service	23
9.	Disposal	23
10.	Distributor in the US/Contact in Canada for product information and complaints	23

1. Applicable to

- For item-specific instructions for use and information on material compatibility, see also the Aesculap Extranet at www.extranet.bbraun.com

2. Intended use

The S⁴ Cervical system is used for posterior occipital cervical and thoracic stabilization and fusion. Indications and contraindications are listed in the instructions for use for the implants (TA011796). The instruments listed in the captions form part of this system. They are used for adapting the S⁴ Cervical implants to the individual patient and to position and insert the implants.

3. Safe handling and preparation

CAUTION

Federal law restricts this device to sale by, or on order of a physician!

- Always follow the instructions for use of the S⁴ Cervical Implants (TA011796) and the respective operating manual (O34202).
- Ensure that the product and its accessories are operated and used only by persons with the requisite training, knowledge, or experience.
- Read, follow, and keep the instructions for use.
- Use the product only in accordance with its intended use, see Intended use.
- Remove the transport packaging and clean the new product, either manually or mechanically, prior to its initial sterilization.
- Store any new or unused products in a dry, clean, and safe place.
- Inspect the product after each cleaning and disinfecting cycle to be sure that it is clean, functioning properly and not damaged.
- Prior to each use, inspect the product for loose, bent, broken, cracked, worn, or fractured components.
- Do not use the product if it is damaged or defective. Set aside the product if it is damaged.
- Replace any damaged components immediately with original spare parts.
- To avoid damage to the working end: Carefully insert the product through the working channel (e.g. trocar).
- Do not use products that are in open or damaged sterile packaging.
- Do not use the product after its use-by date.

Aesculap® S⁴® Cervical

S⁴® Cervical Instruments

Single-use products (FW051SU, FW052SU, FW086SU, FW088SU, FW091SU, FW092SU)

The product has been radiation-sterilized and is supplied in sterile packaging.

The product must not be reused.

- ▶ Ensure that the product and its accessories are operated and used only by persons with the requisite training, knowledge, or experience.
- ▶ Read, follow, and keep the instructions for use.
- ▶ Use the product only in accordance with its intended use, see Intended use.
- ▶ Do not use products from open or damaged sterile packaging.
- ▶ Prior to each use, inspect the product for loose, bent, broken, cracked or fractured components.
- ▶ Do not use the product if it is damaged or defective. Set aside the product if it is damaged.
- ▶ Do not use the product after its use-by date.

4. Safe operation

WARNING

Risk of injury and malfunction!

- ▶ Always carry out a function check prior to using the product.

4.1 Preparing the holes for S⁴ Cervical screws

The holes for self-tapping S⁴ Cervical screws are started using the S⁴ Cervical bone awl 2. The cortical layer of the vertebral body is penetrated in this process. A raised edge is provided on the standard bone awl 2 which provides an indication as to when the ideal depth has been reached.

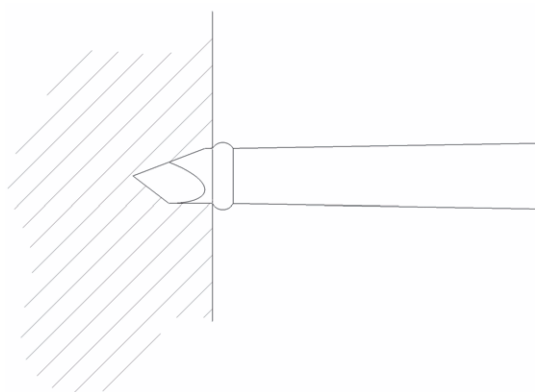
Fig. 1 Awl fully inserted

Opening the cortical layer with the bone awl 2

DANGER

Risk of damage to the spinal cord, nerve roots, adjacent intervertebral space or soft tissue when inserting the bone awl!

- ▶ Never insert the awl beyond the raised edge, see Fig 1.

4.2 Drilling the holes for S⁴ Cervical screws

DANGER

Risk of injury from an incorrectly placed hole or a hole that is too deep!

- ▶ Do not sharpen the drill, as this makes it the reading on the depth gauge inaccurate and falsifies the reading.
- ▶ Replace blunt drills with new ones.

The drill is applied with a drill guide and can be driven manually with twist drill handle 1 or with a motor system and Aesculap intra handpiece (e.g., GD450R/GD456R).

Mounting the drill and handle (only for manual drilling)

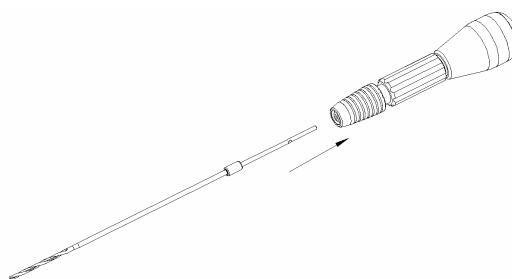
Fig. 2 Mounting the drill

- ▶ Use twist drill handle 1 for manual drilling.

DANGER

Risk of damage to the spinal cord, nerve roots, adjacent intervertebral space or soft tissue when drilling the holes!

- ▶ Always use the correct S⁴ Cervical drilling guides when drilling the holes. Do not combine drills and guides incorrectly.
- ▶ Always check the drill lengths with the guide sleeve being used using a caliper (e.g., AA845R, Caspar instrument for anterior cervical fusion) before starting to drill.

Injury to spinal cord and nerve roots caused by application of a drill that is too long!

- ▶ Use x-ray image to select an appropriate drill length prior to the operation.
- ▶ Alignment and insertion of the drills may only be executed under radiographic monitoring and/or aided by a navigation system.
- ▶ Select a drill of a length equivalent to the intended drill hole depth.

- ▶ To insert drill in twist drill handle 1:
 - Pull back locking sleeve against the spring pressure, in the direction indicated by the arrow and hold it there.
 - Push drill into the adapter of twist drill handle as far as it will go.
 - Slightly twist drill while releasing locking sleeve.

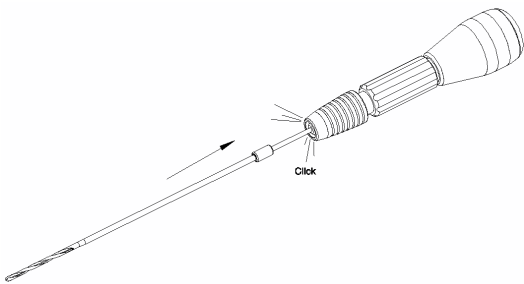
Fig. 3 Drill engages audibly

Risk of injury and damage to the drill if the rotational speed is too high!

- ▶ Be sure to set a rotational speed appropriate for the intended application.

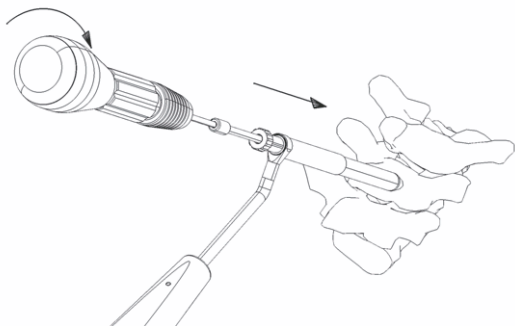
Fig. 4 Drilling with the drill guide (FW053R, FW049R or FW669R)

- ▶ For the variable drill guide 7, set the depth to be drilled on the depth stop by rotating the plastic collar on the guide sleeve. The fixed drill guide (14 mm) 8 has a green handle. The fixed drill guide (12 mm) 8 has a red handle.
- ▶ Check the length to be drilled prior to drilling by inserting the drill in the guide and measuring the cutting depth with a caliper (e.g., AA845R).
- ▶ Drill to the preset depth. The stop unit stops the drill on the guide sleeve.

4.3 Tapping (optional)

S⁴ Cervical screws are self-tapping screws. However, if the bone quality is found to be hard, the operating surgeon can decide to pre-tap the thread with the S⁴ Cervical tap.

Risk of tissue injury when inserting the S⁴ Cervical tap or stripping bone threads!

- ▶ Always check if the sleeve retracts correctly prior to applying the S⁴ Cervical tap.

- ▶ To tap the thread:
 - Slowly drive the S⁴ Cervical tap 4 in steadily until the required depth is reached.
 - The depth is read from the scale as the sleeve retracts during tapping.

Fig. 5 Depth is read from scale

Aesculap® S⁴® Cervical

S⁴® Cervical Instruments

4.4 Positioning and temporary fixation of an S⁴ Cervical screw

Picking up the S⁴ Cervical screw

- Pick up the S⁴ Cervical screw with the S⁴ Cervical self holding polyaxial screwdriver 11.

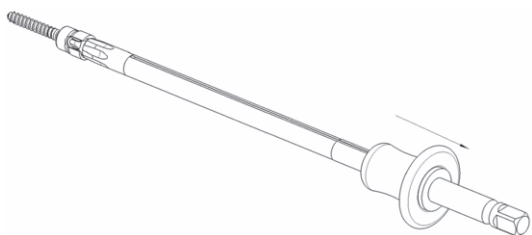
Fig. 6 Picking up the S⁴ Cervical screw

- To prevent the S⁴ Cervical screw from dropping out of the S⁴ Cervical self holding polyaxial screwdriver 11 during transfer to the operating surgeon, lock the instrument by pulling back and releasing the blue knob.
- Make sure that the S⁴ Cervical screw is firmly seated on the S⁴ Cervical self holding polyaxial screwdriver and the polyaxiality of the screw is blocked for insertion.

Using the S⁴ Cervical TH threaded screwdriver 12

- Hold the S⁴ Cervical screw body in one hand.
- Insert the S⁴ Cervical TH threaded screwdriver 12 into the S⁴ Cervical screw hex.
- To affix the S⁴ Cervical TH threaded screwdriver 12 to the screw, rotate the blue handle clockwise until the screwdriver thread engages the screw thread.
- Tighten until resistance is felt. Do not overtighten as this can cause difficulty disengaging the screwdriver after screw placement.
- After insertion of the S⁴ Cervical screw into the patient: Rotate the blue handle counter-clockwise to release the S⁴ Cervical TH threaded screwdriver 12.

4.5 Adapting the S⁴ Cervical rods

The Aesculap S⁴ Cervical rods are supplied straight (excluding pre-bent for occipital plate). Their curvature can be adjusted to the situation or to the required spine lordosis using the S⁴ Cervical bending instruments.

- Use the rod templates 30 for examining the degree of rod bend required.

Note

For long S⁴ Cervical rods, bending should be carried out in several steps in order to avoid excessive or insufficient lordosis.

DANGER

Damage to and breakage of the S⁴ Cervical rods caused by excessive strain on the material!

- Always bend the S⁴ Cervical rods in one direction only. Do not bend back the S⁴ Cervical rods.
- Always use S⁴ Cervical bending instruments (e.g., 31, 32, 33, 34) for bending S⁴ Cervical rods.
- Avoid small bending radii, bending back, notching or scratching of S⁴ Cervical rods.

- To bend the S⁴ Cervical rods close to their ends, use the rod bending plates 31.

DANGER

Risk of damage to patient, personnel or equipment in the operating room due to flying fragment from cut rod!

- Always hold the piece of rod to be cut at both ends.

- To cut the S⁴ Cervical rods use the rod cutter 24 or 25.
- To cut the S⁴ Cervical rods with rod cutter 25:
 - Ensure that the handles of the rod cutter are in the free, relaxed and open position.
 - Rotate the large knob clockwise until both lines match between the body of the device and the rotating knob.
 - Mark the cutting location on the rod.
 - Slide the rod into the center hole of the rod cutter through the rotating knob until the mark on the rod is approximately at the location of the cut line. The rod will be cut at the side of the rod cutter marked with "cut location".
 - Squeeze handles together until they come close together.
 - Release pressure on the handles and allow them to fully open again to the relaxed position.
 - Repeat several times until the rod is fully cut and can easily be removed from the rod cutter.

If desired, sharp edges left after cutting the rod can be filed with a suitable instrument such as the Aesculap durogrip file OL429R.

- Hold the rod during insertion using the S⁴ Cervical rod holding forceps **20**.

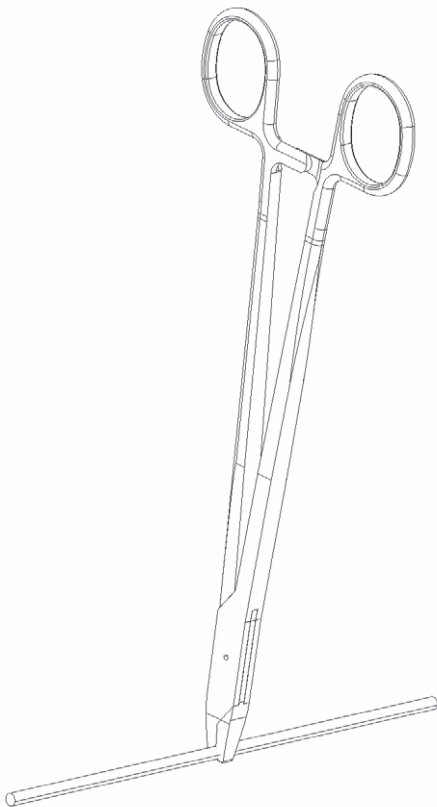
Fig. 7 Holding the rod

- Likewise hold the S⁴ Cervical hooks using the S⁴ Cervical hook holding forceps **23, 24**.

4.6 Rod persuader (26, 27, 28, 29)

Damage to and breakage of the rod persuader!

- Never use the rod persuader to bend the S⁴ Cervical rod.
- Never use the rod persuader to act as a counter torque for the S⁴ Cervical set screws.

- To persuade with rod persuader **28** or **29**:
 - Make sure that the persuader is placed axial to the screw head.
 - Slide the persuader blades over the screw head so that the persuader reaches its operational point.
 - Make sure that the pins of the persuader blades interface with the holes in the head of the polyaxial screw.
 - Make sure that both pins of the persuader blades engage the holes of the polyaxial screw.
 - Make sure to hold the persuasion knob as the passage of the set screw may be obstructed by deflection of the hinge when holding the linear rod persuader **29** by the outer sleeve.

4.7 Insertion of set screw

Damage to the S⁴ Cervical screw caused by incorrectly applying the torque indicating screwdriver 15 when driving in the set screws!

This screwdriver is an indicating instrument. Align the two opposing arrows to apply the proper torque.

- Always use the counter torque instrument **14** when tightening the S⁴ Cervical set screw in the S⁴ Cervical screw body.
- Refer to TA011947 for further information.

Damage to the S⁴ Cervical implants caused by incorrectly applying the occipital flexible set screw starter 49 when driving in the set screws!

- Refer to TA012676 S⁴ Cervical Occiput Mini Set Screw Starter Flexible Shaft.

Aesculap® S⁴® Cervical

S⁴® Cervical Instruments

4.8 Removing the S⁴ Cervical implants

CAUTION

Damage to the S⁴ Cervical Torque Indicating Screwdriver 15 or 16!

- ▶ Never use the S⁴ Cervical Torque Indicating Screwdriver 15 or 16 for extracting screws.
- ▶ Use the S⁴ Cervical set screw removal screwdriver 10 together with the handle 5 and counter torque 14 when extracting set screws.
- ▶ Refer to TA011948 for further information.

CAUTION

Damage to the S⁴ Cervical implants caused by incorrectly applying the occipital torque wrench 50 when driving in the set screws!

- ▶ When tightening the set screw always apply counter torque using the countering instrument for this purpose.
- ▶ Refer to TA012367 S⁴ Cervical Spinal System - Torque wrenches with click function.

4.9 Pseudo rod revision instrument

CAUTION

Damage to and breakage of the pseudo rod revision instrument!

- ▶ Never use the pseudo rod revision instrument to release a locked S⁴ Cervical polyaxial screw tulip.

5. Disassembling

5.1 Twist drill handle (1)

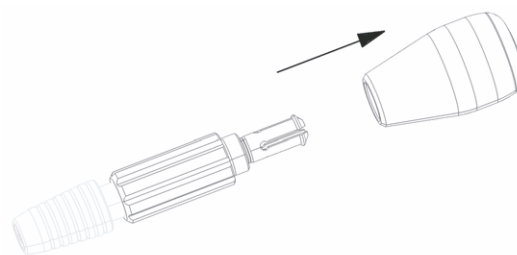
Fig. 8 Disassembling of twist drill handle

- ▶ Pull off twist drill handle 1 from collar in the direction indicated by the arrow.

5.2 S⁴ Cervical tap (4)

Fig. 9 Disassembling of taps

- ▶ Loosen and unscrew nut and remove it from the sleeve in the direction of the arrow.
- ▶ Remove sleeve together with the spring in the direction of arrow.

5.3 S⁴ Cervical self holding polyaxial screwdriver (11)

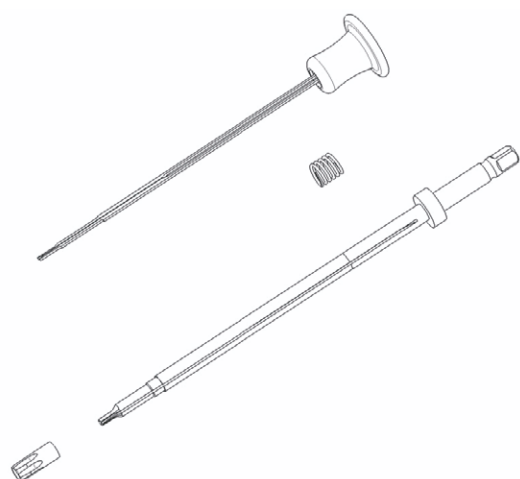
Fig. 10 Exploded view of S⁴ Cervical self holding polyaxial screwdriver

Note

These instructions are also valid for FW070R, FW131R, FW069R, FW132R.

CAUTION

The screwdriver will not function properly if the retaining tongues are bent or crimped!

► **Do not bend or crimp the retaining tongues.**

- Pull back on blue knob (black for FW069R and FW132R).
- Loosen and unscrew the screw sleeve at the working end of the instrument.
- Keep the blue knob (black for FW069R and FW132R) pulled back until the screw sleeve is completely unscrewed.
- Slide the knob and the retaining tongues towards the working end of the instrument and off the screwdriver shaft.
- Remove the spring.

5.4 Rod persuader (26)

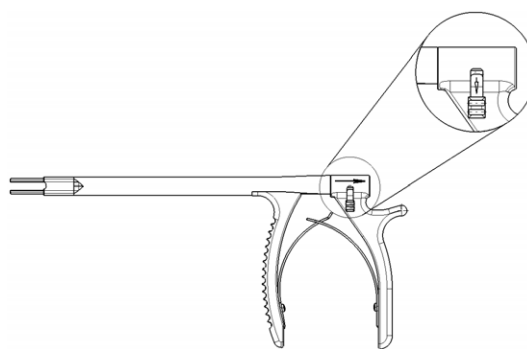
Fig. 11 Disassembling of rod persuader

- Pull the release mechanism on the back handle to the downward position and maintain this position until the rear handle has been removed.
- Holding the outer shaft pull the front handle towards the rear of the persuader until fully disassembled.

5.5 S⁴ Cervical TH Threaded screwdriver (12)

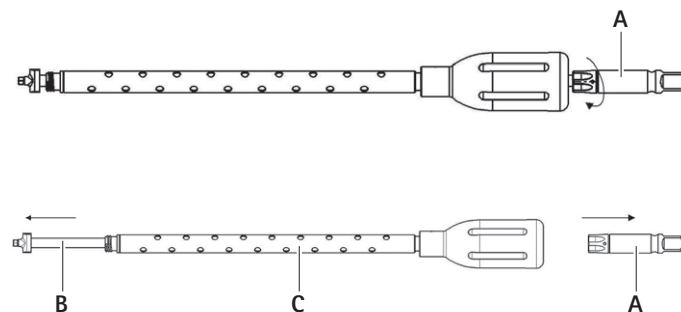
Fig. 12 Disassembling the threaded screwdriver

- Hold the screwdriver 12 at the distal or hex end.
- Turn the adapter for handle A on the back end of instrument counter-clockwise.
- Remove the screwdriver shaft B through the handle-sleeve assembly C.

Aesculap® S⁴® Cervical

S⁴® Cervical Instruments

5.6 Linear rod persuader (5)

- Ensure the proximal knob is in the unlock extended position.
This will retract the locking sleeves from the inner sleeve.

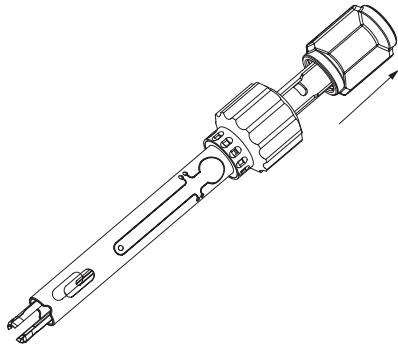
Fig. 13 Proximal knob is in the unlock extended position

- Turn the persuasion nut as if persuading the rod until the screw threads are no longer engaged.

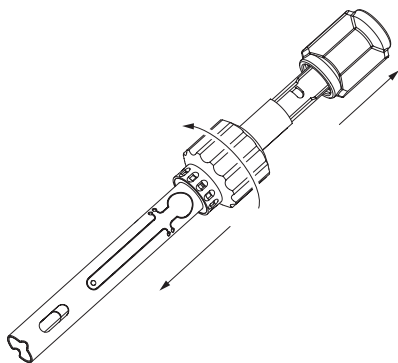
Fig. 14 Turning the persuasion nut

- Press the circular tabs of the hinge.
This will lift the pins out of the inner sleeve, enabling it to be pulled out of the outer sleeve.

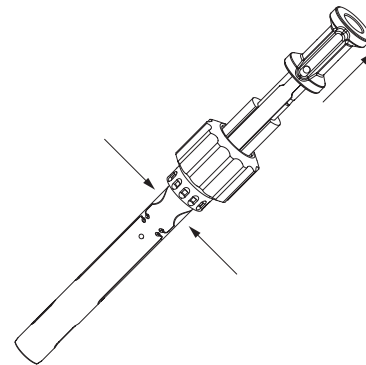
Fig. 15 Pulling out the inner sleeve

- Remove the locking sleeves by rotating them out from under the knob.

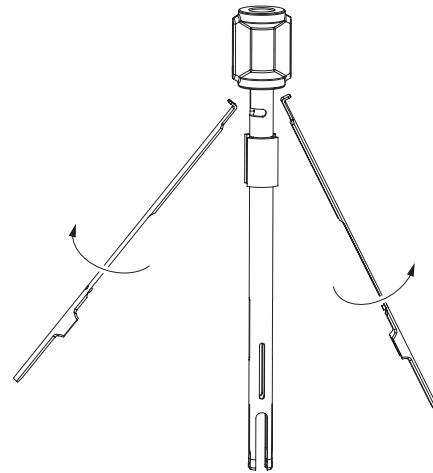
Fig. 16 Removing the locking sleeves

6. Assembling

6.1 Twist drill handle (1)

- Push handle onto collar against the direction indicated by the arrow, see Fig. 8.

6.2 S⁴ Cervical tap (4)

- Slide guide sleeve with spring over the tap against the direction of the arrow and fasten it by turning the nut clockwise, see Fig. 9.

6.3 S⁴ Cervical self holding polyaxial screwdriver (11)

Note

These instructions are also valid for FW070R, FW131R, FW069R, FW132R.

CAUTION

The screwdriver will not function properly if the retaining tongues are bent or crimped!

- Do not bend or crimp the retaining tongues.

- Replace the spring on the shaft.
- Slide the knob with retaining tongues onto the screwdriver shaft in such a way that the retaining tongues engage the grooves of the screwdriver shaft, see Fig. 10.
- Pull back on blue knob (black for FW069R and FW132R).
- Screw the screw sleeve onto the tip of the screwdriver and tighten it down.
- Keep the blue knob (black for FW069R and FW132R) pulled back until the screw sleeve is completely tightened.

6.4 Rod persuader (26)

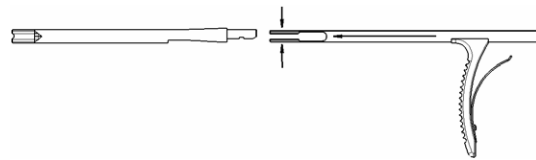
Fig. 17 Assembling of rod persuader

- Align the gripping end with the outer shaft and carefully depress the gripping fingers inward until secured in the outer shaft.
- Position the rear handle concentrically with the outer shaft and apply pressure until an audible click is heard confirming locking.
- Squeeze the handle to verify successful assembly.

6.5 S⁴ Cervical TH Threaded screwdriver (12)

- Slide screwdriver shaft B through the handle-sleeve assembly C, see Fig. 12.
- Insert the adapter for handle A onto screwdriver shaft B.
- Turn the adapter for handle A clockwise and fully tighten thread.
- Make certain the adapter for handle A is threaded onto screwdriver shaft B.

Aesculap® S⁴® Cervical

S⁴® Cervical Instruments

6.6 Linear rod persuader (5)

- Engage the locking sleeves into the proximal nut by rotating them into the groove of the knob.

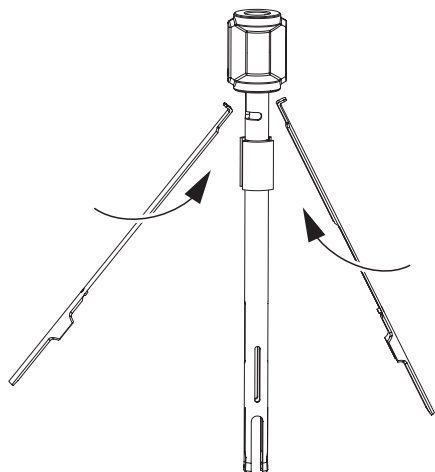
Fig. 18 Engaging the locking sleeves

- Slide the inner assembly into the outer sleeve.

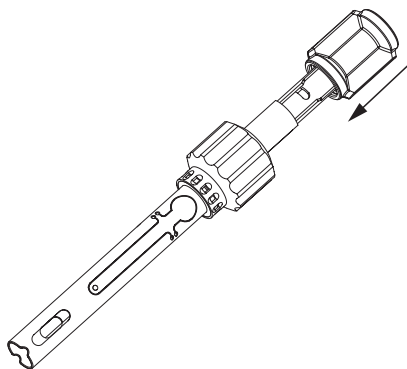
Fig. 19 Sliding the inner assembly into the outer sleeve

- Press the circular tabs of the hinge to lift the pins and slide the inner tube into position.

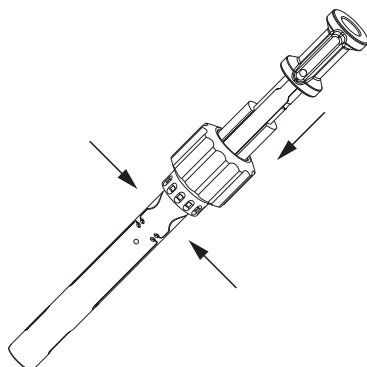
Fig. 20 Sliding the inner tube into the position

- Re-engage the persuasion nut and thread completely to the starting position to prepare for use.

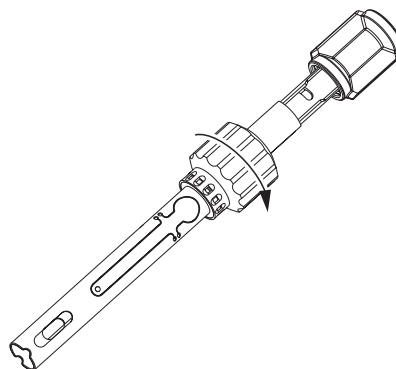
Fig. 21 Re-engaging the persuasion nut

6.7 Straight cut rod cutter (25)

- Apply Aesculap STERILIT® lubricant as needed.

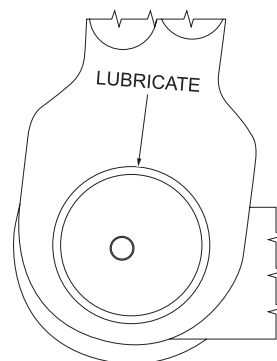
Fig. 22 Applying lubricant

7. Validated reprocessing procedure

7.1 General safety instructions

Note

Adhere to national statutory regulations, national and international standards and directives, and local, clinical hygiene instructions for sterile processing.

Note

For patients with Creutzfeldt-Jakob disease (CJD), suspected CJD or possible variants of CJD, observe the relevant national regulations concerning the reprocessing of products.

Note

Mechanical reprocessing should be favored over manual cleaning as it gives better and more reliable results.

Note

Successful processing of this medical device can only be ensured if the processing method is first validated. The operator/sterile processing technician is responsible for this.

The recommended chemistry was used for validation.

Note

If there is no final sterilization, then a virucidal disinfectant must be used.

Note

For the latest information on reprocessing and material compatibility see also the Aesculap extranet at www.extranet.bbraun.com

The validated steam sterilization procedure was carried out in the Aesculap sterile container system.

7.2 Single-use products (FW051SU, FW052SU, FW086SU, FW088SU, FW091SU, FW092SU)

WARNING

Risk of infection for patients and/or users and impairment of product functionality due to reuse.
Risk of injury, illness or death due to contamination and/or impaired functionality of the product!

► Do not reprocess the product!

7.3 General information

Dried or affixed surgical residues can make cleaning more difficult or ineffective and lead to corrosion. Therefore the time interval between application and processing should not exceed 6 h; also, neither fixating pre-cleaning temperatures >45 °C nor fixating disinfecting agents (active ingredient: aldehydes/alcohols) should be used.

Excessive measures of neutralizing agents or basic cleaners may result in a chemical attack and/or to fading and the laser marking becoming unreadable visually or by machine for stainless steel.

Residues containing chlorine or chlorides e.g. in surgical residues, medicines, saline solutions and in the service water used for cleaning, disinfection and sterilization will cause corrosion damage (pitting, stress corrosion) and result in the destruction of stainless steel products. These must be removed by rinsing thoroughly with demineralized water and then drying.

Additional drying, if necessary.

Only process chemicals that have been tested and approved (e.g. VAH or FDA approval or CE mark) and which are compatible with the product's materials according to the chemical manufacturers' recommendations may be used for processing the product. All the chemical manufacturer's application specifications must be strictly observed. Failure to do so can result in the following problems:

- Optical changes of materials, e.g. fading or discoloration of titanium or aluminum. For aluminum, the application/process solution only needs to be of pH >8 to cause visible surface changes.
- Material damage such as corrosion, cracks, fracturing, premature aging or swelling.
- Do not use metal cleaning brushes or other abrasives that would damage the product surfaces and could cause corrosion.
- Further detailed advice on hygienically safe and material-/value-preserving reprocessing can be found at www.a-k-i.org, link to Publications, Red Brochure – Proper maintenance of instruments.

Aesculap® S⁴® Cervical

S⁴® Cervical Instruments

7.4 Disassembling the product before carrying out the reprocessing procedure

- ▶ Disassemble the product immediately after use, see Disassembling and TA012367.
- ▶ Open up instruments with hinges.

7.5 Preparations at the place of use

- ▶ If applicable, rinse non-visible surfaces preferably with deionized water, with a disposable syringe for example.
- ▶ Remove any visible surgical residues to the extent possible with a damp, lint-free cloth.
- ▶ Transport the dry product in a sealed waste container for cleaning and disinfection within 6 hours.

7.6 Preparation before cleaning

- ▶ Disassemble the product prior to cleaning, see Disassembling.

7.7 Cleaning/disinfection

Product-specific safety notes on the reprocessing procedure

CAUTION

Damage to the product due to inappropriate cleaning/disinfecting agents and excessive temperatures!

- ▶ Use cleaning and disinfecting agents according to the manufacturer's instructions which
 - are approved for e.g. aluminum, plastics, high-grade steel,
 - do not attack softeners (e.g. in silicone).
 - ▶ Observe specifications regarding concentration, temperature and exposure time.
 - ▶ Do not exceed the maximum allowable washing temperature of 60 °C.
- ▶ Mount jaws protection on the product.
 - ▶ Carry out ultrasound cleaning:
 - as an effective mechanical supplement to manual cleaning/disinfection.
 - as a pre-cleaning procedure for products with encrusted residues, in preparation for mechanical cleaning/disinfection.
 - as an integrated mechanical support measure for mechanical cleaning/disinfection.
 - for additional cleaning of products with residues left after mechanical cleaning/disinfection.
 - ▶ Clean and disinfect microsurgical products mechanically if they can be placed securely in the machine or on the positioning aids.

Validated cleaning and disinfection procedure

Validated procedure	Specific requirements	Reference
Manual cleaning with immersion disinfection	■ For FJ968R, FJ983R, FJ984R, FJ988R, FW036R, FW037R, FW041R, FW042R, FW044R, FW045R, FW046R, FW047R, FW049R, FW053R, FW054R, FW058R, FW059R, FW062R, FW063R, FW064R, FW065R, FW066R, FW069R, FW070R, FW071R, FW073R, FW074R, FW076R, FW077R, FW078R, FW080R, FW081R, FW082R, FW084R, FW085R, FW087R, FW089R, FW090R, FW093R, FW094R, FW095R, FW096R, FW097R, FW098R, FW099R, FW101R, FW104R, FW109R, FW131R, FW132R, FW133R, FW134R, FW135R, FW136R, FW213R, FW422R, FW427R, FW428R, FW523R, FW528R, FW669R, FW671R, FW674R, FW675R, FW676R, FW677R, FW678R, FW679R, ME084R ■ Cleaning brush: e.g., TA011327, TA006874 for FW095R, FW096R, FW097R, FW098R ■ Disposable syringe 20 ml ■ Keep working ends open for cleaning. ■ When cleaning products with movable hinges, ensure that these are in an open position and, if applicable, move the joint while cleaning. ■ Drying phase: Use a lint-free cloth or medical compressed air	Chapter Manual cleaning/ disinfection and sub-chapter: ■ Chapter Manual cleaning with immersion disinfection
Manual cleaning with ultrasound and immersion disinfection	■ For FW061R, FW067R, FW103R, FW116R, FW165R, FW672R, FW673R ■ Cleaning brush: e.g. TA011944 ■ Disposable syringe 20 ml ■ Keep working ends open for cleaning. ■ When cleaning products with movable hinges, ensure that these are in an open position and, if applicable, move the joint while cleaning. ■ Drying phase: Use a lint-free cloth or medical compressed air	Chapter Manual cleaning/ disinfection and sub-chapter: ■ Chapter Manual cleaning with ultrasound and immersion disinfection
Mechanical alkaline cleaning and thermal disinfection	■ For FJ968R, FJ983R, FJ984R, FJ988R, FW036R, FW037R, FW041R, FW042R, FW044R, FW045R, FW046R, FW047R, FW049R, FW053R, FW054R, FW058R, FW059R, FW062R, FW063R, FW064R, FW065R, FW066R, FW069R, FW070R, FW071R, FW073R, FW074R, FW076R, FW077R, FW078R, FW080R, FW081R, FW082R, FW084R, FW085R, FW087R, FW089R, FW090R, FW093R, FW094R, FW095R, FW096R, FW097R, FW098R, FW099R, FW101R, FW104R, FW109R, FW131R, FW132R, FW133R, FW134R, FW135R, FW136R, FW213R, FW422R, FW427R, FW428R, FW523R, FW528R, FW669R, FW671R, FW674R, FW675R, FW676R, FW677R, FW678R, FW679R, ME084R ■ Place the product in a tray that is suitable for cleaning (avoiding rinsing blind spots). ■ Connect components with lumens and channels directly to the rinsing port of the injector carriage. ■ Keep working ends open for cleaning. ■ Place products in the tray with their hinges open.	Chapter Mechanical cleaning/disinfecting and sub-chapter: ■ Chapter Mechanical alkaline cleaning and thermal disinfecting

Aesculap® S⁴® Cervical

S⁴® Cervical Instruments

Validated procedure	Specific requirements	Reference
Manual pre-cleaning with brush and subsequent mechanical alkaline cleaning and thermal disinfection	■ For FJ839R, FW061R, FW067R, FW069R, FW070R, FW116R, FW103R, FW131R, FW132R, FW165R, FW669R, FW672R, FW673R ■ Cleaning brush: e.g. TA011327, TA006874 for FW095R, FW096R, FW097R, FW098R ■ Disposable syringe 20 ml ■ Place the product in a tray that is suitable for cleaning (avoiding rinsing blind spots). ■ Connect components with lumens and channels directly to the rinsing port of the injector carriage. ■ Keep working ends open for cleaning. ■ Place products in the tray with their hinges open.	Chapter Mechanical cleaning/disinfection with manual pre-cleaning and sub-chapter: ■ Chapter Manual pre-cleaning with a brush ■ Chapter Mechanical alkaline cleaning and thermal disinfecting
Manual pre-cleaning with ultrasound and brush, and subsequent mechanical alkaline cleaning and thermal disinfection	■ For FJ985R, FW128R ■ Cleaning brush: e.g. TA011944 ■ Disposable syringe 20 ml ■ Mount jaws protection on the product. ■ Place the product in a tray that is suitable for cleaning (avoiding rinsing blind spots). ■ Connect components with lumens and channels directly to the rinsing port of the injector carriage. ■ Keep working ends open for cleaning. ■ Place products in the tray with their hinges open.	Chapter Mechanical cleaning/disinfection with manual pre-cleaning and sub-chapter: ■ Chapter Manual pre-cleaning with ultrasound and brush ■ Chapter Mechanical alkaline cleaning and thermal disinfecting

7.8 Manual cleaning/disinfection

- Prior to manual disinfecting, allow water to drip off for a sufficient length of time to prevent dilution of the disinfecting solution.
- After manual cleaning/disinfection, check visible surfaces visually for residues.
- Repeat the cleaning /disinfection process if necessary.

Manual cleaning with immersion disinfection

Phase	Step	T [°C/°F]	t [min]	Conc. [%]	Water quality	Chemical
I	Disinfecting cleaning	RT (cold)	>15	2	D-W	Aldehyde-free, phenol-free, and QUAT-free concentrate, pH ~ 9*
II	Intermediate rinse	RT (cold)	1	-	D-W	-
III	Disinfection	RT (cold)	15	2	D-W	Aldehyde-free, phenol-free, and QUAT-free concentrate, pH ~ 9*
IV	Final rinse	RT (cold)	1	-	FD-W	-
V	Drying	RT	-	-	-	-

D-W: Drinking water

FD-W: Fully desalinated water (demineralized, low microbiological contamination: drinking water quality at least)

RT: Room temperature

*Recommended: BBraun Stabimed

- Note the information on appropriate cleaning brushes and disposable syringes, see Validated cleaning and disinfection procedure.

Phase I

- Fully immerse the product in the cleaning/disinfectant for at least 15 min. Ensure that all accessible surfaces are moistened.
- Clean the product with a suitable cleaning brush in the solution until all discernible residues have been removed from the surface.
- If applicable, brush through non-visible surfaces with an appropriate cleaning brush for at least 1 min.
- Mobilize non-rigid components, such as set screws, links, etc. during cleaning.
- Thoroughly rinse through these components with the cleaning disinfectant solution (at least five times), using a disposable syringe.

Phase II

- Rinse/flush the product thoroughly (all accessible surfaces) under running water.
- Mobilize non-rigid components, such as set screws, joints, etc. during rinsing.
- Drain any remaining water fully.

Phase III

- Fully immerse the product in the disinfectant solution.
- Mobilize non-rigid components, such as set screws, joints, etc. during rinsing.
- Rinse lumens at least 5 times at the beginning of the exposure time using an appropriate disposable syringe. Ensure that all accessible surfaces are moistened.

Phase IV

- Rinse/flush the product thoroughly (all accessible surfaces).
- Mobilize non-rigid components, such as set screws, joints, etc. during final rinse.
- Rinse lumens with an appropriate disposable syringe at least five times.
- Drain any remaining water fully.

Phase V

- Dry the product in the drying phase with suitable equipment (e.g. cloth, compressed air), see Validated cleaning and disinfection procedure.

Aesculap® S⁴® Cervical

S⁴® Cervical Instruments

Manual cleaning with ultrasound and immersion disinfection

Phase	Step	T [°C/°F]	t [min]	Conc. [%]	Water quality	Chemical
I	Ultrasonic cleaning	RT (cold)	>15	2	D-W	Aldehyde-free, phenol-free, and QUAT-free concentrate, pH ~ 9*
II	Intermediate rinse	RT (cold)	1	-	D-W	-
III	Disinfection	RT (cold)	15	2	D-W	Aldehyde-free, phenol-free, and QUAT-free concentrate, pH ~ 9*
IV	Final rinse	RT (cold)	1	-	FD-W	-
V	Drying	RT	-	-	-	-

D-W: Drinking water

FD-W: Fully desalinated water (demineralized, low microbiological contamination: drinking water quality at least)

RT: Room temperature

*Recommended: BBraun Stabimed

- Note the information on appropriate cleaning brushes and disposable syringes, see Validated cleaning and disinfection procedure.

Phase I

- Clean the product in an ultrasonic cleaning bath (frequency 35 kHz) for at least 15 min. Ensure that all accessible surfaces are immersed and acoustic shadows are avoided.
- Clean the product with a suitable cleaning brush in the solution until all discernible residues have been removed from the surface.
- If applicable, brush through non-visible surfaces with an appropriate cleaning brush for at least 1 min.
- Mobilize non-rigid components, such as set screws, links, etc. during cleaning.
- Thoroughly rinse through these components with the cleaning disinfectant solution (at least five times), using a disposable syringe.

Phase II

- Rinse/flush the product thoroughly (all accessible surfaces) under running water.
- Mobilize non-rigid components, such as set screws, joints, etc. during rinsing.
- Drain any remaining water fully.

Phase III

- Fully immerse the product in the disinfectant solution.
- Mobilize non-rigid components, such as set screws, joints, etc. during rinsing.
- Rinse lumens at least five times at the beginning of the exposure time with an appropriate disposable syringe. Ensure that all accessible surfaces are moistened.

Phase IV

- Rinse/flush the product thoroughly (all accessible surfaces) under running water.
- Mobilize non-rigid components, such as set screws, joints, etc. during final rinse.
- Rinse lumens with an appropriate disposable syringe at least five times.
- Drain any remaining water fully.

Phase V

- Dry the product in the drying phase with suitable equipment (e.g. cloth, compressed air), see Validated cleaning and disinfection procedure.

7.9 Mechanical cleaning/disinfecting

Note

The cleaning and disinfection device must be of tested and approved effectiveness (e.g. FDA approval or CE mark according to DIN EN ISO 15883).

Note

The cleaning and disinfection device used for processing must be serviced and checked at regular intervals.

Mechanical alkaline cleaning and thermal disinfecting

Machine type: single-chamber cleaning/disinfection device without ultrasound

Phase	Step	T [°C/°F]	t [min]	Water quality	Chemical/Note
I	Prerinse	<25/77	3	D-W	-
II	Cleaning	55/131	10	FD-W	■ Concentrate, alkaline: - pH = 13 - <5 % anionic surfactant ■ 0.5 % working solution - pH = 11*
III	Intermediate rinse	>10/50	1	FD-W	-
IV	Thermal disinfecting	90/194	5	FD-W	-
V	Drying	-	-	-	According to the program for cleaning and disinfection device

D-W: Drinking water

FD-W: Fully desalinated water (demineralized, low microbiological contamination: drinking water quality at least)

*Recommended: BBraun Helimatic Cleaner alkaline

► Check visible surfaces for residues after mechanical cleaning/disinfecting.

Aesculap® S⁴® Cervical

S⁴® Cervical Instruments

7.10 Mechanical cleaning/disinfection with manual pre-cleaning

Note

The cleaning and disinfection device must be of tested and approved effectiveness (e.g. FDA approval or CE mark according to DIN EN ISO 15883).

Note

The cleaning and disinfection device used for processing must be serviced and checked at regular intervals.

Manual pre-cleaning with a brush

Phase	Step	T [°C/°F]	t [min]	Conc. [%]	Water quality	Chemical
I	Disinfectant cleaning	RT (cold)	>15	2	D-W	Aldehyde-free, phenol-free, and QUAT-free concentrate, pH ~ 9*
II	Rinsing	RT (cold)	1	-	D-W	-

D-W: Drinking water

RT: Room temperature

*Recommended: BBraun Stabimed

- Note the information on appropriate cleaning brushes and disposable syringes, see Validated cleaning and disinfection procedure.

Phase I

- Fully immerse the product in the cleaning/disinfectant for at least 15 min. Ensure that all accessible surfaces are moistened.
- Clean the product with a suitable cleaning brush in the solution until all discernible residues have been removed from the surface.
- If applicable, brush through non-visible surfaces with an appropriate cleaning brush for at least 1 min.
- Mobilize non-rigid components, such as set screws, links, etc. during cleaning.
- Thoroughly rinse through these components with the cleaning disinfectant solution (at least five times), using a disposable syringe.

Phase II

- Rinse/flush the product thoroughly (all accessible surfaces) under running water.
- Mobilize non-rigid components, such as set screws, joints, etc. during rinsing.

Manual pre-cleaning with ultrasound and brush

Phase	Step	T [°C/°F]	t [min]	Conc. [%]	Water quality	Chemical
I	Ultrasonic cleaning	RT (cold)	>15	2	D-W	Aldehyde-free, phenol-free, and QUAT-free concentrate, pH ~ 9*
II	Rinsing	RT (cold)	1	-	D-W	-

D-W: Drinking water

RT: Room temperature

*Recommended: BBraun Stabimed

- Note the information on appropriate cleaning brushes and disposable syringes, see Validated cleaning and disinfection procedure.

Phase I

- Clean the product in an ultrasonic cleaning bath (frequency 35 kHz) for at least 15 min. Ensure that all accessible surfaces are immersed and acoustic shadows are avoided.
- Clean the product with a suitable cleaning brush in the solution until all discernible residues have been removed from the surface.
- If applicable, brush through non-visible surfaces with an appropriate cleaning brush for at least 1 min.

- Mobilize non-rigid components, such as set screws, links, etc. during cleaning.

- Thoroughly rinse through these components with the cleaning disinfectant solution (at least five times), using a disposable syringe.

Phase II

- Rinse/flush the product thoroughly (all accessible surfaces) under running water.
- Mobilize non-rigid components, such as set screws, joints, etc. during rinsing.

Mechanical alkaline cleaning and thermal disinfecting

Machine type: single-chamber cleaning/disinfection device without ultrasound

Phase	Step	T [°C/°F]	t [min]	Water quality	Chemical
I	Pre-rinse	<25/77	3	D-W	-
II	Cleaning	55/131	10	FD-W	■ Concentrate, alkaline: - pH = 13 - <5 % anionic surfactant ■ 0.5 % working solution - pH = 11*
III	Intermediate rinse	>10/50	1	FD-W	-
IV	Thermal disinfecting	90/194	5	FD-W	-
V	Drying	-	-	-	According to the program for cleaning and disinfection device

D-W: Drinking water

FD-W: Fully desalinated water (demineralized, low microbiological contamination: drinking water quality at least)

*Recommended: BBraun Helimatic Cleaner alkaline

- Check visible surfaces for residues after mechanical cleaning/disinfecting.

Aesculap® S⁴® Cervical

S⁴® Cervical Instruments

7.11 Inspection, maintenance and checks

CAUTION

Damage (metal seizure/friction corrosion) to the product caused by insufficient lubrication!

- Prior to function checks, lubricate moving parts (e.g. joints, pusher components and threaded rods) with maintenance oil suitable for the respective sterilization process (e.g. for steam sterilization: Aesculap STERILIT® I oil spray JG600 or STERILIT® I drip lubricator JG598).

- Allow the product to cool down to room temperature.
- After each complete cleaning, disinfecting and drying cycle, check that the instrument is dry, clean, operational, and free of damage (e.g. broken insulation or corroded, loose, bent, broken, cracked, worn, or fractured components).
- Dry the product if it is wet or damp.
- Repeat cleaning and disinfection of products that still show impurities or contamination.
- Check that the product functions correctly.
- Immediately put aside damaged or inoperative products and send them to Aesculap Technical Service, see Technical Service.
- Assemble dismountable products, see Assembling.
- Check for compatibility with associated products.

7.12 Packaging

- Appropriately protect products with fine working tips.
- Store products with ratchet locks fully opened or locked no further than in the first notch.
- Place the product in its holder or on a suitable tray. Ensure that all cutting edges are protected.
- Pack trays appropriately for the intended sterilization process (e.g. in sterile Aesculap containers).
- Ensure that the packaging provides sufficient protection against recontamination of the product during storage.

7.13 Steam sterilization

Note

The product can be sterilized either in disassembled or in assembled condition.

Note

To avoid breakage due to stress crack corrosion, sterilize the products with the lock fully open or locked no further than on the first ratchet tooth.

- Check to ensure that the sterilizing agent will come into contact with all external and internal surfaces (e.g. by opening any valves and faucets).
- Validated sterilization process
 - Steam sterilization through fractionated vacuum process
 - Steam sterilizer according to DIN EN 285 and validated according to DIN EN ISO 17665
 - Sterilization using fractionated vacuum process at 134 °C/holding time 5 min
- When sterilizing several instruments at the same time in a steam sterilizer, ensure that the maximum load capacity of the steam sterilizer specified by the manufacturer is not exceeded.

7.14 Sterilization for the US market

- Aesculap advises against sterilizing the device by flash sterilization or chemical sterilization.
- Sterilization may be accomplished by a standard prevacuum cycle in a steam autoclave.

To achieve a sterility assurance level of 10⁻⁶, Aesculap recommends the following parameters:

Aesculap Orga Tray/Sterile container (perforated bottom) Minimum cycle parameters*			
Sterilization method	Temp.	Time	Minimum drying time
Prevacuum	270 °F/275 °F	4 min	20 min

*Aesculap has validated the above sterilization cycle and has the data on file. The validation was accomplished in an Aesculap sterile container cleared by FDA for the sterilization and storage of these products. Other sterilization cycles may also be suitable, however individuals or hospitals not using the recommended method are advised to validate any alternative method using appropriate laboratory techniques. Use an FDA cleared accessory to maintain sterility after processing, such as a wrap, pouch, etc.

WARNING for the US market

If this device is/was used in a patient with, or suspected of having Creutzfeldt-Jakob Disease (CJD), the device cannot be reused and must be destroyed due to the inability to reprocess or sterilize to eliminate the risk of crosscontamination.

7.15 Storage

- ▶ Store sterile products in germ-proof packaging, protected from dust, in a dry, dark, temperature-controlled area.
- ▶ Store sterile packed single-use products dust-protected in a dry, dark and temperature-controlled room.

8. Technical Service

Risk of injury and malfunction!

- ▶ **Do not modify the product.**

- ▶ For service and repairs, please contact your national B. Braun/Aesculap agency.

Modifications carried out on medical technical equipment may result in loss of guarantee/warranty rights and forfeiture of applicable licenses.

Service addresses

Aesculap Technischer Service

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78532 Tuttlingen / Germany

Phone: +49 (7461) 95-1602

Fax: +49 (7461) 16-5621

E-Mail: ats@aesculap.de

Or in the US:

Aesculap Implant Systems LLC

Attn. Aesculap Technical Services

615 Lambert Pointe Drive

Hazelwood

MO, 63042

Aesculap Repair Hotline

Phone: +1 (800) 214-3392

Fax: +1 (314) 895-4420

Other service addresses can be obtained from the address indicated above.

9. Disposal

- ▶ Adhere to national regulations when disposing of or recycling the product, its components and its packaging!

10. Distributor in the US/Contact in Canada for product information and complaints

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