



POWER SYSTEMS

AESCULAP® Acculan 4 Sternum Saws

LOOK. FEEL. EXPERIENCE.

AESCULAP® Acculan 4 Sternum Saws

LOOK. FEEL. EXPERIENCE.



Acculan

The Acculan 4 Sternum saws are powerful and simple to use. They combine reliability, longevity and ergonomics.

THE DIFFERENCE IS IN THE DETAIL

AESCULAP® was a pioneer in the introduction of non-sterile batteries. This avoids sterilization not only of the temperature-sensitive battery cells, but also of the equally vulnerable integrated control electronics. This, in turn, has a positive effect on the service life of the batteries and, therefore, on the system's operating costs. A concept that has been proven for years now. The opportunity to replace the battery at a fixed price and the further use of the electronics and housing helps to preserve the environment as well as to deliver financial benefits.

MECHANICAL PROCESSING


Thanks to the titanium housing and the ECCOS storage solution Acculan 4 is suitable for mechanical cleaning, even with alkaline cleaning media.

SERVICE LABEL

Regular servicing and preventive maintenance is essential to ensure the reliability of motor systems over a long service life.

The service label on each Acculan 4 machine clearly shows when the next preventive maintenance is due.

Acculan 4

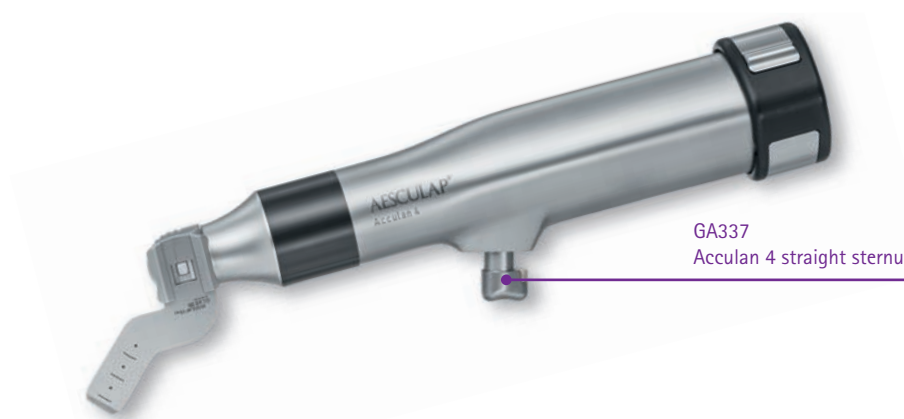

GA337
Acculan 4 straight sternum saw

- Accessories, such as the ECCOS storage solution, can be found in the Power Systems main catalogue 084002.
- The tools can be found in the Burrs & Blades catalogue 017599.

AESCULAP® – a B. Braun brand

Aesculap AG | Am Aesculap-Platz | 78532 Tuttlingen | Germany
Phone +49 7461 95-0 | Fax +49 7461 95-2600 | www.aesculap.com

The main product trademark "Aesculap" is a registered trademark of Aesculap AG.

Subject to technical changes. All rights reserved. This brochure may only be used for the exclusive purpose of obtaining information about our products. Reproduction in any form partial or otherwise is not permitted.