

Aesculap Neurosurgery

YASARGIL Aneurysm Clip System

Double Bayonet Clip Applier



Double Bayonet Clip Applier

Product features

- The applier has a double bayonet design, enabling an excellent view of the clip jaws.
- The applier arms are very slim, occupying significantly less space.
- A new feature is the uniform colour coding by means of a coloured button near the handle.
- The appliers are available either with or without latch.

For appliers with latch

- The lock has been relocated to the rear part of the handle, reducing the danger of accidental actuation by the neurosurgeon during clip application.
- The latch can be actively engaged by pressing a button, allowing the surgeon to decide, whether or not he wishes to use it.





■ for YASARGIL aneurysm clips


Titanium

Phynox

Double Bayonet Clip Appliers

| | | |
|-------------|-------------|---------------|
| ● FT501T | ● FE301K | with latch |
| FT505T | FE305K | without latch |

- Standard
- Total length 225 mm

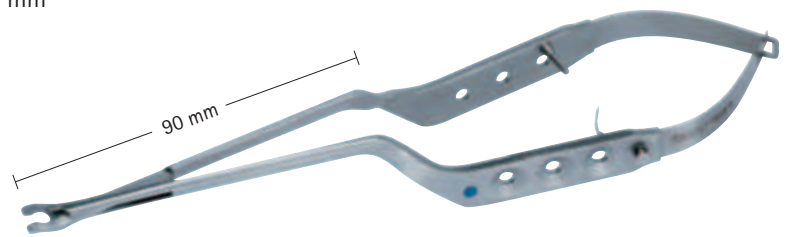


| | | |
|-------------|-------------|---------------|
| ● FT502T | ● FE302K | with latch |
| FT506T | FE306K | without latch |

- Standard
- Total length 245 mm




| | | |
|-------------|-------------|---------------|
| ● FT511T | ● FE311K | with latch |
| FT515T | FE315K | without latch |

- Miniature
- Total length 225 mm

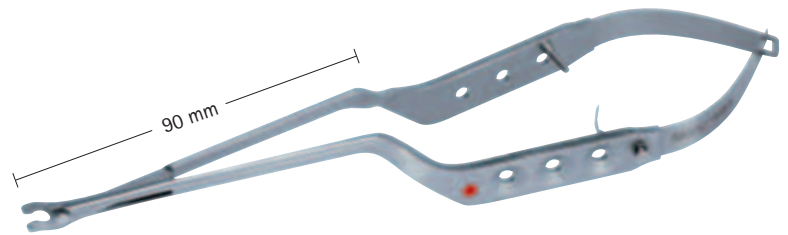

| | | |
|-------------|-------------|---------------|
| ● FT512T | ● FE312K | with latch |
| FT516T | FE316K | without latch |

- Miniature
- Total length 245 mm

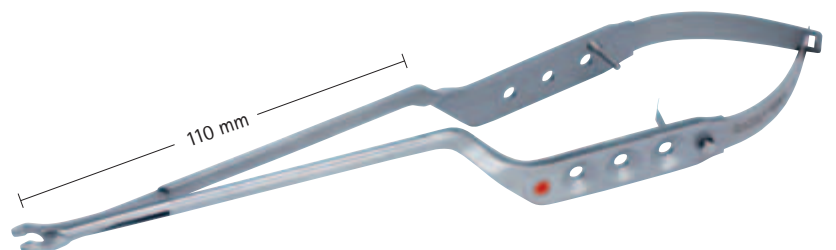



AESCULAP®

B | BRAUN
SHARING EXPERTISE

Aesculap AG

Am Aesculap-Platz
78532 Tuttlingen
Germany

Phone +49 7461 95-0
Fax +49 7461 95-2600

www.aesculap.de

All rights reserved. Technical alterations are possible. This leaflet may be used for no other purposes than offering, buying and selling of our products. No part may be copied or reproduced in any form. In the case of misuse we retain the rights to recall our catalogues and pricelists and to take legal actions.