


# Haemofix

Pressure Monitoring-Set for  
Closed Blood Sampling

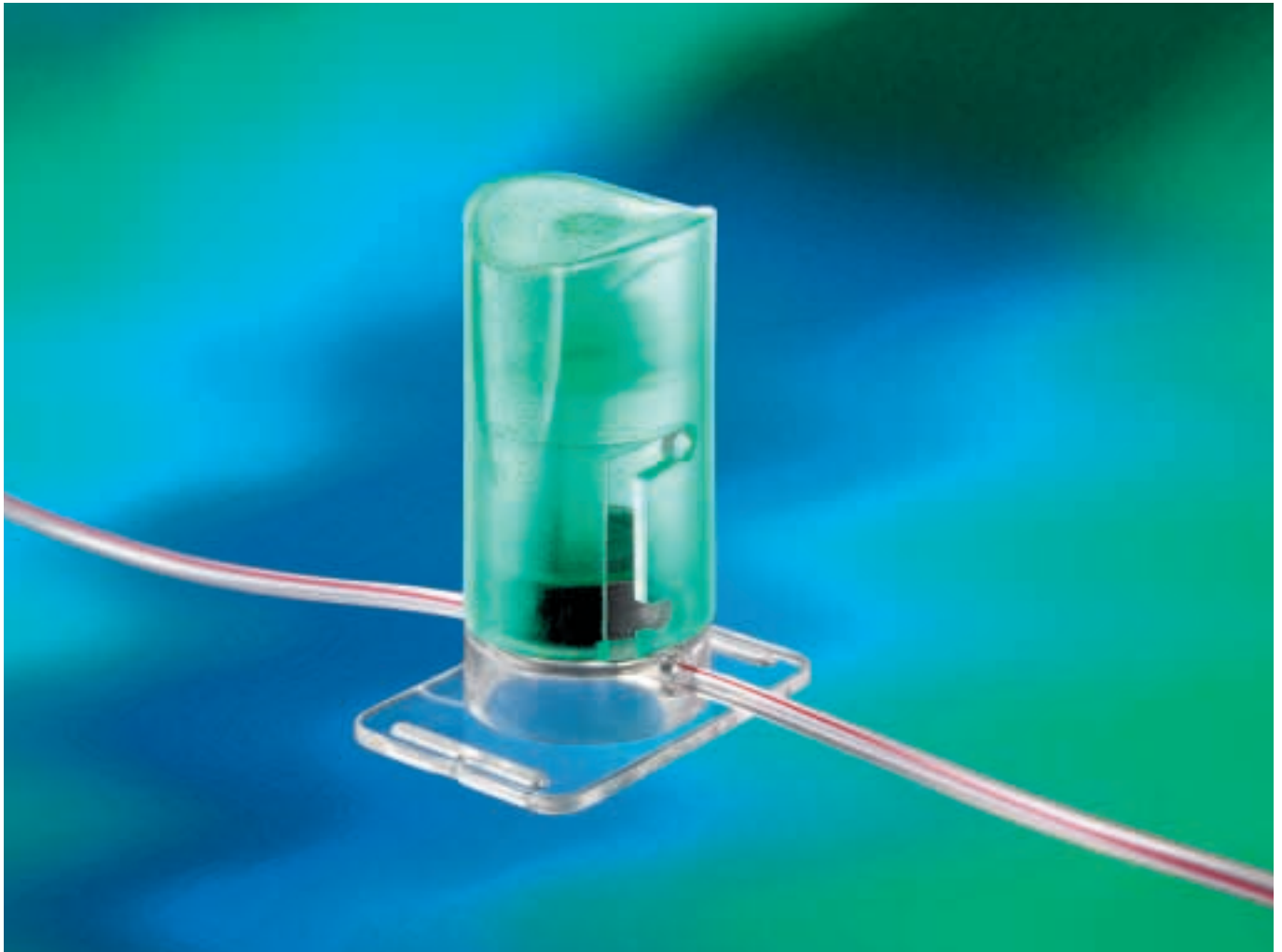

## Pressure Monitoring-Set for Closed Blood Sampling

### Product Features



#### Safety and accuracy in Blood Gas Analysis

- The reservoir is designed for safe blood sampling through arterial and right-heart catheters.
- The integrated spring mechanism allows automatic blood aspiration without affecting the quality of blood gas analysis. The system can also be used with patients with low blood pressure.
- External contamination of the aspirated volume is avoided by a labyrinth-seal between the cap and the housing.
- The design of the fluid channel allows easy air bubble-free filling of the Monitoring-Set.

#### Ease of use

- The reservoir can be handled comfortably with just one hand.
- Accurate blood gas analysis is achieved without blood wastage
- The reservoir can be filled and emptied completely free of dead space.
- It is not necessary to handle On/Off taps due to the integrated lock mechanism in the reservoir.

■ The reservoir and the port have no noticeable effect on resonance frequency.

■ The combination of the blood sampling port together with adapters designed to meet your needs provides a complete system for closed blood sampling.

### Product Summary

Description	Article	Sales Unit
Haemofix-Combitrans Monitoring-Set arterial (150 cm)	05213505	20
Haemofix-Exadyn Monitoring-Set arterial and venous (150 cm)	05213516	20
Haemofix-Combitrans Monitoring-Set double (150 cm)	05213527	20
Sampling Adapter Luer-Lock	05213940	50
Sampling Adapter Vacutainer	05213835	50
Combitrans attachment plate (4-fold)	05203597	1
IV-Pole-Mount	04086007	1

Custom devices on request. Please ask your clinical specialist!

CE 0123

**B | BRAUN**  
HOSPITAL CARE

Cardio-Medical

**B. Braun Melsungen AG**  
Sieversufer 8  
D-12359 Berlin

Tel (0 30) 68 98 97-0  
Fax (0 30) 68 98 97-30  
[www.bbraun.com](http://www.bbraun.com)